€ 19.-

The Guide to Imaging Technology and Informatics in Europe | Columbia | Colum



NEVER STOP

ADVANCING HEALTHCARE

When you're a company who has spent over 80 years striving to move the world forward, you never stop. We'll never stop innovating in the medical areas of prevention, diagnosis and treatment.

With technologies like AI powered imaging solutions that are Designed to help improve the accuracy of diagnoses and surgery. Applying our scientific expertise to advance regenerative medicine to meet unmet medical needs.

Creating nano-particle based drug delivery technologies to help in the fight against cancer, or continuing to pioneer endoscopic and x-ray imaging technologies.

We'll never stop improving the future, creating value from innovation that will change the world of healthcare forever.

FUJ!FILM

Value from Innovation

FUJIFILM and Fujifilm Value from Innovation are trademarks of FUJIFILM Corporation. © 2018 FUJIFILM Corporation. All rights reserved.



HITACHI Inspire the Next

Dear Reader,

it was a bit more than three years ago, in November 2016, that Geoff Hinton at an Al conference in Toronto demanded we stop training radiologists since quite obviously the profession would be obsolete in five or maybe ten years' time and after all there were plenty of radiologists around today.

Now, three years later, there are not plenty of radiologists around – quite the contrary – and Al has not pushed a single radiologist out of a job. Moreover, I dare say that Geoff Hinton would not trust a diagnosis that was made 100 percent by machines. Patients continue to prefer placing their trust in human physicians whose treatment decisions are based on a mix of experience and current research and whose diagnoses are supported by state-of-the-art technology.

Al, no doubt, is gaining ground and machine learning algorithms will be used to automatize error-prone and repetitive tasks.

These algorithms are particularly suited to optimise radiology workflows, to automatically select protocols and settings, analyse images, suggest diagnoses and prepare structured reports. The responsibility for well-founded diagnoses, however, will remain with the radiologist for years to come. To trust systems that self-validate their results and deduct findings from these results might just be a tad naïve.

That's why we will continue to need radiologists! Physicians who – aware of their knowledge and their experience – face the new technology with self-confidence will open up new perspectives in medicine.

Your editorial team

Daniela Zimmermann and Guido Gebhardt

the next level in humancentered design

Titachi Medical Systems Europe Holding AG, Switzerland www.hitachi-medical-systems.com





www.VarexImaging.com | www.Veolity.com



Products & Solutions

Computed Tomography	6
Dual Source CT	7
Volume CT	7
20 to 64 Slices	12
2 to 16 Slices	13
Mammo CT	15
Conebeam CT	15
Oncology CT	15
Accessories / Complementary Systems	16
Magnetic Resonance Imaging	24
PET MR	25
7 Tesla	25
3 Tesla	25
1.5 Tesla	26
Oncology	28
Open	28
MRI Coils	29
Accessories / Complementary Systems	30
Injectors	32
Accessories / Complementary Systems	36
Interventional Systems	37
Interventional Systems Multi Modality Suites	
· · · · · · · · · · · · · · · · · · ·	38
Multi Modality Suites Bi-Plane Single Plane	38 38 39
Multi Modality Suites Bi-Plane	38 38 39
Multi Modality Suites Bi-Plane Single Plane	
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms	38 38 39 42 44
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms	38 38 39 42 44
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories / Complementary Systems	38 38 39 42 44 48
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories / Complementary Systems Artificial Intelligence IT Systems	38 38 39 42 44 48 49
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories/Complementary Systems Artificial Intelligence IT Systems RIS/PACS	38 38 39 42 44 48 49 52
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories / Complementary Systems Artificial Intelligence IT Systems	38 38 39 42 44 48 49 52
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories/Complementary Systems Artificial Intelligence IT Systems RIS/PACS. IT-Solutions Overview	38 38 39 42 44 48 49 52 53 54 59
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories / Complementary Systems Artificial Intelligence IT Systems RIS / PACS IT-Solutions Overview Business Intelligence	38 38 39 42 44 48 49 52 53 54 59 60
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms. Surgical II-C-Arms. Accessories / Complementary Systems Artificial Intelligence IT Systems RIS / PACS. IT-Solutions Overview Business Intelligence. Mammo Workstations.	38 38 39 42 44 48 49 52 53 54 59 60 61
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms. Surgical II-C-Arms. Accessories / Complementary Systems Artificial Intelligence IT Systems RIS / PACS. IT-Solutions Overview Business Intelligence. Mammo Workstations. Pathology.	38 38 39 42 44 48 49 52 53 54 59 60 61
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms Surgical II-C-Arms Accessories / Complementary Systems Artificial Intelligence IT Systems RIS / PACS IT-Solutions Overview Business Intelligence Mammo Workstations Pathology Portal Solutions	38 38 39 42 44 48 49 52 53 54 60 61 61 62
Multi Modality Suites Bi-Plane Single Plane Surgical Flat Panel C-Arms. Surgical II-C-Arms. Accessories / Complementary Systems Artificial Intelligence IT Systems RIS / PACS IT-Solutions Overview Business Intelligence. Mammo Workstations. Pathology. Portal Solutions Mobile RIS / PACS Viewers.	38 38 39 42 44 48 49 52 53 54 60 61 61 62 63


Mammography	66
Mammo CT	67
Tomosynthesis	
Digital Mammography	
Mammo DR	
Biopsy Units	70
Film-Screen Mammography	70
Accessories / Complementary Systems	70
R/F Film-Screen	72
Bucky	73
Fluoroscopy	
Mobile X-Ray	
Accessories / Complementary Systems	78
R/F Digital	80
CR	82
DR	83
Portable DR	
DR Detectors	96
Mobile DR	101
Flatpanel Fluoro	107
DXA	111
Accessories / Complementary Systems	112
Molecular Imaging	115
PET-MR	116
PET-CT	116
SPECT-CT	117
SPECT	117
Displays / Printers	118
Displays – Mammo	119
Displays – Color	
Displays – Grayscale	
DVD Import	
DVD Burner	
Printers	121
Ultrasound	123
Accessories / Complementary Systems	138
Testing Devices	139


5

Computed Tomography



Dual Source CT

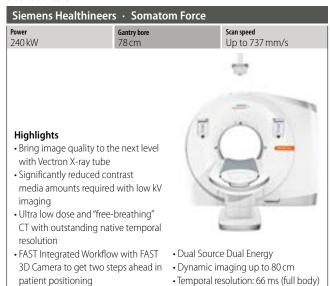

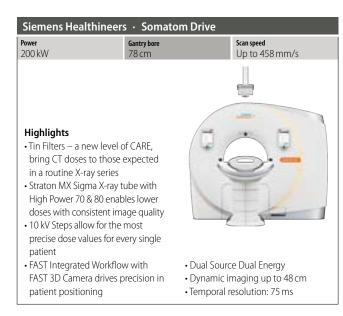

Volume CT

• 0.17 mm spatial resolution

• SURE Position patient centring

• 640 slices/rotation

• 0.17 mm spatial resolution



Isophasic organ perfusion

CT Fluoroscopy

Energy applications


Isophasic organ perfusion






Canon · Aquilion One Genesis Edition 350 Gantry bore Scan range 150/200 cm Power 100 kW Highlights Advanced Intelligent Clear-IQ Engine (AiCE) Artificial intelligence in scanning • SURE Position patient centring and processing · Lateral table movement • PUREVISION Optics imaging chain • 300 kg patient table • 0.35 s rotation • SEMAR (Metal Artefact Reduction) • 16 cm coverage • Wide range of Cardiac and Dual · 640 slices/rotation **Energy applications** • 0.17 mm spatial resolution Isophasic organ perfusion

Volume CT

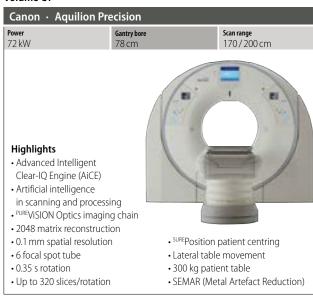

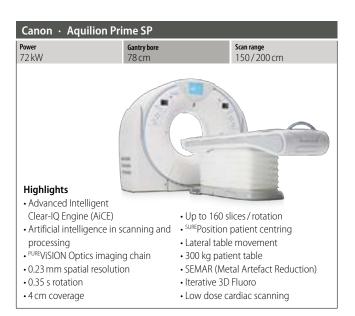







(total: 200 mm) • 650 mm wide patient table with weight limit of 250 kg • Slices per rotation 64/128 • Dual Energy Scan

Scan range

200 cm

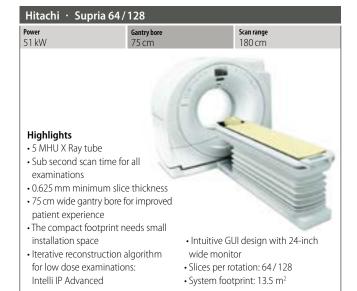

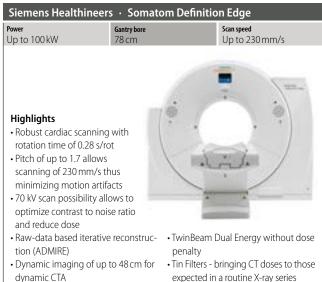

Dunlee is your CT tube partner now and in the future

Dunlee, a supplier of reliable medical imaging components, announced at RSNA 2019 its investment in the development of a new CT replacement tube with liquid metal bearing. With this investment, customers can be confident that Dunlee will offer high-quality alternatives to OEM tubes for current and future generations of GE scanners.

Customer testimonials on early installations planned for RSNA 2020

Since RSNA, Dunlee has reached the next level of product development, and plans to share the first customer testimonial cases of the new tube at RSNA 2020. "We've built our business on helping our customers secure a prosperous future," says Alexander Eitel, Head of Marketing & Business Development., Because we have loyal, long-term relationships with them, as well as many new prospects, we are continually investing in new products and anticipating future needs so we can provide a robust portfolio of alternatives to OEM tubes."

Dunlee's full line of CT replacement tubes offers reliability, quality and economic value, which helps minimize CT scanner operational costs and downtime. Dunlee also provides 24/7/365 technical support to help customers efficiently replace, install and calibrate tubes.

Dunlee is a leading provider of CT tubes for the replacement market, as well as CT, X-ray and MR solutions for the OEM market.

www.dunlee.com





Intelligent CT Imaging

How can intelligent and personalized CT scanning ease radiology workflow?

Increasingly, hospitals are aiming for more than just precision, quality, and speed in radiology. They want to go one step further – a step the University Hospitals of Erlangen and Navarra have now taken.

Francisco Cubo was desperate: The 52-year-old Spaniard was diagnosed with cancer of the kidney and malignant lung nodules. The real estate broker from Madrid sought out the best possible medical care and found the Clínica Universidad de Navarra. This clinic in Pamplona, northern Spain, is among the country's most prestigious referral hospitals for treatment of cancer and heart diseases. His kidney cancer was successfully dealt with. But recently, new metastases have formed in his lungs, which are now being treated

and checked regularly. The journey from Madrid is long, but for him, it's worth the effort.

At 10:30 a.m., the computed tomography (CT) scan is performed. By 11:20, Francisco has already received the good news: The metastases have not spread any further. The treatment in the test phase seems to be effective. "In a regular hospital, I would have waited a week for the result. That's an eternity when you're in fear of your life. Here, I get the result within 50 minutes. That's priceless," says Francisco.



Seeking quality and personalized treatment

Professional excellence combined with short waiting times and fast, personalized treatment are the hallmarks of this university clinic, which attracts patients from all over Spain, other European countries and even Africa.

"To be able to work both rapidly and with high precision, we need the latest scanning technology," explains Gorka Bastarrika, MD, Head of Radiology.

Radiology nurse Begoña Sara finds that especially corpulent or less mobile patients can be positioned more easily thanks to the bigger gantry opening.



"SOMATOM X.cite is not just one of the most versatile devices on the market, but it is also the first to facilitate personalized, individualized scanning. As far as I'm concerned, this machine marks the beginning of a new era in computed tomography," says Bastarrika.

Matthias May, MD, Assistant Professor at the radiology department of University Hospital Erlangen, Germany, which also took part in the SOMATOM X.cite pilot study, agrees: "The scanner is definitely the first step toward personalized radiology."

"By using the predefined decision-making criteria, even less experienced staff are able to scan rapidly and flawlessly with the optimal image quality, X-ray dose, and contrast medium dose," says the radiologist. This, he believes, is the upcoming trend in CT imaging. "We envision a future where you won't need to be an expert to operate a CT scanner. This will be crucial not only in regions lacking a radiologist, but even here in our hospital. For instance, if less experienced personnel find themselves in stressful emergency situations or alone during their night shift."

Since August, Bastarrika's radiology team has been using the new SOMATOM® X.cite from Siemens Healthineers.

"Thanks to intelligent support from the new user concept myExam Companion as well as the automated capture of imaging and reconstruction settings, less experienced personnel can also carry out even the most complex scans quickly and accurately," says the radiology chief. He notes that especially with cardiology patients suffering from arrhythmias, arteriosclerosis, or high heart rates, it is now much easier to standardize settings through intelligent decision trees.

A new era in computed tomography

The new CT scanner also offers advantages for patients. The bigger gantry opening means larger and less mobile patients, or those attached to medical devices, can be positioned more easily, says radiology nurse Begoña Sara. For patient Francisco Cubo, the biggest benefit is that he can understand the breathing instructions without his hearing aid. "The new visual color codes indicate when you must breathe in or hold your breath, which makes it much easier for me," he explains.


Huge time savings in positioning


On her tablet, radiology technologist Laura Schwarzfaerber clicks through the simple questions that drive the decision criteria of myExam Companion together with the patient input. Automated algorithms read the patient's ECG and feed it directly into the user system, which calculates the appropriate protocols and doses of radiation adjusted for ECG data, weight, size, and age. Schwarzfaerber stands beside the patient, who is already lying on the scanner table. The fully automated positioning program, supported by a 3D camera, places the patient in the optimal position. "We save loads of time at this stage, since all of these steps had to be done manually before," says Schwarzfaerber. With the advanced tablet functionalities, she can now remain at the patient's side even longer – "and that makes them much more calm and relaxed, which also assists our scanning process."

siemens-healthineers.com/somatom-xcite

The statements by Siemens Healthineers customers described herein are based on results that were achieved in the customer's unique setting. Since there is no "typical" hospital and many variables exist (e.g., hospital size, case mix, level of IT adoption) there can be no guarantee that other customers will achieve the same results.

Volume CT

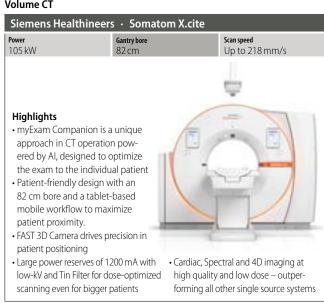


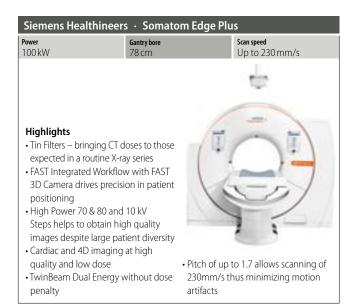







20 to 64 Slices

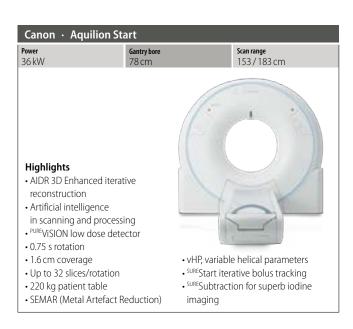





20 to 64 Slices

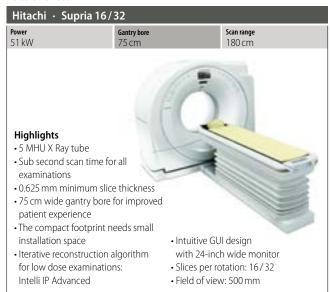

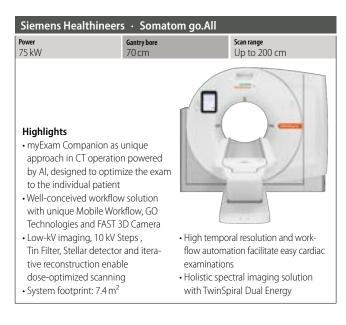

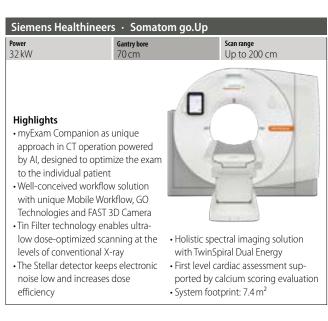

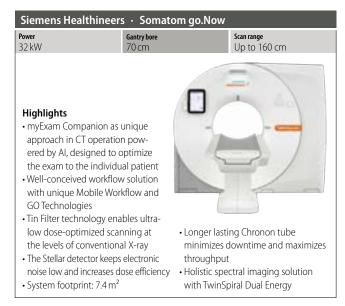

2 to 16 Slices

• Compact footprint to maximize

easiness of installation




• Slices per rotation: 16/32

• System Footprint: 13.5 m²

SternMed · Cytom 16 Power 60 kW Gantry bore 70 cm Scan range 160 cm Highlights · Ultra-fast, Non-contact PowerLink gantry technology

- · Slices per rotation: 16
- · Slice thickness: 0.625, 1.25, 2.5, 5.0, 10 mm
- · Lower dose, lower consumption
- Ultra-fast rare earth ceramic detector
- 24 rows detector, one-disc rotation collects 16 slices images
- 21,504 detector units
- Coverage of each single scan is 20 mm
- Multiple X-ray dose reduction technology
- Intelligent scattering line automatic elimination technology • 120 second non-stop scanning
- 13 RADBOOK 2020

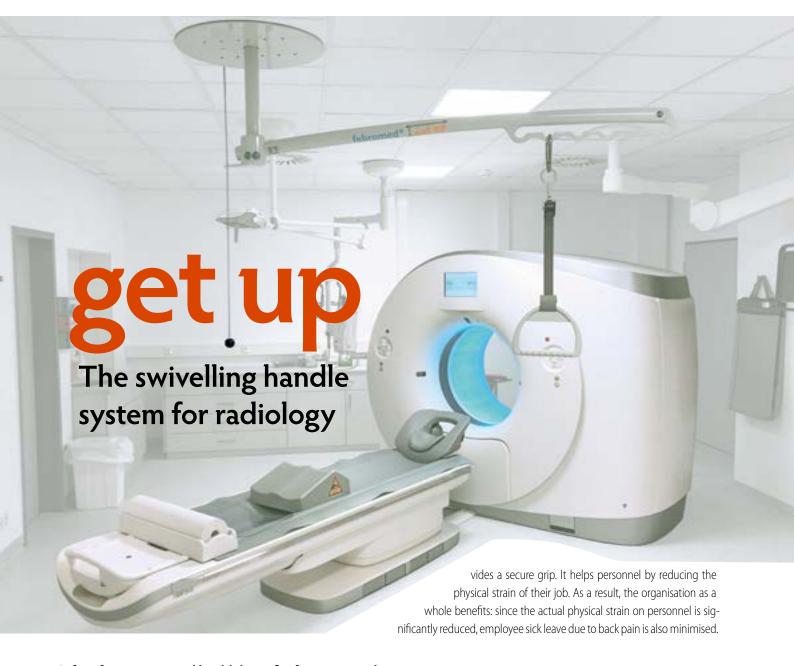

Safety for patients and health benefits for personnel

Febromed GmbH & Co. KG, the expert in delivery room equipment and medical accessories from Oelde, Germany, has developed "get up", an innovative handle system for radiology. The new swivelling system was installed for the first time in a state-of-the-art CT scan room at the Institute of Diagnostic and Interventional Radiology and Neuroradiology at Essen University Hospital.

For a secure grip

Many patients find getting onto the examination table for a CT scan difficult. In particular, restricted mobility leads to uncertainty as the patient is positioned and arranged, thus placing increased physical strain on care personnel, predominately in the back area. The new "get up" handle system from Febromed offers a solution: this swivelling system helps patients get onto the table before their scan and stand up again safely and comfortably afterwards. It minimises the risk of falling and pro-

Positive experiences

After installing the handle system in May 2017, the Institute of Diagnostic and Interventional Radiology and Neuroradiology at Essen University Hospital has consistently had positive experiences. As Anton S. Quinsten, Ltd. MTRA, reports, "We are really happy with the "get up" system from Febromed. The first few months have shown that the handle system is considered a real asset by both patients and personnel."

Space-saving and durable

The 'get up' handle system is designed for space-saving mounting on the ceiling and can be swivelled by 360°. The structure can be locked in 15° increments so that the system is always in the optimal position for the patient. This purely mechanical construction ensures easy handling and extended durability.

www.febromed.com

Mammo CT

grams

compression

· Low dose in the range of mammo-

• Excellent patient comfort without



mastodynia, implants

• Unique single photon counting

detector for highest sensitivity, accuracy and speed (direct conversion)





Conebeam CT


Cefla · NewTom 5G XL		
FOV 6×6cm-21×19cm	Scan time 18 s	Pixel size 100 µm
Highlights Cone Beam CT with oper and supine position. Backside access available High definition volumetr bone tissues, non-overla and fewer artifacts Safe Beam: Automatic de necessary Dose, Pulsed e Extensive range of discip Otorhinolaringology, Oral and Maxillofacial sur "Cine X" Dinamic acquisit	e: ic images of pping sections etection minimal mission lines in Orthopaedics,	isition



Oncology CT



Oncology CT

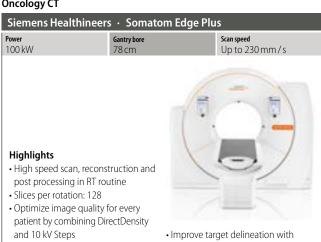

- Precise reduction of metal artifacts TwinBeam and Dual Spiral Dual with iMAR Energy scans
 - Flexible CT imaging with Edge Plus Sliding Gantry

Siemens Healthineers · Somatom Confidence RT Pro		
	Gantry bore 80 cm	Scan speed Up to 96 mm/s
HighlightsDirectDensity allows image any kV value to be used in		
TPS without the need of s calibration curves	everal	

- Enable visualization up to 80 cm bore size with HD FoV Pro
- · Precise reduction of metal artifacts with iMAR
- Improve visualization with Dual Spiral **Dual Energy**
- Leverage your ability to be more accurate with iterative reconstruction

Scan speed

Up to 200 mm/s

- · Slices per rotation: 20/64
- Detector: Stellar RT
- · Dual Spiral Dual Energy

Siemens Healthineers · Somatom go. Open Pro Power 75 kW Gantry bore Scan speed Up to 200 mm/s Highlights

• Direct i4D: First 4D CT scan mode to adapt to breathing patterns in real time for dramatic motion artifact reduction

· Comprehensive motion manage-

ment solutions for moving tumors

- 4 cm detector coverage and 0.35 s rotation times for deep inspiration breath-hold scanning
- DirectORGANS: Al powered organs-at risk contouring directly at the CT console for advanced contouring results
- TwinSpiral Dual Energy scanning and Tin filter for less variability in target contouring



- Direct Laser: Seamless integration of patient marking lasers and laser QA for time saving and error avoidance
- · Large bore of 85cm with 60cm true scan Field of View; recon. slices per rotation: 128

Siemens Healthineers · Somatom go.Sim

Gantry bore

Highlights

Power 75 kW

- DirectORGANS: Al powered organs-at risk contouring directly at the CT console for consistent results
- · Direct Laser: Seamless integration of patient marking lasers and laser QA for time saving and error avoidance
- · Mobile Workflow: Re-designed workflows with mobile tablet and SIM&GO technologies to increase efficiency and patient satisfaction
- · DirectDensity: Personalized scanning and scan kV from 70 – 140kV (in steps of 10)



- 227 / 307 kg Patient table (fully TG-66 compliant) with new flat table top
- Large bore of 85 cm with 60 cm true scan Field of View; recon. slices per rotation: 64



Scan speed Power 32 kW Gantry bore Up to 42 mm/s

Highlights

- The mobile workflow brings an unparalleled level of flexibility and mobility to daily CT routines
- GO technologies help standardize and simplify the simulation process
- Stellar detector keeps electronic noise low, increases dose efficiency, and improves spatial resolution
- Tin Filter (Sn) cuts out lower energies to reduce dose and optimizes contrast at the interface between soft tissue and air



- The dedicated Patient Table RT is fully compliant with AAPM TG-66, compatible
- · Guide&GO is an intuitive tablet-based solution for CT-guided interventions

Accessories / Complementary Systems

Canon Electron Tubes & Devices · LM-CT Tube



Highlights

- For CT systems (2-MHU to 4-MHU)
- · Uses a liquid metal bearing
- Supports 0.5 s full scans
- Our unique liquid metal bearing technology uses an all-metal target, enabling high anode heat dissipation with low noise and long bearing life.

Accessories / Complementary Systems

Dunlee · CT Replacement Tubes



Highlights

Dunlee's CT replacement tubes:

- Meticulously engineered to be compatible with a variety of popular GE scanners
- · Offer excellent quality
- Include 24/7/365 technical support
- Tube stocks at major airport hubs in the United States, Europe and Asia

Dunlee · Components Bundle for CT



Highlights

- The Xpert bundle with CT8000 meets your needs for an ideal integration into high-end CT systems
- Consists of an X-ray tube, generator, cooling unit and cables
- For a fast integration into your development process to speed up the time to market
- Extraordinary image quality combined with wide coverage

NEWTOM 5G XL EXTRA.VISION

www.newtom.it

CUTTING EDGE PERFORMANCE FOR ADVANCED DIAGNOSTICS.

Advanced diagnostics with 5G XL, the only CBCT device with lying down patient positioning that offers excellent stabilisation and a broad range of FOVs for very high quality 3D and 2D images.

5G XL is the first device with native FOV 21 x 19 cm for highly detailed investigations. The ideal device that produces clear, high definition images for orthopaedics, otorhinolaryngology, maxillofacial surgery and dentistry applications.

SUPERIOR DIAGNOSTIC QUALITY Very high resolution 2D and 3D images with a broad range of FOV for an extensive selection of clinical applications.

OPTIMAL LYING DOWN POSITION The only CBCT system with patient lying down positioning, a motor-driven patient table and an open gantry. Perfect patient stabilisation considerably reduces any

movement-induced artifacts.

LOW X-RAY DOSE The ECO Scan mode and SafeBeam[™] technology offered by

and SafeBeamTM technology offered by 5G XL further reduce the dose radiated to the patient, compared to examinations performed with CT technology.

SPECIALIST SOFTWARE The adaptive user interface displays images and provides access to innovative 3D and 2D analysis functions for rapid and accurate diagnoses and optimal workflow.

ECR 2020 · Vienna · AUSTRIA 15-19 JULY EXPO X1 BOOTH N. 111

New Tome what's next

The power of the new Aquilion ONE / PRISM Edition

Since its global launch in December 2019, the Aquilion ONE / PRISM Edition has been making waves all around the world particularly with its Advanced intelligent Clear-IQ Engine (AiCE), Deep Learning Spectral and CT fluoroscopy capabilities. To find out more about what our customers think, we sat down with leading clinicians from two-prestigious institutions, Radboud University Medical Center (the Netherlands) and Strasbourg University Hospital (France), for an in-depth review of the CT system that's designed for deep intelligence.

One of the most exciting features of the PRISM is our Advanced intelligent Clear-IQ Engine (AiCE), a cutting-edge Deep Learning Reconstruction technology that's been trained to reduce noise and boost signal to deliver sharp, clear and distinct images at speed.

Can you tell us a bit about how your teams have adopted AiCE?

Prof. Ohana: We weren't sure in the beginning whether to use AiCE for specific applications or as an additional reconstruction, so we tried it with the first patients. The difference with AIDR 3D was so immense that we decided, ok, let's use it for every patient.





"AiCE lets us obtain thinner images at high quality, which helps us make a more confident diagnosis."

Prof. Mathias Prokop

Head of Radiology, Radboud University Medical Center, the Netherlands

Prof. Prokop: Our people accepted AiCE extremely quickly and have used it immediately for most applications.

What's the most significant difference that AiCE has made to your practice?

Prof. Prokop: AiCE lets us obtain thinner images at high quality, which helps us make a more confident diagnosis. When you're looking at complex cases in which you have to focus on detail, your confidence goes up.



Dr. Smit: I really like AiCE because of its low contrast detectability, especially with liver lesions or peritoneal metastases. Also, the contours are much sharper, which helps with the delineation of bowels, to assess fat stranding, or to see if a tumour is infiltrating another tissue.

What about dose, have you seen any benefits there?

Dr. Smit: Because AiCE in integrated in the automatic exposure controle, there's no guessing on how low the dose should be. It's just about switching on the button for dose reduction. I generally check the dose report for every scan, and I'm often surprised at how low the dose is.

Prof. Prokop: We wanted to move on from using morphological information, like in the past, to using more functional information, and we wanted to do that using the lowest possible dose and in a repeatable, reproducible fashion. AiCE lets us do that.

"I really like AiCE because of its low contrast detectability and the much sharper contours."

Dr. Ewoud Smit

Radiologist, Radboud University Medical Center, the Netherlands



"The image quality is excellent with the new PRISM Edition."

Prof. Catherine Roy

Radiologist, Strasbourg University Hospital, France

Spectral Imaging: Improved quality and workflow

In addition to AiCE, the new Aquilion ONE / PRISM Edition features advanced spectral imaging to improve workflow and help reduce patient iodine load. Images can be analysed directly at the workstation with a range of new Vitrea applications.

Have you seen an overall reduction in iodine load?

Prof. Roy: That's the main application in our department and all around the world. Spectral imaging means that we're now able to use a contrast agent of 270 mg I/mL instead of staying in the 400 mg I/mL range. Plus, you can still see everything perfectly, with all the tiny details. This is a user-point that still amazes me.

Have you noticed an improvement in tissue classification?

Prof. Ohana: In oncology imaging, we sometimes have subtle lesions. With spectral imaging, we are more confident and can see if there is or isn't contrast uptake in these lesions. Spectral imaging is also more specific, mostly for gout or kidney stones. We are able to characterize the tissue and give more specific answers to the clinician.

Radiology Team Strasbourg University Hospital, France

What about workflow, is that better?

Prof. Ohana: Acquisition time is the same as a single energy examination; there's no difference in terms of timing. But also, speed of reconstruction is almost the same.

CT Fluoroscopy: More control where you need it

The new Aquilion ONE / PRISM Edition features a hybrid CT Fluoroscopy interface that enables one-handed operation with ergonomically designed controls and a versatile touchscreen tablet.

Has the new CT Fluoroscopy panel made things easier for your team?

Dr. Jenniskens: I really like working with the new console because it's faster and I have everything under my control. I have an overview of the whole process and can stay more focused for a more efficient and safe procedure.

Prof. Prokop: It's important to control the system from inside to do what you want quickly and without too much hassle. The Aquilion ONE PRISM Edition offers more flexibility for improved control and safety during the procedure.

Are there any benefits for the patient?

Dr. Jenniskens: Absolutely. You can now stay close to the patient, which is more comfortable for everyone. Patients also benefit from the new console because the procedure can be expedited, and they can be out of the room earlier.

www.canon-europe.com/medical



"I really like working with the new console. It's faster and I can stay more focused for a more efficient and safe procedure."

Dr. Sjoerd Jenniskens

Radiologist, Radboud University Medical Center, the Netherlands.

Accessories / Complementary Systems

Febromed · Get Up

The swivelling handle system for radiology offers the following highlights:

For the patients

- Independent moving
- Safe support in any position
- Safe motion for seniors and disabled people

For the staff

- Ergonomic working
- Reduced physical workload
- Fast changing of the sling



For the facility


- Various combination with existing systems
- Small space requirement
- Mounting on wall, floor or ceiling on customer request

DUNLEE IS YOUR PARTNER NOW AND IN THE FUTURE

BEYOND SOLUTIONS

Dunlee announces a new product develop ment: a CT replacement tube with liquid metal bearing. The new tube will join the popular DA200Ultra and DA200P40 tubes in the DA200 Series product portfolio. Now and in the future, you can rely on Dunlee's CT replacement tubes as an excellent, cost-effective solution and as an alternative to OEM tubes.

GCTechnology GmbH · CIRS Phantoms

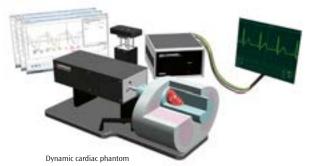

Highlights

- Electron density phantom for calibration
- Dynamic lung phantom
- CT dose phantoms
- Bone analysis CT simulator
- Plastic water and tissue equivalent materials
- Spiral / helical CT phantom
- AAPM CT performance phantom
- 3D sectional torso phantom
- Head phantom



DA200ULTRA



DA200P40



Cutted Model: Insert with Liquid Metal Bearing

Visit dunlee.com to learn more



I.A.E. · RTC 165



Accessories / Complementary Systems

Leoni · Cable Systems



Highlights

Leoni cable harnesses and systems for medical devices integrate a wide variety of functions: besides energy supply and control, they can also transmit light and data up to including so-called assist functions, such as the cooling of device components. Leoni provides wiring services for CTs and can collaborate in the development phase.

As a systems supplier, Leoni provides consultation and analysis at local level, development and design of the optimum cabling and the supply of prototypes. In addition, Leoni offers individual and system verification testing, as well as customer-specific logistics solutions.

PTW · Nomex Dosemeter



Highlights

- Diagnostic dosemeter,
 CE marked, class IIb certified, fully compliant with IEC 61674
- Suitable for CTDI measurements acc. to IEC 60601-2-44 using a 100 or 300 mm CT ion chamber
- Data and waveform export to Excel via USB or Bluetooth
- Accessories:
 CTDI head and / or body phantoms,
 CE marked, class I certified

Ultrasound Technologies · MedicCO₂LON



Highlights

Colonic Insufflator for CT Colonography. The MedicCO₂LON provides automated colonic distension with CO₂ gas for CT colonography procedures, providing reliable colon distension while improving patient comfort.

- State of the art design allowing ease of operation
- Near silent operation
- Large, colour touchscreen LCD
- LED backlight and wide view angle
- Compact, lightweight design
- Multilingual interface
- · Locking connectors

Varex Imaging · Cardinal CT Tube



Highlights

The Cardinal (Computed Tomography) CT tube is being designed into new OEM equipment and is also a direct replacement for the Stargate/CTR-2150 tube used in Philips Brilliance 6 and 16 CT scanners. The Cardinal has a high heat capacity with excellent image quality and throughput allowing for quicker imaging.

Varex Imaging · CT Replacement Tube - MCS 6074

Highlights

The MCS 6074 Computed Tomography (CT) tube is designed as a replacement tube for the GE Light-Speed and BrightSpeed families of CT scanners. The tube has a 200 mm (7.9") 140 kV, 4.7 MJ (6.3 MHU) maximum anode heat content, rotating anode insert. The insert features a 7° tungsten-rhenium facing on molybdenum with a graphite backed target. The MCS 6074 calibrates like the original, has long life bearings with 0.5 second full scans and offers a full 12-month warranty.



Varex Imaging · CT Replacement Tube – MCS 8064



Highlights

The MCS 8064 is an anode end grounded (AEG) Computed Tomography (CT) tube designed as a replacement tube for Lightspeed VCT series scanners. The tube has a 240 mm (9.4") 140 kV, 5.7 MJ (8.0 MHU) maximum anode heat content, rotating anode insert. The insert features a 7° tungsten-rhenium facing on molybdenum with a graphite backed target. The MCS 8064 offers lower life cycle costs and boasts a full 12-month warranty. The MCS-8064 installs and calibrates on the LightSpeed VCT like the original OEM tube.

Varex Imaging

Service Solutions Group






"The Varex Imaging Service Solutions Group provides world class imaging components and services. We deliver cost effective service components and support, enabling Varex customers to lower their cost to service X-ray Diagnostic imaging systems and optimize asset life cycles. Leveraging our global presence and service centers, Varex Imaging supports our customers servicing X-ray imaging systems around the world."

DAVID HURLOCK, Services Solutions Group

At Varex Imaging, we aren't just a supplier to our customers. We want to be an extension of their teams; a partner in their success; and to provide solutions to their challenges. We facilitate our customers' ability to provide world class service. Our knowledge, our people, our innovation and our 65+ year history of dedication to the X-ray imaging industry, make us who we are. At Varex Imaging, we have Solutions in Sight.

Varex Imaging Corporation's Service Solutions team account managers and support engineers focus on customers who service X-ray imaging equipment. Recognizing our customers' challenges, Varex Imaging provides cost effective components, technical support, and training to enable our customers to provide value in today's competitive healthcare service market.

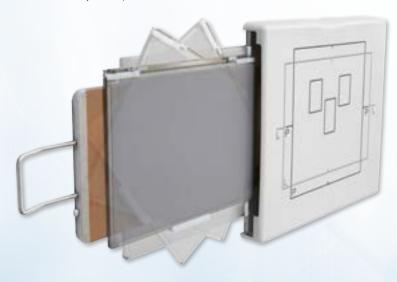



Varex Imaging Corporation recently acquired Reliable Healthcare Imaging Technology to broaden our range of X-ray tubes, and deepen our technical support resources. This acquisition adds expertise in preowned X-ray tubes that have been refurbished, system tested, and include a detailed system performance report.

The Service Solutions Group provides support to enhance our customers' ability to service systems with Varex components:

- A wide range of new X-ray tubes for Rad, R/F and CT applications
- Refurbished certified preowned X-ray tubes
- High voltage components for CT
- X-ray tube installation or installation support
- inventory in regional stocking locations
- Industry leading digital radiography upgrade packages and DR flat panel detectors

www.vareximaging.com

Magnetic Resonance Imaging



PET MR

Siemens Healthineers · Biograph mMR Channels Gradient 200 T/m/s¹ 45 mT/m¹ Up to 102 × 32



Highlights

- · Largest customer base of installed PET-MR systems worldwide
- State-of-the-art 3T MRI with 2nd order shim
- Comprehensive set of surface coils available for full range of MR-only exams
- Not only simultaneous, but synergistic PET-MR: MR-based motion compensation of PET images
- Whole-body MR-based PET attenuation correction including major bones
- Up to 10 bed positions with PET-MR
- · Latest applications available with syngo MR E11 software

7 Tesla

Gradient 80 mT/m ¹	Slewrate 200 T/m/s ¹	Channels Up to 64×64
[†] Maximum gradient amplitude and	slewrate can be applied simultaneously	
Highlights World's first 7T MR for clinical use Dual Mode – secu research and clinic 50 percent lighter nology** for easie clinical environme	re switch between cal operation* 7T magnet tech- r integration into	

- Double SNR for more precision***
- XR 80/200 gradients; Up to 16 channel parallel transmit in research mode
- Submillimeter BOLD fMRI precision for pre-surgical evaluation
- Latest applications available with syngo MR E12 software
- Additional metabolic information with ²³Na imaging and ³¹P spectroscopy
- *Research mode as part of dual mode is available as an option and not intended for clinical use
- **Compared to previous 7T generation ***Compared to 3T systems

3 Tesla

Canon · Vantage Galan 3T Gradient Slewrate 200 mT/m/s Channels 33/45 mT/m Highlights • Patient friendly 71cm wide bore and • EasyTech solutions silent scanning with Pianissimo Zen for automated scan-planning and • Fully integrated Deep Learning increased productivity

- Reconstruction: Advanced intelligent • Next generation scan techniques
 - · Enhanced throughput with Compressed SPEEDER, dockable table and live planning tool ForeSee View

Siemens Healthineers · Magnetom Prisma

80 mT/m¹

Channels Up to 204 × 128



Highlights

- A unique MR design driving innovation in research applications
- Unique scannner technology in one package: benchmark 3T magnet homogeneity; highest gradient amplitude

and performance with XR 80 / 200 gradients; advanced parallel transmit technology for zoomed imaging and ultra-high coil element density with Tim 4G designed for maximum SNR and extreme iPAT performance

- Driving the largest and most active MRI research network
- Latest applications available with syngo MR E11 software

Siemens Healthineers · Magnetom Vida with BioMatrix

200 T/m/s¹ Up to 60 mT/m¹ Up to 228 × 128 Maximum gradient amplitude and slewrate can be applied simultaneously

Highlights

•The first MRI scanner with BioMatrix technology

Clear-IQ Engine (AiCE). AiCE produces

exceptionally detailed MR images

with high SNR

- An all-new, 3T magnet with a large Field-of-View of $55 \times 55 \times 50 \text{ cm}^3$
- Up to 60 / 200 XT gradients for up to 25 percent higher SNR for DWI
- Up to 50 percent faster routine exams with Turbo Suite*
- Perform free-breathing liver dynamics and cardiac MRI with Compressed Sensing GRASP-VIBE and Compressed Sensing Cardiac Cine
- Reliable and predictable whole-body MRI exams in just 25 minutes
- Explore new diagnostic frontiers based quantitative information with MR Fingerprinting

*Data on file

Siemens Healthineers · Magnetom Lumina with BioMatrix

Gradient Channels 36 mT/m¹ 200 T/m/s1 180×32

Maximum gradient amplitude and slewrate can be applied simultaneousl

Highlights

- New 3T magnet with 70 cm open bore and large $55 \times 55 \times 50 \text{ cm}^3$ FoV
- Unique BioMatrix technology automatically adjusts to patient biovariability
- · New Turbo Suite acceleration packages enable up to 50 percent* faster clinical routine examinations
- GO technologies powered by artificial intelligence boost patient throughput
- Innovision** the revolutionary patient infotainment solution
- syngo Virtual Cockpit provides game-changing remote scanning support for standardized results across your MR system fleets
- ** Innovision is still under development and not yet commercially available. Its future availability cannot be guaranteed.

3 Tesla



· Low operating cost through low power consumption and Zero Helium boil off

• Fast break-even due to unmatched financial performance

• Latest applications available with syngo MR E11 software


1.5 Tesla









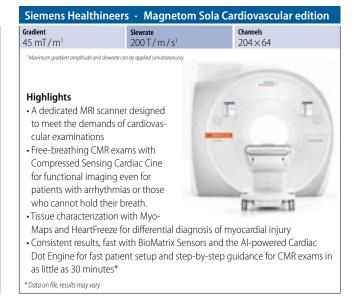


1.5 Tesla



 * Data on file, results may vary. ** Cardiac Triggering is still under development and not commercially available yet. Its


Compressed Sensing acceleration

future availability cannot be ensured.

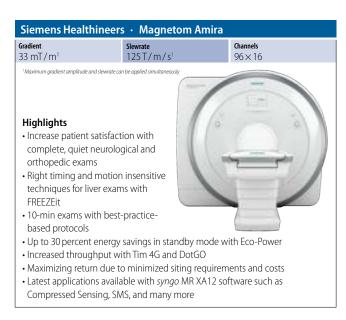




1.5 Tesla

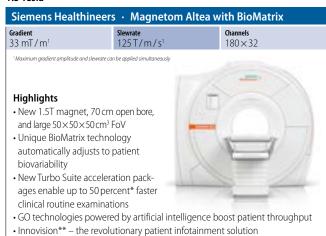


• syngo Virtual Cockpit provides game-changing remote scanning support for

** Innovision is still under development and not yet commercially available. Itsfuture availability cannot be guaranteed


standardized results across your MR system fleets

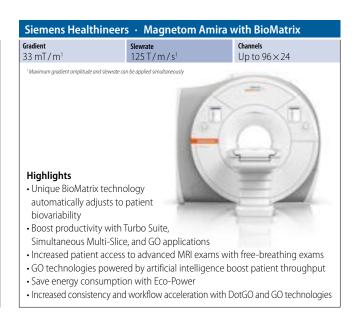

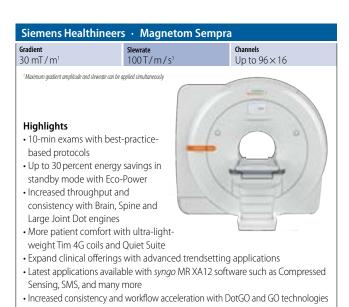




Oncology

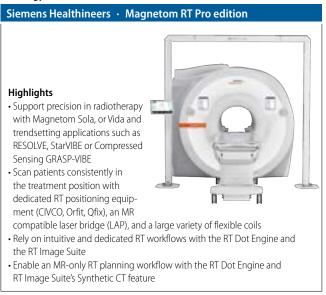

Open

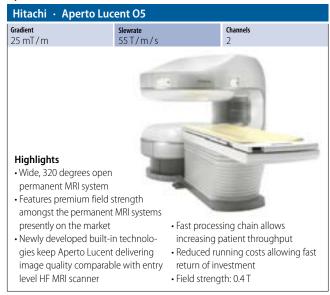

Open




• Field strength: 0.3 T

· High magnetic field homogeneity





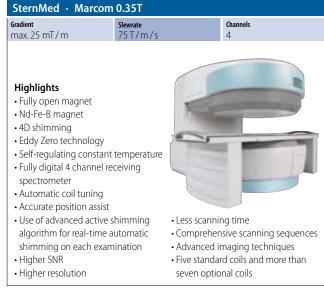



MRI Coils




MRI Coils



- nation with a radiotherapy mask · MRI control before and after stereo-
- tactic radiotherapy · High-resolution MR diagnostics of head and neck
- Transversal, sagittal, coronal and tilted images possible
- Removable double mirror for claustrophobic patients
- *in development




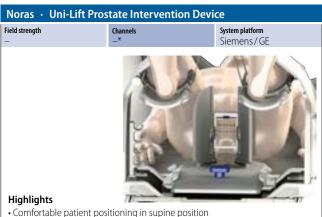

- Optimal utilization of the patient space in 70 cm Bore MR Systems
- Transperineal access to the prostate
- · MR-guided biopsy and therapy
- · Focal biopsy and therapy instead of "blind" biopsy
- *Compatible with standard MR coil portfolio



• Brainlab AIR Support

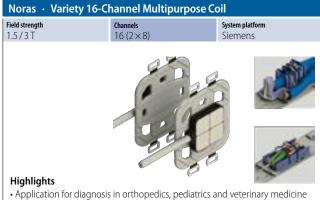

- High signal quality based on a design with 8+8 array elements with high coil element density
- High resolution examinations of even small body regions with reduced scan
- Slim design and optional dedicated positioning aids enable coil placement close to anatomy of interest

Accessories / Complementary Systems



- Enabling fetal cardiac MRI, smart-sync enables a secondary imaging modality for prenatal diagnosis
- Improving adult cardiac MRI, smart-sync is easy to apply and independent of field strength

Accessories / Complementary Systems

allMRI GmbH · MRI safe metal free cleaning tool set



Highlights

- Velcro mop holder with microfiber mop 40 cm
- Glass fiber handle 145 cm long with mop press wringer, optional with complete cleaning cart available

allMRI GmbH · Mobile MRI procedure lamp



allMRI GmbH · Foldable MRI wheelchair



- Including the ball bearing
- •Two swing out adjustable footrests and armrests
- Solid rubber tires

phantom

Febromed · Get Up

The swivelling handle system for radiology offers the following highlights:

For the patients

- Independent moving
- Safe support in any position
- Safe motion for seniors and disabled people

For the staff

- Ergonomic working
- Reduced physical workload
- Fast changing of the sling



For the facility

- Various combination with existing systems
- · Small space requirement
- Mounting on wall, floor or ceiling on customer request

GCTechnology · CIRS Phantoms



Leoni · Cable Systems




Highlights

Leoni cable harnesses and systems for medical devices integrate a wide variety of functions: besides energy supply and control, they can also transmit light and data up to including so-called assist functions, such as the cooling of device components. Leoni provides wiring services for MRTs and can collaborate in the development phase.

As a systems supplier, Leoni provides consultation and analysis at local level, development and design of the optimum cabling and the supply of prototypes. In addition, Leoni offers individual and system verification testing, as well as customer-specific logistics solutions.

RADBOOK 2020 31

MR head distortion Phantom for SRS

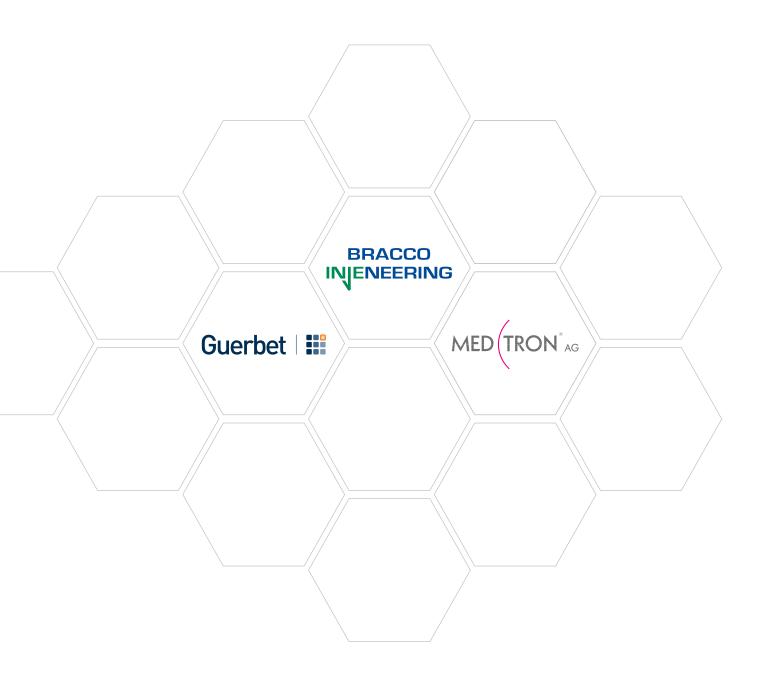

Bracco · CT Exprès FlowRate Application Pressure 0.5 - 9.9 ml/s *

Highlights

- Syringeless injector
- Direct injection from contrast media bottles
- Air and occlusion detection on fluid channels
- Unidirectional flow of fluid
- · Locking and automatic filling
- Digital interface, dual touch screens
- DiluJect (optional): contrast media and saline are injected in rapidly alternating flow through the injector
- Day Set III HP designed for 24 hours
- Pre-warmed contrast media maintained at 37°C
- * in steps of of 0.1 ml/s



Bracco · EmpowerCTA+ FlowRate 40 – 325 psi $0.1 - 10.0 \text{ ml/s}^*$ Highlights • Syringe: 200 ml (CM), 200 ml (NaCl) • Tilt sensor / lockout · Arming at the injector • Independently rot0ating and very compact injector head (270 degrees) • Integrated electroluminescent display

- Modular flexibility of components and Windows based software allow optimal serviceability and enhanced expandability
- Touch-screen color LCD display and intuitive software
- Multi-patient tubing kit for CT: Must kit PLUS CT
- * 0.1-10.0ml/s in user-specified increments of 0.1 ml/s

Guerbet · OptiVantage single use



Bracco · EmpowerMR

Application MR	Pressure	FlowRate
MR	40 – 300 psi	0.1 – 10.0 ml/s *

Highlights

- Syringe: 100 ml (CM), 100 ml (NaCl)
- Hydraulic injector system
- MRI compatible through the use of polymers and non-ferromagnetic metals
- Little contrast media waste due to the very short distance between injector head and patient
- Very lightweight injector head
- No active components in the shielded room (no battery)
- · Unique multi-patient tubing kit for MR: Must kit PLUS MR
- * 0.1-10.0ml/s in user-specified increments of 0.1 ml/s



Dual head CT contrast deliverry system

- · Scan delay, phase delay, auto-fill, auto purge
- Timing bolus, inject delay, patency check
- Fully programmable touchscreen powerhead
- Scanner interface to CAN Open Class 4*
- · OptiBolus bolus shaping software extends the window of imaging opportunity
- · Configurations:

Highlights

- Pedestal and ceiling mount options
- Loading, filling & priming: Automatic/manual • Simultaneous injection: 10% - 90% (5% steps)
- Heater: 37° ± 3°
- Connectivity with Contrast&Care (optional)

Components and consumables certified by the manufacturer

FlowRate 0.1-10 ml/s

* dependent on scanner manufacturer

Guerbet · Illumena Néo

Application 5.2 – 82.7 bar¹ / 5.2 – 21 bar² 0.1 – 40 ml/s¹/0.1 – 10 ml/s² CT/Angio/Cardio

Highlights

Multi-mode contrast delivery system

- · High visibility screen
- One finger operation fill bar
- Single or multi-injection procedures
- · Switch between operating modes
- Hand switch and foot switches available
- Air Detection Aid & Warning System (ADAWS) identifies blood, empty syringes and air bolus
- · Configurations: Pedestal, ceiling or table mount
- Heater: 37° ± 3°
- · Connectivity with Contrast&Care (optional)

Components and consumables certified by the manufacturer

1 Anaio mode/2 CT mode

Guerbet · FlowSens

FlowRate 0.3 - 10 ml/s 21 bar

Highlights

Syringeless CT contrast delivery system

- · Advanced touchscreen interface
- Only few seconds between patients
- 12H manyFlow
- (closed pre-connected day set)
- Secufill patient line
- (scientific study on demand) • All available media containers
- Check-valves (no backflow)
- 4 Air sensors
- Temperature maintenance
- · Pressure monitoring: graphical and numerical


 Simultaneous injection (optional) 20 - 80 %

- Vein test · Loading, filling &
- priming: Automatic · Configurations:
- pedestal & ceiling · Connectivity with Contrast&Care
- (optional)

Components and consumables certified by the manufacturer






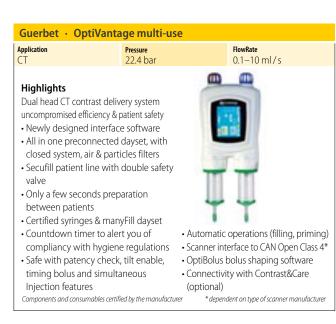





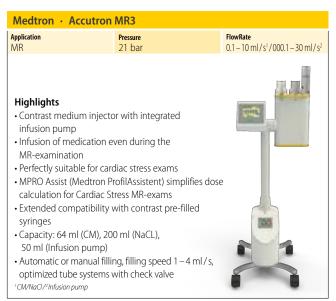



Guerbet Diagnostic Imaging has designed a portfolio of interconnected contrast imaging solutions to enhance your decision-making at each point of the patient journey from diagnosis, to treatment, to follow-up, so you can focus on what matters most, efficiently improving patient outcomes. This is UNIK.



Medtron · Accutron CT-D		
Application	Pressure 21 bar	FlowRate 0.1 – 10 ml/s
At-a-glance visualisation status Aluminium housing – wa CANopen, Class 4 or Phil Capacity: 200 ml (CM), 20 Syringe: Automatic or ma	e holder for Easy Loading note control (as an option) contrast extravasation gement and liquid proof de of injection parameters an all or ceiling suspension ips interfaces (as an option	esign d



Accessories / Complementary Systems

Bracco · Nexo [Dose]



Highlights

- Nexo [DOSE] supports compliance with imminent European Directive (2013/59/EURATOM)
- Single-server, fully automated system enables enterprise-wide data acquisition
- Multi-modality, vendor-neutral software minimizes installation time and costs
- Customized e-mail alerts help improve control and implement the ALARA* principle
- Global Dose Registry for study dose comparison
- * ALARA (as low as reasonably achievable)

Guerbet · Contrast&Care



Highlights

Contrast&Care is a solution dedicated to contrast dose management. It connects to all Guerbet injectors and Hospital Information Systems (RIS, PACS, EMR...) and collects all relevant data about contrast media usage, patient history, and injector activity. Contrast&Care facilitates the traceability of contrast media and provides several tools that help imaging centers optimize contrast media consumption.

Guerbet · Dose&Care




Highlights

Dose&Care is a state of the art radiation dose monitoring solution, which allows documenting patient exams, understanding the reasons for excessive exposure and monitoring activities related to patient exposure. It provides the means to remain compliant with an ever-evolving regulation while improving the workflow and ensuring patient safety.

Interventional Systems



Multi Modality Suites

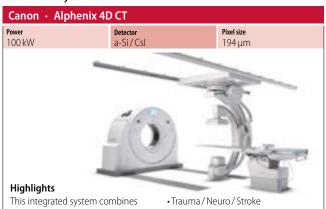

premium CT and ceiling-mounted angiography technology. The perfect diagnostic and treatment set-up for

high-risk procedures in various inter-

ventional segments such as: · Interventional Oncology

- General Vascular
- · Additional or Backup CT
- Detector size: $30 \times 30 \,\text{cm} / 30 \times 40 \,\text{cm}$ The system is available with two different CT configurations: Aquilion One Genesis and Aquilion Prime

Siemens Healthineers · nexaris Angio-CT Pixel size Detector 100 kW 160 µm

Highlights

- First hybrid suite with a common coordinate system that fuses images instantly
- Direct access to angio and CT with Quick Switching
- Efficient multi-room configurations to share imaging equipment
- Enabling combined CT and angio guidance in one session

Siemens Healthineers · nexaris Angio-MR-CT


Pixel size 100 kW 160 µm

Highlights

and digital acquisition system.


- Unprecedented imaging capabilities with all the advantages of angiography, MR, and CT imaging
- Patient transfer without repositioning for barrier-free intraoperative imaging

Bi-Plane

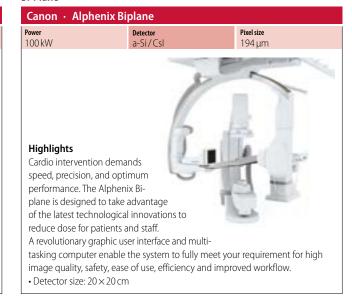

• Innovative treatments that combine multiple imaging modalities in a single procedure



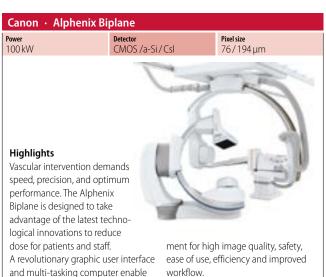

• Detector size: 30 × 30 cm / 30 × 40 cm

38 RADBOOK 2020

• Detector size: 20 × 20 cm / 30 × 40 cm

the system to fully meet your require-

Bi-Plane



- SCORE Navi / Navi+Plus
- · SMART design concept
- Comprehensive dose management package

Siemens Hea	Althineers · Artis zee, Ar	tis Q, Artis Q.zen
Power 100 kW	Detector a-Si/Csl	Pixel size 184 µm
Highlights		Ψ
Biplane system	for interventional	-
offers high per	rtis biplane system formance in inter-	
ventional imag	ing combined with	// 50
J ,	ane imaging position	

- · Single plane operation with extended position flexibility enabled by rotated table
- · Ergonomic system controls for smooth table-side operation
- 3D acquisition rate up to 75 f/s

Single Plane

ing technology

• SCORE Pro Advance image process-

• Unique pioneering imaging technology:

motion-tolerant SCORE RSM

Canon · Alphenix Hy	brid+	
	Detector a-Si/Csl	Pixel size 194 µm
	TAP	
Highlights	1	
The combination of the Alphenix ceiling mounted	1/23	
C-Arm with a fully integrate dedicated surgical table, e.		
acarcatea sargical table, c.	9.	

Maguet Magnus or Trumpf Trusystem 7500, perfectly meets the requirements for the rapidly growing number of hybrid procedures.

With its unique double sliding C-Arm the Alphenix Hybrid+ provides ultrafast whole body 3D coverage, free head access and a unique lateral C-Arm stroke for better ergonomics, improved productivity and stunning 3D images from head to toe.


• Detector size: 30 × 40 cm

Canon · Alphenix Hybrid Power 100 kW Pixel size 194 µm

Highlights

The combination of the Alphenix ceiling mounted C-Arm with a fully integrated dedicated surgical table, e.g. Maquet Magnus or Trumpf Trusystem 7500, perfectly meets the requirements for the rapidly growing number of hybrid procedures. Its flexibly designed ceiling rail system allows perfect patient access in any situation.

• Detector size: $30 \times 40 \,\text{cm} / 30 \times 30 \,\text{cm} / 20 \times 20 \,\text{cm}$

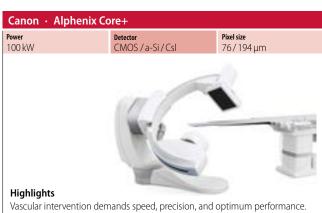

The flexible Alphenix Core+ is designed to take advantage of the latest technological innovations to reduce dose for patients and staff.

A revolutionary graphic user interface and multi-tasking computer enable the system to fully meet your requirement for high image quality, safety, ease of use, high efficiency and improved workflow.

• Detector size: 30 × 30 cm / 30 × 40 cm


Canon · Alphenix Core+ Power 100 kW Detector Pixel size 194 µm Highlights


Cardio intervention demands speed, precision, and optimum performance. The flexible Alphenix Core+ is designed to take advantage of the latest technological innovations to reduce dose for patients and staff.

A revolutionary graphic user interface and multi-tasking computer enable the system to fully meet your requirement for high image quality, safety, ease of use, advanced efficiency and improved workflow.


• Detector size: 20 × 20 cm

Single Plane

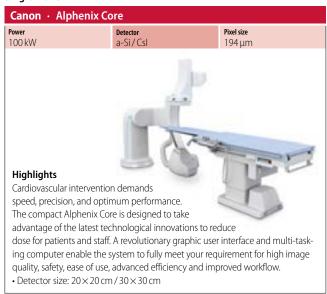

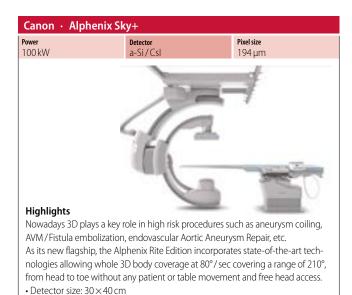




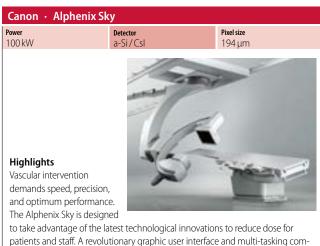

to take advantage of the latest technological innovations to reduce dose for patients and staff. A revolutionary graphic user interface and multi-tasking computer enable the system to fully meet your requirement for high image quality, safety, ease of use, high efficiency and improved workflow.

• Detector size: $30 \times 30 \, \text{cm} / 30 \times 40 \, \text{cm}$

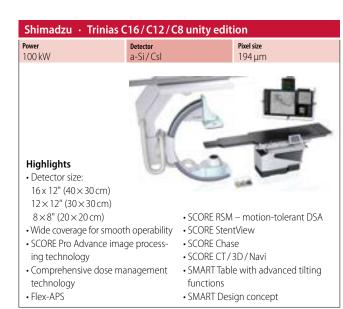



Single Plane



- Detector size:
- 12×12 " (30×30 cm)/ 8×8 " (20×20 cm) • Wide coverage for smooth operability
- · SCORE Pro Advance image processing technology
- Comprehensive dose management technology
- SCORE RSM motion-tolerant DSA
- SCORE StentView
- SCORF Chase
- SCORE CT/3D/Navi
- SMART Table with advanced tilting functions
- · SMART Design concept

Siemens Healthineers · Artis pheno Pixel size Detector 100 kW 160 μm Highlights Artis pheno delivers the right images for preprocedural planning, intraoperative guidance, and immediate assessment Detector: zen40HDR, hi-res cristalline silicon / CsI, $30 \times 40 (2,496 \times 1,856 \text{ px})$, 160 µm


- Optimal treatment of any patient virtually regardless of patient size, condition, or positioning needs
- · Lower complication rates and improved outcomes even in most complex procedures
- Designed to support infection control measures to maintain highest standards for infection control

Siemens Healthineers · Artis zee, Artis Q, Artis Q.zen Power 100 kW Detector Pixel size 154/160/180 µm Highlights The Artis ceiling-mounted system enables clinicians to care with greater ease, precision and flexibility. Detector: $20 \times 20 (1,024 \times 1,024 \text{ px}), 184 \mu\text{m}$ $30 \times 40 (1,920 \times 2,480 \text{ px}), 154 \mu\text{m}$ $30 \times 40 (1,920 \times 2,480 \text{ px}), 154 \mu\text{m}$ • Positioning flexibility that supports any angle • Ergonomic system controls for smooth table-side operation • 3D acquisition rate up to 75 f/s • Complete 3D-portfolio including cross-sectional imaging with syngo DynaCT and syngo 3D Roadmap

Siemens Healthineers · Artis zee, Artis Q, Artis Q.zen Pixel size 100 kW 184 µm Highlights The Artis floor-mounted system enables clinicians to care with greater ease, precision and flexibility for small rooms. • Small footprint of 29 gm • Slim-line design for easy patient access • Ergonomic system controls for smooth table-side operation • 3D acquisition rate up to 75 f/s · Complete 3D-portfolio including cross-sectional imaging with syngo DynaCT and syngo 3D Roadmap

Siemens Healthineers · Artis zee Power 100 kW Pixel size 154 µm Highlights Artis zee multi-purpose is designed to meet the demands of interventional radiology and fluoroscopy. The system left suspension meets the needs of endoscopic applications in gastroenterology • Detector: 30×40 (1,920 \times 2,480 pixels), 154 μm • Ergonomic system controls for smooth table-side operation

• 2k imaging with highly practical and user-friendly handling features

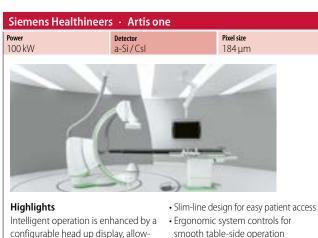

a completely new, intuitive way. • Detector: as30, (1,560 × 1,420 pixels),


ing you to interact with the system in

Small footprint of 25 gm²

- smooth table-side operation
- Full patient coverage imaging up
- Integrated 3D-Imaging and review with acquisition rate up to 66 f/s

• 3D acquisition rate up to 75 f/s

Surgical Flat Panel C-Arms

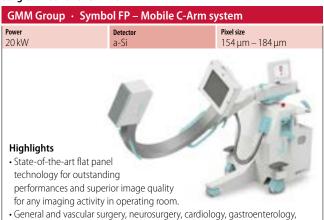


• Easy patient positioning thanks to the wide C-Arm opening.

management of the operating parameters.

• Detector Size: 26 × 30 cm


• Exclusive user interface with LCD touch screen display ensuring complete

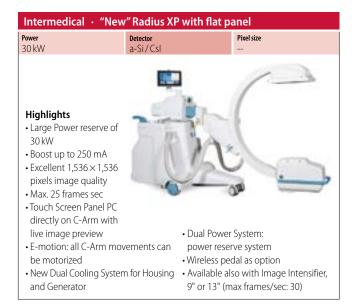



positioning

- Designed for surgery
- 13s true 360° 3D scan Fully mobile
- Flexible intra-operative 2D- and 3D-imaging
- · Large 2D-image size and large 3D scan volume up to 40 cm width
- Seamless integration in OR workflow
- simple control panel
- Position memory remembers four scan positions
- Easy draping of the breakable gantry
- · Seamless integrating with Stealth-Station Navigation
- Detector size: 30 × 40 cm

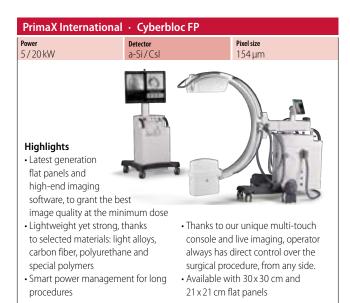



Siemens Healthineers · Cios Spin Pixel size **Power** 12/25 kW $30 \times 30 \, \text{cm}$ 152 µm Highlights • More certainty in demanding cases with precise intraoperative quality control based on 3D technology • More efficiency in intraoperative 3D with Easy 3D package

• More cost-effectiveness in surgery through intraoperative corrections based

42 RADBOOK 2020

on 3D images

Surgical Flat Panel C-Arms



- Up to 25 percent more coverage* even during image rotation with Full View FD
- See and do more with a powerful 25 kW mobile C-Arm
- Effortless operability full table-side control and single-touch positioning
- * Compared to today's conventional 33 cm image intensifiers





- · Orthopaedic, head, spine, thorax, abdomen, vascular, cardiac
- · Large C-Arm depth and wide orbital rotation
- Adjustable height & angle of medical displays
- Dynamic FPD with high DQE and MTF
- Removable grid
- Advanced functions : APR, post-processings, DSA
- DICOM connectivity
- Detector size: 21 × 21 cm/30 × 30 cm

Stephanix · Omniscop DReam S Power 4.2/5 kW Pixel size 21×21 cm/30×30 cm 194/205 μm Highlights Orthopaedic/Urology/Cerebral/ Thoracic/Pain therapy/Peripheral vascular using DSA function • Single unit system, all components integrated into the C-arm stand Very small footprint • 4 Mpixel 27" medical monitor on an · Advanced functions:

- articulated arm, adjustable height and angle
- Dynamic FPD with high DQE and MTF
- APR, post-processings, DSA
- DICOM connectivity
- Detector size: 21 × 21 cm/30 × 30 cm




- 30×30 or 21×21 cm
- Digital memories: 1.5k×1.5k, with acquisition up to 25 fps $1k \times 1k$
- 30 kW H.V. generator
- Dual cooling systems for immediate and effective heat removal
- Dual power system: Power reserve system
- E-Motion system: C-arm movements, fully motorized (optional)
- 215 mm horizontal run
- (17 mm in motorized version) • Wide orbital movement 150°
- C-arm lateral rotation: ± 180°
- 10" touch screen control console



for pediatric application Full DICOM connectivity

Surgical Flat Panel C-Arms

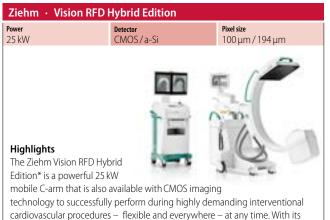

cardiovascular procedures – flexible and everywhere – at any time. With its zero room preparation, the comprehensive mobile hybrid solution easily takes your OR to the next level. Plug in your system and start your hybrid procedure.

- Detector size: 31 × 31 cm / 20.5 × 20.5 cm (CMOS) · 30 × 30 cm (a-Si)
- * Ziehm Vision RFD Hybrid Edition represents a group of optional hardware and software that creates an option package on the device named Ziehm Vision RFD.

Ziehm · Vision RFD	BD	
Power 25 kW	Detector CMOS/a-Si	Pixel size 100 µm / 194 µm
C Highlights		
_	, ,	aoperative control, the Ziehm e CT scans and costly corrective
	th Ziehm iterative rec	onstruction to minimize fan and

• Detector size: 31 × 31 cm (CMOS) · 30 x 30 cm (a-Si)



and demanding procedures and the intuitive Ziehm usability concept* helps surgeons ensure consistently high clinical standards. This impressive feature lineup make the systems ideal for challenging interventions.

- Detector size: 31 × 31 cm / 20.5 × 20.5 cm (CMOS) · 30 × 30 cm (a-Si)
- *The Usability Concept includes a variety of hard- and software features. Due to regulatory reasons the availability of each feature may vary. Please contact your local Ziehm Imaging sales representative for detailed information.

Power 2.4 kW	Detector CMOS / a-Si	Pixel size 100 μm / 150 μm
Highlights Now in the upg line*, the Ziehm features an enha	Vision FD	
thanks to the Ad In addition, fine patient outcom now also availal allows to cover • Detector size: 2	dvanced Active Cooling – is de y tuned workflows and new s es and further increase produ ole with a 31 × 31 cm/a-Si flat	oftware features help to optimize ctivity. The Ziehm Vision FD is -panel. The bigger detector size rthopedic and vascular surgery. 1 cm (a-Si)

Surgical II-C-Arms



Del Medical · **Vision R*** by Ziehm CCD-matrix Power 20 kW 23 cm/31 cm Highlights · Feature rich mobile surgical fluoroscopy system with high frequency pulsing generator • 9 inch or optional 12-inch image intensifier · Advanced Active Cooling - the gen-

- Powerful 20 kW monoblock generator with rotating anode
- Up to 30 pulses per second for minimizing dose exposure while maximizing image quality
- * Available for US market only
- erator remains at the ideal operating temperature for uninterrupted use
- Dual 19-inch flat screen stand-alone monitor cart and adjustable C-arm mounted touchscreen display

and powerful performance

tion image and lower dose

pulsed fluoroscopy modes

• Dose reduction through various

· High performance with high resolu-





• Wider SID, wiser operation


• Two control panel on both side

• C-arm movement fully counter

• 19" dual monitor (43" monitor is option)

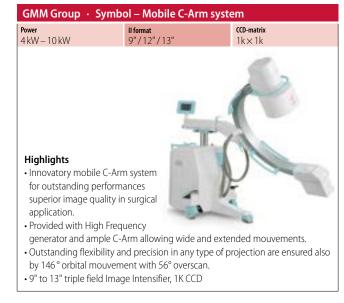

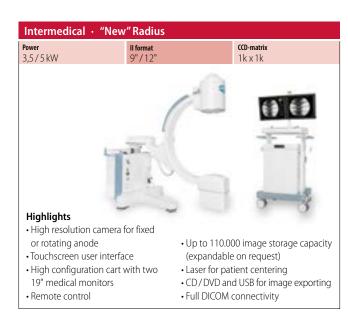

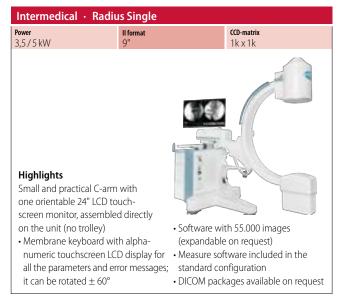





availability



• Optimally balance image quality and dose - with IDEAL

• Increase asset utilization – with preventive maintenance and high system




• Seamless integration with 3D Navigation – through fully digital NaviLink

(Option) with major vendors

- Stephanix · Omniscop Series CCD-matrix Up to 15 kW Highlights Surgery, traumatology, orthopedics, vascular... · Wide range of movements, large orbital rotation, small footprint • Touch screen user interface • High resolution CCD camera coupled • Post-processing software highlight with Thales Image Intensifier tiny details • Collimator with motorised and rotat-· Advanced functions: APR, DSA, ing iris, continuously adjustable DICOM connectivity
- SternMed · Xenox C200 Power 5 kW CCD-matrix Highlights • DSA 25 frames/ses • 2 × 19" active matrix LCD, anti-glare and hard coated (standard) • 2 × 19" TFT with 2000:1 contrast ratio (optional) • 1 TB hard disks up to 440,000 images • 270° on each side arm rotation storage capacity Dose meter (DAP meter) • Touch screen keyboard with 5.7" LCD • Full Dicom (optional) • 200 mm horizontal C-arm run display with the possibility to rotate · Orbital movement 125° of $\pm 60^{\circ}$ • 1k CCD camera delivers sharp and • 12° on each side C-arm swiveling measure & cine software detail-rich images



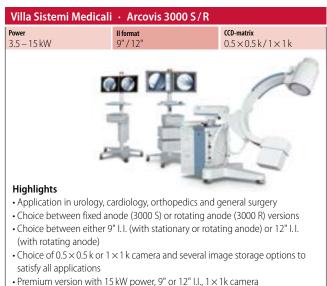



dose levels. The outstanding power reserves make this C-arm particularly

cardiology, including AAA, PTCA and coronary angioplasty.

suited for demanding procedures in cardiovascular surgery and interventional







Accessories/Complementary Systems

Canon Electron Tubes & Devices · X-ray Image Intensifier



Highlights

- Suitable for mobile C-Arms
- Smart design with smooth surfaces
- Excellent performance and high reliability
- Advanced simulation technologies used in development and production
- Our unique technologies provide a high Gx value, reducing radiation exposure to the patient.
- · Environmentally friendly
- · Compliant with the RoHS directive
- Free from hazardous substances such as hexavalent chromium and cadmium
- Detector: Xray Image Intensifier
- Size: Field size 9 inch, 9/6/4.5 inch Output image size Ø 20mm, Ø 25mm
- Design: For C-Arm

Canon Electron Tubes & Devices · LM-Angio Tube Capacity 2.1 MHU(Anode Heat Capacity) 100 kW



- Highlights • For angiography systems
- (2.1 MHU) · Uses a liquid metal bearing
- · Our unique liquid metal bearing technology
- Compact Housing provides a long tube life, quiet operation, high stability, and excellent reliability.

Canon Electron Tubes & Devices · Angio Tube assembly



- (3 MHU)
- Uses a liquid metal bearing
- · Our unique liquid metal bearing technology
- Compact Housing provides a long tube life, quiet operation, high stability, and excellent reliability.

GCTechnology · CIRS Phantoms



- Multi modality abdominal biopsy phantom (for CT, US, MRI)
- Biopsy breast phantom
- Thyroid training phantom
- Prostate training phantoms family
- Kidney training phantom
- · Vascular access training phantom

I.A.E. · C30-RTM 70



Highlights


- Rotating anode X-Ray tube unit designed for mobile c-arm equipment
- Lead lined single piece aluminium body, internal pump for oil circulation, to improve thermal exchange
- Choice of HT cable socket: Parker or Claymount mini
- Optional remote water-air heat exchanger increases heat dissipation to 500W continuous for demanding interventional applications
- · Water cooling can be mounted or upgraded on field

Artificial Intelligence



Artificial Intelligence

AlgoMedica · PixelShine CT

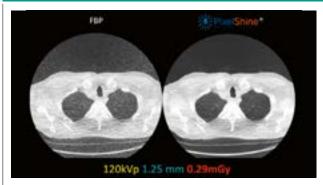

Highlights

PixelShine is an Artificial Neural Networking (ANN) which processes ultralow-dose images improving quality to emulate high dose images. It is vendor neutral, providing crystal clear, diagnostic quality images across the entire radiology. Improve diagnostic confidence, increase patient safety, substantially lower radiation dose, Integrate without changing radiology workflow.

Canon · Advanced intelligent Clear-IQ Engine for MR



Highlights

The power of AI is brought to routine MR imaging by Canon Medical's deep learning reconstruction technology: Advanced intelligent Clear-IQ Engine (AiCE). AiCE is the world's first fully integrated DLR technology for MRI and produces exceptionally detailed MR images. AiCE intelligently removes the noise from the images, which results in higher SNR enabling increased resolution or decreased scan time.


Fujifilm · REiLI



Highlights

Under the REiLI brand, Fujifilm is developing AI technologies that strongly support diagnostic imaging workflow, leveraging the combination of deep learning and Fujifilm's image processing heritage. Working with both Fujifilm developed algorithms and market leading specialist vendors, the REiLI platform can automate alerts and send critical information directly to the relevant clinician, increasing both the speed and accuracy of diagnosis augmenting the decision making process.

Fujifilm · Lunit nano

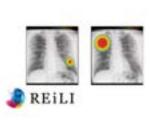

Highlights

- · Ground breaking compact and lightweight cart. By utilizing the ultrahigh sensitive D-EVO II and virtual grid technology, resulting in a compact mobile X-ray cart only weighing 90 kg
- · Integrated AI CAD software Lunit insight CXR-MCA provides abnormality score by heatmap
- · Major chest abnormalities including nodule, consolidation, pneumothorax, atelectasis, fibrosis, pleural effusion, pneumoperitoneum and mediastinal widening supported
- · Providing an advanced workflow and improved patient care pathway

mediaire · mdbrain



Please visit us at

www.healthcare-in-europe.com



Highlights

mdbrain is an easy-to-use, platform independent & economic software solution for brain volumetry and the characterization of white matter lesions. mdbrain uses a combination of several deep learning methods to derive meaningful, well-visualized quantitative reports from MR images. The seamless integration into customers' local PACS allows for fast processing and the best possible data protection.

Artificial Intelligence

Mindray · ME



Highlights

ME is the first laptop ultrasound system powered by ZST+ platform which will be newly launched at ECR 2020. Its Al empowered analysis tools, smart cardiopulmonary assessment solution, help to quickly assess patient status under critical care environment.

Find Mindray at ECR 2020, Expo X5/543 to explore more.

ScreenPoint Medical · Transpara



Highlights

ScreenPoint Medical is the leading developer of AI driven image analysis technology for automated reading of 2D and 3D mammograms. With proven accuracy matching experienced radiologists, Transpara is the most advanced commercially available multivendor AI solution (CE marked for 2D and 3D, FDA cleared for 2D and pending 3D). Find us at ECR 2020, stand #112, X1.

Siemens Healthineers · Al-Rad Companion



Highlights

The Al-Rad Companion, a family of multi-modality (MRI, CT, X-ray) cloud-based augmented workflow solutions may increase your diagnostic precision when interpreting medical images. Its solutions provide automatic post-processing of imaging datasets through our Al-powered algorithms. The automation of routine workflows with repetitive tasks and high case volumes helps you to ease your daily workflow – so that you can focus on more critical issues.

SIUI · Apogee 6500 Color Doppler Ultrasound Imaging System

Highlights

Apogee 6500 is the best proof of SIUI's perseverant endeavor in ultrasonic innovations throughout the decades. Special features are:

- OB/GYN Auto OB is a deep learning technology based on deep neural network. It supports intelligent analysis and precise measurement for obstetric images with high efficiency, applicable for the most common measurement items.
- Cardiology Auto EF automatically measures cardiac volumes based on speckle tracking technology. The system identifies the endocardium and establishes a ROI using an anatomical database.

LV EF, ESV, EDV and GLS are calculated at the same time. It provides an effortless and efficient way to acquire accurate 2D EF and volumes, which is a decision supporting tool to cardiac function diagnosis.

$\textbf{Therapixel} \; \cdot \; \textbf{MammoScreen}$



Highlights

MammoScreen: an Al-based concurrent reading aid for interpreting screening mammograms. It detects and characterizes suspicious regions of the breasts and generates a summary report containing the recommended action for the mammogram level and characterization of lesions and breasts suspiciousness.

- Trustable characterization score, clear recommended actions
- More confident decision making
- Workflow untouched
- Optimal human / machine complementarity

CE mark pending



Please visit us at

www.healthcare-in-europe.com

IT Systems



Agfa HealthCare · Enterprise Imaging VNA



Highlights

A robust solution for enterprise archiving of DICOM and non-DICOM data. As part of the Enterprise Imaging solution, the VNA consolidates all your imaging data, from multiple systems, departments, facilities and vendors, into a central clinical data foundation. Your data ownership, migration and storage costs are reduced, while management is simplified

Agfa HealthCare · Enterprise Imaging for radiology



Highlights

Agfa HealthCare Enterprise Imaging for Radiology is a unified imaging management platform that provides PACS, reporting, advanced image processing capabilities and integration of clinical information. The solution offers diagnostic tools and powerful task-based workflow, designed to achieve gains in clinical productivity.

Del Medical · EvoView PACS



Highlights

- Highly efficient and effective diagnostic image management and storage
- Virtually limitless user customization and accessibility
- Android, Apple, and Windows compatible for image viewing and control on any stationary or mobile workstation device
- Healthcare IT friendly with centralized maintenance control
- Powerful tools and features such as built in RIS, EMR and referring physician portal, permitting easy and secure image sharing
- Supports multiple languages

Examion · X-AQS




Highlights

Universal software platform for radiological image acquisition and management of all medical image data.

- High quality images in a few clicks
- Intuitive GUI with clear menu structure and icons
- Modular architecture, adaptable to all needs
- Certified diagnostic viewer with comprehensive measurement functions

Fujifilm · Synapse 3D

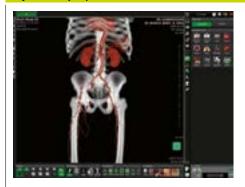


Highlights

Synapse 3D is Fujifilm's vendor-neutral advanced visualization platform with more than 50 clinical modules.

The advanced pre-surgical planning modules allow surgeons & clinicians to plan the most efficient, least invasive surgical activities. Ensuring that patient outcomes are maximised

Fujifilm · Synapse 5

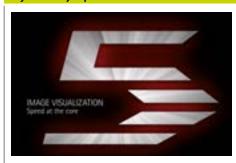

Highlights

Fully scalable, Synapse PACS can fit small installations, designed to handle few modalities, to the complete management of a regional diagnostic imaging solution. PACS provides one-touch access to all relevant patient information. Synapse makes the most of your infrastructure with Fujifilm's unique patented compression technology. This technology combines high compression ratios without sacrificing image quality. Giving clinicians easy access to all of the relevant tools to provide the most accurate and efficient diagnosis.

	Manag	gement		iving		Image Dis		
	RIS	Business Intelligence	Enterprise PACS	VNA	Inhouse	Teleradiology	Portal Solution	Cloud Computing
Agfa HealthCare	•	•	•	•	•	•	•	•
AlgoMedica								
BMS Informationstechnologie								
Bracco Injeneering								
Canon Medical Systems	•	•	•	•	•	•	•	
Del Medical			•		•	•	•	
Examion	•		•	•	•		•	
Fujifilm Medical Systems	•	•	•	•	•	•	•	•
Guerbet								
Image Information Systems	•		•	•	•	•	•	•
i-Solutions Health	•	•						
ITZ Medicom	•	•	•	•	•	•	•	•
mediaire								
medigration	•		•		•	•	•	•
Mindray								
Nexus / Chili	•		•	•	•	•	•	
Protec			•					
ScreenPoint Medical								
Siemens Healthineers	•	•	•	•	•	•	•	•
SIUI								•
Therapixel								
Visaris	•		•					

Workstations

Multimodality	Advanced Visualization	Mammography	Clinical Applications	CAD	Dose Management	Artificial Intelligence	
•	•	•	•	•	•	•	AGFA 🐠 HealthCare
•	•				•	•	& Algo Medica
					•		BMS INFORMATIONSTECHNOLOGIE® GMBH
					•		BRACCO INJENEERING
•	•	•	•	•	•	•	Canon
•			•				DEL MEDICAL
•		•	•				EXAMION' x-Ray Systems - Digital Imaging - Service
•	•	•	•	•	•	•	FUJ¦FILM
					•		Guerbet !!!
•			•				IMAGE
					•		N-SOLUTIONS HEALTH
•	•	•	•	•	•	•	itz-medi,com
						•	medi ai re
•	•	•	•		•		dender gruppe medigration
•			•	•		•	mindray
•	•	•			•		nexus chili inaging & radiology solutions
•							PROTEG TRANSPORT JAMEN
						•	SCREENPOINT Medical
•	•	•	•	•	•	•	SIEMENS Healthineers
		•	•			•	SIUI
						•	THERAPIXEL
•							VISARIS

Fujifilm · Synapse VNA



Highlights

Synapse VNA is core building block of an enterprise imaging architecture, with a vendor neutral best of breed approach; it is a secure, scalable, standard-based application allowing clinicians and caregivers to access any relevant clinical object. Focused on all medical data, DICOM and native non-DICOM objects. Synapse VNA can also enable easy upload of content from the desktop or mobile device, helping Healthcare providers to reduce silos and ensure all data is available when a clinician needs it.

Image Information Systems · iQ-Web RIS



Highlights

iQ-Web RIS is a flexible and modular web-based radiology information system, which can easily be customized to meet the needs of any hospital or imaging center.

- HTML5
- · Advanced web-scheduling
- Flexible cross-site reporting
- Extensive statistics
- Secured access

Image Information Systems · iQ-System PACS



HighlightsiQ-System PACS is an easily configurable, highly scalable picture archiving and communication system. It is installed in more than 10,000 facilities ranging from small, individual, imaging centers to large multi-modality, multi-site hospital installations across 118 countries. It is full-featured, state-of-the-art, robust and reliable, and available in most major world languages.

Image Information Systems · iQ-View



Highlights

iQ-View is the vendor neutral easy-to-use multimodality reading station that has been designed by radiologists for imaging specialists. A unique previous study management using artificial intelligence accelerates the diagnostic process by automatically presenting relevant previous studies of any modalities. iQ-View Pro automatically merges different patient identities from any PACS.

i-Solutions Health $\,\cdot\,$ RadCentre Cockpit & Speech Integration



Highlights

RadCentre is a comprehensive process and data management solution for radiology, nuclear medicine and radiotherapy. Based on latest technologies it offers high usability with an innovative user interface (Cockpit) and most efficient reporting with integrated speech recognition.

- Integration of received reports (specification depends on cooperating system)
- Fast and efficient creation of reports for treatment without delay


ITZ Medicom · ITZ Hyper.PACS with archiving-system Hyper.ARC



Highlights

- Fast, stable and safe PACS
- Reliable software solutions
- One frontend for the data
- One database for all data
- Fireproof longterm-archive
- ITZ-parallel-archiving-concept; no archiving of errors like with backup-principle
- Easy to support
- Compliant to RöV and MDD/MDR Class IIb

ITZ Medicom · ITZ Hyper.PACS – Reporting & advanced visualization

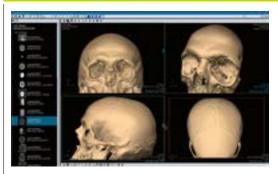

Highlights

- · Solution for all purposes with special, easy hanging protocols
- Selection of postprocessing software for Radiology and Cardiology
- Real-time viewing. LVA, QCA and 3D-high-end-postprocessing
- One surface for viewing, diagnosis and telemedicine
- Viewing-history, session-parking, MRT-space-time-presentation
- Unlimited and automatic lists for demo, science and presentations

medigration · ImageVision

- Mammo MR
- Screening Calcium scoring
- CFA
- Curonaries / heart
- \times Lung ☐ EP planning ☑ Functional Ir
- Functional Imaging
- Vessel measurement ☐ Virtual colonoscopy



Highlights

- Easy to use, high performance examination and analysis system for radiological routines
- · Access to all images (including previous images) within seconds
- Unique and hierarchical data compression without any loss
- Individually configurable hanging protocols
- Independent individual scaling of your interfaces

medigration · RIS/PACS



Highlights

Our RIS / PACS solutions are designed for multisite and manufacturer-independent networks. The WinRadiolog RIS product portfolio implies the whole patient management for your medical institution. Our PACS product portfolio comprises a proven DICOM archive, an intuitive operating reporting 3D ImageVision workstation, teleimaging and mobile solutions, patient CD system dosemanagement software solution..

Nexus/Chili · PACS



Highlights

- Multimedia PACS
- One viewer for all areas
- Scalable (practice to enterprise)
- Multitenancy
- Fail over and load balancing
- · Archiving in existing systems
- Interfaces and synchronisation with HIS/RIS
- Supports multiple IHE workflows
- Referring physician access
- Teleconferencing
- Consultation
- · Portal functionality

Nexus/Chili · Web



Highlights

- · Multi-media (DICOM, JPEG, MPEG, PDF, ...)
- · Very well suited for teleradiology
- · Referring physician access
- User concept with roles and rights
- · Central user administration (LDAP, AD)
- · Security measures
- Data compression (lossy & lossless)
- Suited for reporting (MPG class IIb)
- · Works with any PACS

Nexus/Chili · Nexus/RIS



Highlights

- Modern and intuitive user interface
- Scheduling and resource management
- Seamless integration with all our radiology products, e.g. PACS and portals
- · Context-sensitive integration of 3rd party solutions, e.g. speech recogni-
- tion, structured reporting and dose management
- Integration server for the management and monitoring of DICOM or HL7 interfaces
- Business intelligence tools


Nexus/Chili · Diagnost



Highlights

- Independent of modality
- CT, MR, CR, DR, PET, PET-CT, US, XA, . . .
- Mammography
- · Radio therapy
- · Powerful hanging protocols
- Independent of OS
- Integrated teleradiology
- Extensible by other applications
- HIS / RIS integration
- Consultation functionalities
- Teleconferencing

Nexus / Chili · Teleradiology Gateway

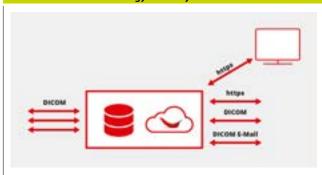

Highlights

- · Vendor-independent protocols
- DICOM, DICOM E-Mail, HTTPS,
- Rule-based autorouting
- Automatic recovery after interruption
- Comprehensive security measures
- Lossy and lossless compression
- Data encryption
- Audit trails
- Diagnostic web-viewer
- · Web-based administration
- Compliant to German RöV and DIN 6868-159
- · Works with any PACS

Nexus/Chili · Import PACS



Highlights

- PACS for external data from CD/teleradiology
- Temporary archive in addition to regular PACS
- · Manual web-based import
- · Automatic import with import robot
- Data reconciliation with own IDs (IHE compliant)
- Delivery to regular PACS
- Adjustable automatic data removal
- DICOM Q/R capable
- Works with any other PACS

Protec · Propaxx



Highlights

- Administrative and assisting functions, e.g. the integrated interface for reporting the clinical findings or synchronic viewing images
- Detailed 10-bit display of the X-ray images
- Support of US images /sequences (US,CT,MRI in Propax Vet)
- Configurable menu with guide access
- Individual system size: single or multiple workstations
- Individual system size as multi-user/ multi-client PACS solution
- · Integrated backup function

Protec · Conaxx 2



Highlights

User-friendly and intuitively operable software for the acquisition of X-ray images and operation of several DR-modalities and X-ray generators.

- Three clicks only to get your X-ray image
- Professional image tuning function
- Automatic image optimisation
- Image diagnose directly in CONAXX two possible (optional/single workstation solution)
- Compatible with any DICOM PACS
- Extraordinary workflow efficiency

Siemens Healthineers · syngo.via



Highlights

syngo.via is the intelligent, integrated imaging software, which offers multi-modality and fast 3D reading, innovative and Al-powered applications. It speeds up your routine and provides actionable imaging based results to enhance care delivery and outcomes.

- Simplifying routine streamlined reading and reporting with powerful tools and integrated reporting solutions
- Empowering innovation latest technologies and syngo.via open apps provide a gateway to innovations and boost your clinical capabilities
- Adapting to you integrating seamlessly into your IT environment and growing with all your medical and operational needs from workstation to multi-site

Siemens Healthineers · syngo.plaza



Highlights

syngo.plaza is the smart PACS for reading and reporting a large variety of cases – from routine to complex.

- Centerpiece robust performance, intuitive operation and intelligent reading tools
- Smart PACS 3D technology, powerful storage capacities and vendor-neutral archiving even enterprise-wide
- Lasting investment highly scalable long-term solution growing with your plans

Siemens Healthineers · syngo Dynamics



Highlights

syngo Dynamics is a single, enterprise wide, multi-modality intelligent reading and structured reporting platform to streamline data transfer and workflow. It helps enable high quality outcomes, efficient workflows, and improved operational efficiency

- Efficient workflows: Enables healthcare teams to quickly and easily access study data
- Consistent data: Quickly connect to the right data and avoid missing or conflicting data, to enable high-quality outcomes and faster reimbursement
- Simplify operations: A single platform that helps coordinate care across the continuum

Siemens Healthineers · syngo Virtual Cockpit



Highlights

syngo Virtual Cockpit, a software for remote scanning assistance, lets you make the most of your imaging devices. Medical staff can use this software solution to connect remotely to scanner workplaces to assist personnel at a different location, especially where more sophisticated examinations are required.

- Boost confidence by sharing in-house expertise
- Enhance patient satisfaction by improving availability
- Relieve cost pressure by enhancing flexibility

Siemens Healthineers · next generation VNA



Highlights

Enterprise-wide clinical image data management (IDM) made easy: Your data strategy to lead the future. IDM is the universal solution that meets the requirements of a powerful enterprise data management solution for managing, sharing and archiving clinical data independent of format and origin (DICOM and Non-DICOM). Scalable storage capacities allow data management across departments.

- Patient-centric storage
- Single point of integration
- Cost-saving data management
- Virtual patient CD app
- Universal zero-footprint enterprise viewer

Visaris · Avanse 360



Highlights

- Image archive and communication servers (PACS)
- Web-based radiology workflow system
- Diagnostic radiology workstations (reading stations)
- Advanced DICOM services including modality scheduling (MWL) and exam status monitoring (MPPS) servers



• Range of fully scalable server and viewing station extensions that can be tailored to specific needs of enterprise level healthcare providers at any time

Business Intelligence

Agfa HealthCare · Enterprise Imaging Business Intelligence



Highlights

Easy access to the information you need through standard and customizable reports. Your Enterprise Imaging solution contains a wealth of information about your healthcare enterprise and its operations. Agfa HeatlhCare Business Intelligence reports are a cornerstone in better understanding operational realities, identifying areas for focused improvement and help build efficiency gains.

Business Intelligence

Canon · Vitrea Advanced Visualization



Highlights

Vitrea Intelligence allows you to analyze the DICOM data generated directly by modalities.

- Gain accurate visibility into all modality utilization including CT dose tracking
- Discover and resolve unwanted variations in acquisition workflow
- · Better manage sites, departments and resources

i-Solutions Health · RadCentre Analytics



Highlights

RadCentre Analytics offers an integrated solution for specific data analysis and interactive reporting to increase performance in radiology.

- Predefined and high performant processing of operating figures
- Unlimited analysis options for optimisation of business outcomes
- Integrated data warehouse solution
- Visualization of radiation exposure extracted from PACS

Mammo Workstations

Fujifilm · Amulet Bellus II



Highlights

- Multi modality diagnostic workstation
- Tomosynthesis reconstruction for time saving image transfer
- Customizable GUI and workflow
- Report functionality and 3rd party report integration
- Can be integrated into existing environments
- Up to five clients

Hologic · Intelligent 2D imaging technology



Intelligent 2D Reveal more details in high-resolution 3D exams.



• Lower dose. Generate a new 2D image using advanced machine-learning algorithms and high-resolution

3D data to produce a robust, natural-looking image.

Improved patient experience.
 3.7-second scans reduce compression time to minimize motion and improve patient comfort.

* Compared to 2D.

Hologic · Quantra 2.2 Breast Density Assessment



Highlights

A New Standard for Automated Breast Density Assessment Discover the Difference with Quantra 2.2.

Now with texture and pattern analysis, Quantra 2.2

offers more consistent, more reliable scoring. It delivers the accurate information clinicians need to confidently customize personalized patient care decisions and provide the highest quality assessments.

- Breast Density Estimation Algorithm
- Texture and Other Features Calculations
- BI-RADS 5th Edition-like Breast Density Scores

medigration · MammoView

- Default display protocol
- Hi-Res displays
- or mixed setups

 Digital dictation
- integration
- □ Dedicated keypad
- ✓ WebClient



Highlights

- Extremely easy to use and manage
- Direct findings in the image
- CAD support (optional) and a second view area to examine US and MRT images
- Hanging protocols can be configured individually to automate your routine workflow
- Outstanding image quality (2,048 greyscale)

60

Mammo Workstations

Siemens Healthineers · syngo.Breast Care



Highlights

syngo. Breast Care is the advanced reading and reporting solution with powerful tools for efficient screening and comprehensive multimodality diagnostics.

- Choose the most suitable solution from a stand-alone workstation to a multiple-user server
- Customize your automated reading workflow to your personal preferences
- Easily include multimodality and 3D ultrasound reading, synthetic views, contrast enhanced mammography, breast density and CAD information
- Integrated CAD solutions with interactive decision support based on highly trained Al-based algorithms

Pathology

Fujifilm · Dynamyx



Hiahliahts

A vendor neutral digital pathology application which can be integrated with any lab information system or digital slide scanner. Supporting LEAN workflow and collaboration (including online sharing). It allows pathology departments to move to digital at their own pace and allows the integration of any scanner or Al vendor via an open API throughout the life of the solution. The mature platform was designed by pathologists for pathologists and brings all of the tools to enable a pathology department to digitise and introduce LEAN working with minimal disruption.

Portal Solutions

Canon · Vitrea Advanced Visualization



Highlights

- Suite of advanced applications provide full-powered solutions for 2D, 3D and 4D advanced visualization used to process and analyze clinical data from multiple modalities –MRI, CT, CR, DX, RG, RF, US, XA, NM, PET, PET/CT and SPECT
- Modular viewing platform that provides a broad range of clinical applications for Cardiology, Neurology, Oncology, Women's Health and MSK
- With multi-vendor support, Vitrea's broad range of clinical applications can be used to read data for all the major vendors' equipment


medigration · MultiPortal



Highlights

- To connect your referring practices · Efficient and encoded transferral of
- image data
- · Secure, user-defined access control
- No elaborate VPN neccessary
- Fast display of images and findings as PDF or SR
- For PC / MAC: Intuitive, web-based tool, to be launched without any installation via any standard browser

medigration · webConnect

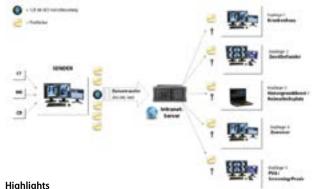

Highlights

- Uncomplicated exchange of image data via the internet
- Highly cost effective since only the actual transferred data is calculated
- · No VPN connection necessary
- Images and results can be called up within seconds due to intelligent data compression
- Total security by means of 256 bit AES encryption

Nexus/Chili · Telemedicine Record



Highlights

Web-based platform for the exchange of multimedia documents, e.g. diagnoses, lab results, DICOM images

- Capture, display and administration of patient data
- Forwarding to referring doctors
- Upload and download of DICOM and other images
- Inter-sector exchange of multimedia patient data
- · Multicentre studies with DICOM images

Portal Solutions

Siemens Healthineers · eHealth Solutions



Highlights

eHealth Solutions provides solutions and products to network health institutions, for centralized access to patient-related information and different facilities. With eHealth Solutions, our customers build an IHE-compliant infrastructure for exchanging medical information in cross-institutional, regional, and national eHealth structures. Built for sharing and collaboration, eHealth Solutions connects physicians, patients, hospital employees, referrers, and other parties and organizations and summarizes relevant data needed for reliable treatment decisions.

Siemens Healthineers · teamplay



Hiahliahts

teamplay is the departmental performance management solution for radiology. It runs on our awarded cloud-based platform with the goal to enable healthcare providers secured processing of data within or outside a hospital's network and serves as the entry point to access a variety of digital applications. The teamplay performance management applications unlock the insights of operations data and hereby help to make quick and well-informed decisions to streamline operations in radiology*.

* Please check if teamplay is available in your country.

Mobile RIS/PACS Viewers

Agfa HealthCare · Enterprise Imaging



Highlights

By seamlessly creating a comprehensive medical imaging record, and providing you with the tools to collaborate, exchange, view and manage it, Agfa HealthCare Enterprise Imaging supports you to build a system that will bring you clinical value all along the care continuum.

Agfa HealthCare · Enterprise Imaging Universal Viewer



Highlights

Patient-centric image access from across all specialties in the enterprise, with enhanced viewing, collaboration and sharing, on a single web viewer. XERO Viewer provides secure access to imaging data from different departments and multiple sources, in one view, to anyone who needs it. With the mobile device support, you can truly work on the go, capturing and uploading images wherever you are.

Agfa HealthCare · Image Exchange



Highlights

Fast, secure, reliable transfer of patient studies between hospitals, with no CDs or DVDs. With unlimited inbound and outbound uploading and downloading of images and a web-based way to share images with patients, referring physicians and other hospitals, Agfa HealthCare Image Exchange solutions provide the enhanced image sharing you need to improve the delivery of care while decreasing costs.

Image Information Systems · iQ-3DView



Highlights

iQ-3DView is a zero-footprint viewer for web-based 3D visualization from anywhere. It runs on both tablet and desktop computers without requiring client installation. Feature highlights include curved MPR, volume rendering including cropping and virtual endoscopic view.

Mobile RIS/PACS Viewers

Image Information Systems · iQ-4View



Highlights

iQ-4View is a ground-breaking diagnostic multimodality zero-footprint viewer, suitable for virtually all browsers and operating systems. It runs on almost any device (desktop computer, tablet PC or smartphone) and requires no installation on the client. iQ-4View allows reading, viewing or reviewing any kind of images, structured reports and Encapsulated PDFs.

ITZ Medicom · ITZ Hyper.mView – Mobile solution



Highlights

- Receive your images wherever you are with high image quality
- ITZ Hyper.mView supports all mobile devices and tablet-PC
- The solution is scaleable to your needs and budgets
- Secure by encryption and anonymized transmission
- Different functionalities from viewing up to diagnosis
- Administration from any location
- Several helpful measurements

medigration · MultiPortal





Highlights

- To connect your referring practices
- Efficient and encoded transferral of image data
- Secure, user-defined access control
- Fast display of images and findings as PDF or SR
- · No elaborate VPN neccessary
- For tablets & smartphones: Installation and updates easily via AppStore

Nexus/Chili · WebViewer^{NG}




Highlights

- · Mobile image viewer
- Teleradiology
- PACS administration
- Easy integration into any other system, such as HIS/RIS/PACS/EPR
- Works without an app store
- Independent of operating system (iOS, Android, . . .)
- Device independent (Apple, Google, . . .)
- No app but HTML5!
- Works with any PACS

Dose Management Systems

Agfa HealthCare · Dose

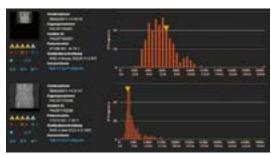

Highlights

The Web-based Dose Monitoring platform, integrates directly into existing picture archiving and communication system's (PACS) environment. Collecting the dose and metadata information already there, it can create patient radiation dose analyses at the study, patient, device, modality or institution level. It also provides all the tools you need for Root Cause Analysis, to help you understand and solve potential problems. The Dose Monitoring solution gives the tools needed to manage, analyze and balance the organization's radiation dose management.

BMS Informationstechnologie · EasyDoseQM



Highlight

EasyDose^{QM} liberates care professionals from most time consuming manual tasks: acquisition, documentation, analysis and archiving. It utilizes DICOM, HL7 and integrates seamlessly within existing HIS / RIS and PACS systems. Dose information about individual patients, modalities and departments can be obtained without complicating search mechanisms with a few mouse clicks. EasyDose^{QM} has been developed with the end-user in mind.

Dose Management Systems

BMS Informationstechnologie · EasyDose^{QM} – Options



Highlights

EasyDose^{QM} complementary options:

- · Integration of measuring stations and column scales
- RFID tracking of mobile devices, e.g. C-arms
- GPU based Monte Carlo Simulation
- · HIS/RIS/PACS integration

Bracco · Nexo [Dose]




Highlights

- Nexo [DOSE] supports compliance with imminent European Directive (2013/59/EURATOM)
- $\bullet \ {\sf Single-server}, fully \ {\sf automated} \ {\sf system} \ {\sf enables} \ {\sf enterprise-wide} \ {\sf data} \ {\sf acquisition}$
- Multi-modality, vendor-neutral software minimizes installation time and costs
- Customized e-mail alerts help improve control and implement the ALARA* principle
- Global Dose Registry for study dose comparison
- *ALARA (as low as reasonably achievable)

Fujifilm · Syncro-Dose

Highlights

Syncro-Dose is a comprehensive system for monitoring and managing patient radiation exposure across different imaging modalities. It is a support for the optimization of radiological procedures and acquisition protocols, a tool for supporting clinical audit and it provides

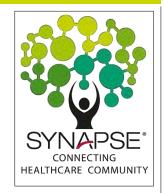

a comprehensive patient dosimetric history. General and specific dashboards track key performance indicators (KPI) to measure productivity, to achieve quality assurance and to support quality of care. SyncroDose is the radiation dose index monitoring system developed by Fujifilm, compliant with the directive 2013/59 / EURATOM of the European Union.

Guerbet · Contrast&Care



Highlights

Contrast&Care is a solution dedicated to contrast dose management. It connects to all Guerbet injectors and Hospital Information Systems (RIS, PACS, EMR...) and collects all relevant data about contrast media usage, patient history, and injector activity. Contrast&Care facilitates the traceability of contrast media and provides several tools that help imaging centers optimize contrast media consumption.

Guerbet · Dose&Care



Highlights

Dose&Care is a state of the art radiation dose monitoring solution, which allows documenting patient exams, understanding the reasons for excessive exposure and monitoring activities related to patient exposure. It provides the means to remain compliant with an ever-evolving regulation while improving the workflow and ensuring patient safety.

i-Solutions Health · RadCentre Dose View



Highlights

RadCentre Dose View is a stand-alone and RIS-independent dose management system to assess patient exposures due to ionizing radiation. The system is able to meet legal requirements (i.e. EU-Directive EURATOM 2013/59 and related national regulations for radiation protection) by offering consistent standards to increase the quality of radiological examinations.

Dose Management Systems

Medigration · Domako



Highlights

Domako. Simple software solution for dose management (DM). Collects, classifies and evaluates dose data; graphs them. Efficiently control DM process. Optimize protocols of modalities purposefully. Observes dose guidelines of BfS. Holistic/detailed, be it in terms of individ. protocols, pat. groups or individuals. Fulfils function of an autom. X-ray book. Enables to react proactively to deviations. Web-based on-premises system. Can be integrated into other software systems.

Siemens Healthineers · teamplay Dose



Highlights

teamplay Dose is an enter-

prise-wide radiation dose management solution that provides access to dose data and supports quality assurance processes for monitoring imaging radiation dosage - vendor and modality agnostic.

- · Analysis: Patient-centric dose monitoring
- Insights: Size-specific dose estimate (SSDE) and effective dose calculation
- Comparison: Benchmark against peers and national and institutional reference
- Report: Dose events and storage for future quality assurance

Accessories / Complementary Systems

Canon · Advanced Edge Enhancement



Highlights

Enhanced visibility of catheters, fine structures and bones

- Better visualization of foreign structures in the image
- Enhanced display of fine structures Better definition of the structures in
- soft tissue and low dose area's
- · Obtain enhanced images suitable for measurement or other applications
- Catheter, Small Structure and bone settings depending on the specific application
- Improved visibility of bone contours for easier measurement of length and angles

Canon · Scatter Correction



Highlights

Excellent image contrast without a grid. Canon's new image processing software Scatter Correction could reduce radiation dose by up to 60 percent on your radiographic examinations.

Where a grid physically reduces scatter

and thereby increases the image contrast, the software mimics this process virtually. The software works by creating a scatter model, which is subsequently subtracted from the image. The result is an image with reduced scatter and increased contrast. The software is available for Canon FPD imaging systems.

i-Solutions Health · RadCentre Technician Profile



RadCentre Technician Profile visualizes requested or performed examinations and reports at a glance and supports a fast and modality based workflow.

- · Specific icons show examination status or patient information
- Images of prior examinations via integrated PACS viewer
- Interactive icons to change information or workflow status
- Scanned document files and laboratory results

i-Solutions Health · RadCentre Quality Manager



RadCentre Quality Manager supports the justification and documentation process. It increases quality assurance, patient safety and efficiency of examinations and offers quick overview of information for doctors to initiate the justification.

- Integrated justification process
- Overview of non-validated examinations
- · Easy planning of examinations and specific information for technologists

Mammography



Mammo CT

compression



Tomosynthesis

Canon · Artemis		
Pixel size 83 µm	Scan angle 30°	Scan time 15 s
Highlights The Artemis mammograp provides both 2D (Full Fie 3D (Digital Breast Tomosy The system system support additional work-up exami procedures.	ld Digital) and nthesis) imaging. rts both screening and	

Fujifilm · Amulet Innovality					
Pixel output 50 μm / 100 μm / 150 μm	Scan angle 15°/40°	Scan time 4s/9s			
	- 1				
Highlights	- 1	1			
Unique Fujifilm develope using hexagonal pixels for of different resolutions					

 Dual angle tomosynthesis for dose efficient with maximum diagnostic performance
 Ergonomic design for user and patients
 Dynamic image processing with advanced options like fine structure correction FSC

• Tomosynthesis biopsy, vertical and lateral

and dynamic visualization II

for biopsy procedures
• DICOM 3.0 connectivity

with CAD software

• Optional diagnostic workstation available

approach

* Breast CT is an alternative breast imaging modality. Classified as "mammography" here for editorial purposes only

racy and speed (direct conversion)

11/11

IMS Giotto – GMM Group · Giotto Class						
Pixel size 85 µm	Scan angle 38.6°	Scan time 13 s				
ου μιτι	30.0	133				
		6				
Highlights						
Giotto Class is an advance	d and innovative					
three dimensional breast	imaging technology	18				
able to perform	100	4				
Digital mammography e						
Breast Tomosynthesis (36 Conthesis of 3D income.)						
Synthesized 2D image godataset	enerated from 3D					
Combo: Tomosynthesis 8	& digital mammo-					
graphy						
Stereotactic biopsy in pro	one or upright					
position						
	tral Mammography (CESM	1)				
IMS Giotto is a company o	ot GMM Group					

Intermedica	l · Wom	am		
Pixel size 85 µm		Detector size 24 × 30 cm	Detector type a-Se/a-Si	
of effective b	ersions: ana d 3D digital ction al modality: reast densit ompressed	logic,		

Planmed Oy · Clarity 3D						
Pixel size 83 µm	Scan angle 15°	Scan time 13 S				
imaging, diagnostic ima and Digital Breast Tomo • Continuous Sync-and-Si	noot tomosynthesis imagi construction and Tomo- nable sharp and artifact	5				

Tomosynthesis

tomosynthesis

Siemens Healthineers · 50° Wide-Angle T Scan angle 50° 85 µm Highlights Advance screening and diagnostic results with high accuracy. 50° Wide-Angle Tomosynthesis has proven an increase in cancer detection rate of 41.5 % for invasive cancer with a one-view tomo scan only. · Highest depth resolution with 50° Wide-Angle Tomosynthesis • Gain a fast overview – with our synthetic visualization Insight 2D • 40 % dose reduction as opposed to FFDM as an adjunct to tomosynthesis · Decrease tomo reading time - with our unique, synthetic visualization Insight 3D a unique, rotating 3D display in breast

Tomosynthesis	Siemens Healthineers · Mammomat Revelation				
Scan time 25 s	Pixel size 85 μm	Scan angle 50°	Scan time 25 s		
	Highlights State-of-the-art digital m system for screening and • Make anatomical details with our unique 50° Wic in Tomosynthesis and Bic • Automated breast densi right at the acquisition of instant risk stratificat • InSpect – our integrated facilitates the immediate biopsy directly at the system of the syst	diagnostics s clearly visible de-Angle – opsy ty measurement vorkstation allows ion d specimen scanner e control of the stem information fast with			

Villa Sistemi Medicali · Melody IIID TS **Scan angle** 15°/24°/36° Scan time Pixel size 85 µm 5s/6s/9s Highlights • Digital mammography system with tomosynthesis function · Special anti-scatter grid for tomosynthesis allowing superior image quality • Collimator with automatic recognition of compressor paddle • User-friendly interface with touch screen display on each side of C-arm • DICOM 3.0 connectivity

Digital Mammography **DMS Imaging** · Serenys DR Bym Detector size Detector type 18×24cm/24×30cm 85 µm Highlights • The Serenys DR Bym, with the added advantage of an isocentric C-arm including stereotactic biopsy • The isocentric C-arm can be fully motorized and permits all breast projections without moving the patient and without adjusting the height of the C-arm, making exams faster and more comfortable • The newest version dbt offers the possibility for tomosynthesis • The device is also available in analogic or in compact version. • kV Range: 20 - 35 kV

IMS Giotto – GMM Group · Giotto Class Smartfinder					
Pixel size 85 µm	Detector size 24 × 30 cm	Detector type a-Se			
Highlights Giotto Class is a patented tomosynthesis system offa multitude of diagnostic interventional solutions, in greater to the system or upright position using specific prone table access	ering and nclud- prone the	A			

• Optional diagnostic workstation available with CAD software

• Detector size: 24 × 30 cm

• High precision tomo guided

technique and tomo biopsy

• Combination of traditional stereo

· Integration with accessory for real-

time acquisition of biopsy cores

biospy

imaging

The compact design allow the operator to use the system in the same room for both diagnostic and interventional procedures

IMS Giotto is a company of GMM Group

IMS Giotto – GMM Group · Giotto Class 40000					
Pixel size	Detector size	Detector type			
85 μm	24 × 30 cm	a-26			
Highlights The system is designed to drastically improve the screening and diagnostic throughput thanks to an high rotation speed and an improved vertical run speed The gantry is ergonomically designed to give patients a natural and more relaxed positioning The operating and interventional modalities include: • Digital mammography examinations (2D) • Breast Tomosynthesis (3D) • Synthesized 2D image generated from 3D dataset • Combo: Tomosynthesis & digital mammography • High precision tomo guided or stereotactic biopsy • Contrast-Enhanced Spectral Mammography (CESM) IMS Giotto is a company of GMM Group					


Digital Mammography







Siemens Healthineers · Mammomat Fusion Detector size Detector type 83 µm Highlights Premium mammography system to enhance everyday screening and diagnostics · Help your patients to relax with the Mood-Light option • Stereotactic biopsy option for fast seamless procedures • New generation Csl detector technology for higher spatial resolution at low dose • Refined workflow to perform complex tasks at the click of a button • Personalized OpComp and OpDose · Focus on total cost of ownership including operating costs and service

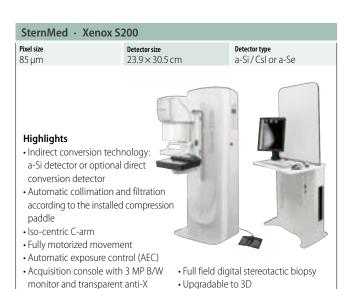

Stereotactic biopsy device

protective barrier



Mammo DR

Examion · X-Mammo DR M Detector size Detector type 76 µm $18 \times 24 \text{ cm}$ Highlights

Retrofit analogue and CR mammography systems. Get the advantages of a digital detector without replacing your existing mammography equipment.

- Pixel size of 76 µm for high image quality
- Image processing developed specifically for digital mammography
- Efficient work flow and high operating comfort

Biopsy Units

IMS Giotto – GMM Group · Giotto Flexitable				
Pixel size 85 µm	Detector size 24 × 30 cm	Detector type a-Se		

Highlights

Flexitable is an accessory which, in combination with the Giotto CLASS system and the Smartfinder biopsy kit, enables interventional prone biopsy procedures

- · High manoeuvrability, thanks to its reduced weight, the handle and the special wheels
- Excellent ergonomics for the patient thanks to the possibility of adapting the position of the breast and inclining or raising the front end of the table to compensate for bending and come into closer contact with the chest



Film-Screen Mammography Planmed Oy · Sophie Classic S Filter 20 – 35 kV Mo/Rh Highlights • An optimum performance analogue system with ease-of-use · High technology and smart solutions enable efficient workflow at smaller clinics and busy screening centers. Also diagnostic mammograms are available by upgrading the unit with a magnification platform and stereotactic biopsy device · Automated features, such as the Automatic Exposure Control and automatic selection of the collimator ensure correct parameters for imaging and help to avoid additional retakes

Siemens Healthineers · Mammomat Select

23 – 35 kV	Mo Mo	Mo/Mo or Mo/Rh
		in the
Highlights		

An analogue system that is easy to use, provides images at the right dose and is cost-effective to offer women the standard of care they need

- Easy touch screen control for streamlined workflow
- · Easy to dose right with AEC control
- · Easy to invest with flexible service and upgrades

Villa Sistemi Medicali · Melody III

Optional CR interface available



• Available with $18 \times 24/24 \times 30$ cm bucky or special potter accepting both

Accessories / Complementary Systems

GCTechnology · CIRS Phantoms



BR3D tomosynthesis breast phantom Mammographic accreditation phantom

Highlights

- · Multi-modality breast biopsy and sonographic trainer (CT, US, MR)
- Stereotactic needle breast phantom
- · Mammography test tools
- · Mammography phototimer consistency testing slabs
- Digital mammography phantoms
- Mammoview markers

cassette sizes

Accessories / Complementary Systems

I.A.E. · C340



Hiahliahts

- Water cooled mammography tube unit for beam scanning mammography equipments, high patients throughput screening applications
- Brass body lead free X-ray shielding internal pump for oil circulation improves oil to casing thermal Exchange
- · Water cooled jacket avoids remote oil circulation
- Compact lightweight structure
- 800 W continuous dissipation for high energy techniques, high patients throughput

I.A.E. · XK1016T-400W



- Highlights
- Rotating anode mammography X-ray tube, with special bi-angled target, for optimal performances with all techniques
- Two separate focal tracks, small focus on 10° and large focus on 16°, optimal resolution performances
- Reduced thermal stress on the bearings improves tube life duration
- Severe tests during conditioning assure best performances
- Compact light weight structure

Leoni · Cable systems



Highlights

Leoni cable harnesses and systems for medical devices integrate a wide variety of functions: besides energy supply and control, they can also transmit light and data up to including so-called assist functions, such as the cooling of device components. Leoni provides wiring services for mammographs and can collaborate in the development phase.

As a systems supplier, Leoni provides consultation and analysis at local level, development and design of the optimum cabling and the supply of prototypes. In addition, Leoni offers individual and system verification testing, as well as customer-specific logistics solutions.

PTW · Normi MAM – Digital X-Ray test object



Highlights

- Checks all relevant parameters of digital Modularly composed CE marked mammographic X-ray installations
- Fully complies with DIN 6868-162 and DIN 6868-14
- class I certified test object
- Incl. different absorbers and test elements

PTW · Nomex Multimeter



Highlights

- Dosimetry system, CE marked, class Ilb certified, acc. to IEC 61674
- Single exposure captures all dose values, kVp, time, TF, HVL, frequency, pulses and waveforms
- Angular independent for positioning within the beam
- Fully automatic adjustment
- Data and waveform export to Excel
- Ideal for tomosynthesis measurements
- · Accessories: Normi MAM test objects

Varex Imaging · B-121 – Mammography housing



Highlights

The B-121 is an air cooled mammography housing that fits a standard threeinch X-ray tube insert meant for digital and tomography applications. The housing has two shroud configurations; with and without quiet D/C fans. The B-121 offers 300 watts of continuous heat dissipation with fans, which is approximately 200% greater than standard mammography housings.

R/F Film-Screen



Bucky



general-purpose radiography, being highly accurate and efficient. It is possible to expand the original system to meet the particular clinical requirements of the user. When the system is combined with a portable FPD (35 \times 43 cm) and digital processor, a wide range of applications can be performed.













Bucky



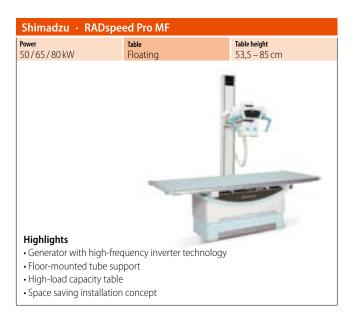





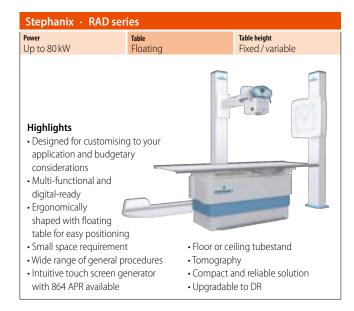



Fluoroscopy



Canon Plessart VIVO is a remote control R/F system comprising an R/F diagnostic table with an over-table X-ray tube configuration, an X-ray high-voltage generator, and a digital imaging system. This system is intended for use as a general-purpose system for abdominal angiography, general abdominal radiography, general skeletal radiography, support of endoscopic procedures, etc.

Canon · Plessart EX8 II format CCD-matrix 80 kW $1k \times 1k$

Highlights


The Canon Plessart EX8 is a digital remote control R/F system comprising a R/F diagnostic table with an over-table X-ray tube configuration, an X-ray high-voltage generator,



and a digital imaging system. This system is intended for use as a generalpurpose system for abdominal angiography, general abdominal radiography, general skeletal radiography, support of endoscopic procedures, etc.





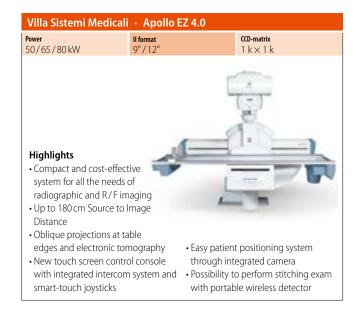

Makila V D

Mobile X-Ray			
DRgem · Jade			
Power 4 kW	Operation Mains	Motorized No	
		-	
Highlights		===	
,	pt: Portable radiography syster	m	
	powerful design nd intuitive operation		
	C (Free voltage) input	112	
• 40 ~120 kV, 1			
• Includes man		2015	
 Three-way co (Main body company) 	ntroi ontrol console, and remote cor	ntrol)	
,	ed APR data and user-program		
, ,	sible mobile stand with extern	40-	
	nterface, with bluetooth or DR	interface	
options			

Mobile X-Ray













Mobile X-Ray



Power 32 kW	Operation Battery	Motorized Yes
	Î	5-10
powered • Exposures are power supply	oile unit, battery cossible without e unit to an external ture and flexible	
• ± 320° rotating scopic arm	fι g adjustment through • k\	rontal bumper with anti-collision unction V Range: 40 – 125 As Range: 0.1 – 320

Villa Sistemi N	ledicali · Visitor T3	OR .
Power 32 kW	Operation Mains	Motorized No
care units as we pediatric or surg • Compact design verability of the • ± 90° arm rotati	on for increased by tube positioning mode	High performance generator and double focal spot (0.8/1.3 mm) tubehead kV Range: 40 – 125 • mAs Range: 0.1 – 220

Index of Advertisers

Canon	back cover / 156
Cefla/NewTom	17
Dunlee	21
ECR	79
Fujifilm	inside front cover / 2
Guerbet	35
Hitachi	3
Konica Minolta	93
Medtron	27
Mindray	129
Protec	47
Siemens Healthineers	inside back cover / 155
SIUI	135
Varex	4
Villa Sistemi	95

Published by:

Imprint

EUROPEAN HOSPITAL Verlags GmbH

Theodor-Althoff-Straße 45 · 45133 Essen · Germany tel +49 201 871 26-851 · www.healthcare-in-europe.com

 $\textbf{Radiologie Magazin} \cdot \textbf{Guido Gebhardt}$

Adalbert-Stifter-Weg 2b · 85661 Forstinning · Germany tel +49 151 15 67 48 33 · www.radiologiemagazin.de

Editors-in-Chief: Guido Gebhardt · Daniela Zimmermann

Executive Director: Daniela Zimmermann

Advertising:

Ralf Mateblowski (D, A, CH) · Eric Jund (E, F, I) Simon Kramer (B, GB, L, NL) · Gavin Hua (China) Jane Park (Korea) · Hanna Politis (U.S.A, Canada) Charles Yang (Taiwan)

Art Director: Christoph Muschiol

Agentur für Werbung, Beratung & Organisation \cdot Kiel \cdot Germany www.muschiol-online.de

 $\textbf{Printing:} \ \mathsf{Schaffrath} \ \mathsf{Medien} \cdot \mathsf{Geldern} \cdot \mathsf{Germany}$

Subscription:

 $Dorothea\ Fleischer \cdot fleischer @european-hospital.com$

Subscription rate: \in 19.– plus postage

The information and opinions expressed in articles and product entries published in RADBOOK are solely those of the manufacturers / companies, their authors and contributors, for which the publisher holds no responsibility. Errors and omissions excepted.


Disclaimer:

All company, brand and product names in this publication are the property of their respective holders. Not all products are available in all European countries

Please note that the manufacturers' websites may contain further product disclaimers.

© 2020 by EUROPEAN HOSPITAL Verlags GmbH. All rights reserved.

Mobile X-Ray

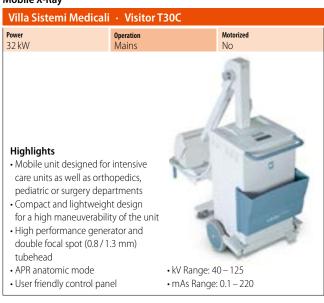



Accessories / Complementary Systems




I.A.E. · C20



Highlights

- · A new compact lightweight housing, specifically designed for mobile equipment.
- A low weight, less than 8.5 kg, combined with compact dimensions,
- 116 mm diameter and 342 mm length, allows significant reductions in the equipment supporting structures.
- A range of tube inserts up to 54 kW peak radiographic power at high rotation speed is available for this unit.

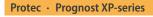

Model 903 radiography fluoroscopy QA pantom



- Mobile patient table to position the patient directly above the corresponding image receptor
- For digital DR detectors or with bucky tray integrated
- Fixed table height or elevating with
- floating carbon fibre table top
- · Elevating versions with line connection or battery powered
- Power: Line or battery
- Table: Fixed or adjustable height (optional), carbon fiber table top

Protec · Bucky series




Highlights

- Outstanding compatibility with X-ray tables, wall stands and U-arm of various brands
- High cost effectiveness due to continuation of use of existing grids and AEC chambers
- · All established detector types are supported
- Suitable for cassettes / detectors of different dimensions
- Perfectly prepared for simple realisation when upgrading an existing analogue system to a fully digital DR

EUROPEAN CONGRESS OF RADIOLOGY

ECR 2020 Vienna, July 15–19

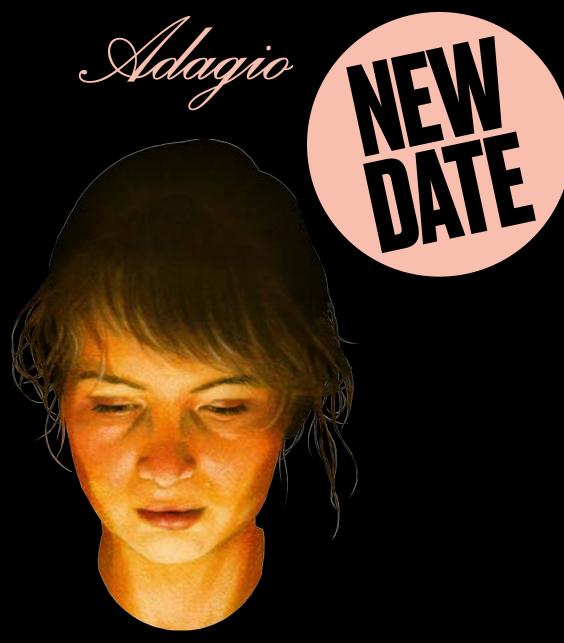

REGISTER NOW!

myESR.org/registration



R/F Digital





CR







- Customer-chosen optimal workflow
- Robust, yet easy to install and maintain
- Suited for mobile applications
- Networking capabilities deliver seamless integration
- Cassette size: From 15×30 cm to 35×43 cm, incl. mammography



Agfa HealthCare · DX-M*

Agfa HealthCare · DX-M* Slots Capacity Pixel size 6,7 – 20 Pixel/mm

Highlights

DX-M: Mixed to perfection

- Next-generation CR digitizer
- NIP and PIP detectors for general radiography and mammography
- Superb image quality and potential for dose reduction
- Five cassette drop-and-go buffer
- $\bullet \, \mathsf{Small} \, \, \mathsf{footprint} \,$
- Musica Image Processing
- Cassette size: From 15×30 cm to 35×43 cm, incl. mammography
- * DX-M with CR Mammography application is not available in the US





CR









Agfa HealthCare · I	DR 400 (floorm	nounted)	
Power 40/50/65/80kW	Detector type Csl/GOS	Pixel size <150 µm	
Highlights Cassette size bucky car from landscape to port Build-in Dose Area prod	rait	Tout Tr	

options for most needs • Supports CR and DR integration


• Flexible configurations and

• Scalable, flexible and affordable

meter (optional)

modality

- Requires limited space (4×2 m)
- "Musica processing" provides superior contrast detail and consistent, exam-independent image quality
- NX acquisition workstation offers comprehensive functionality for integrated workflow

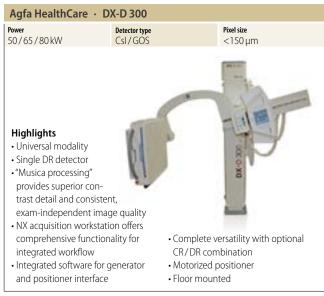


• Motorized Z-movement for easy

positioning

· Floating tabletop

DR



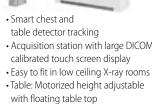








Canon · Easy DR

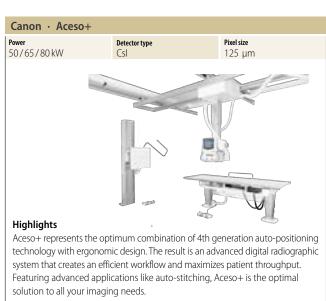



Our unique ceiling wagon provides unsurpassed usable stroke for high ceiling


heights while the CUBE solution allows the system to be installed in rooms


with ceilings as low as 2,5 m without requiring ceiling support infrastructure.


Pixel size Detector type Power 50/65/80 kW 125 µm Highlights Versatile solution for multipurpose examinations Multipurpose floor mounted X-Ray system • Suitable for mobile installations (i. e. truck or container) • Retractable anti-scatter grid · Vertical and horizontal positioning of the U-arm • Acquisition station with DICOM calibrated touch screen display • Asymmetrical diaphragm, specially designed for Thorax examinations • Optional low-power consumption X-Ray Generator






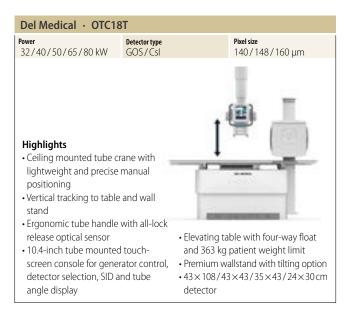

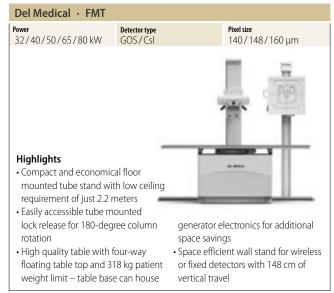

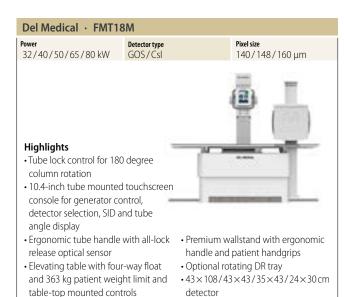


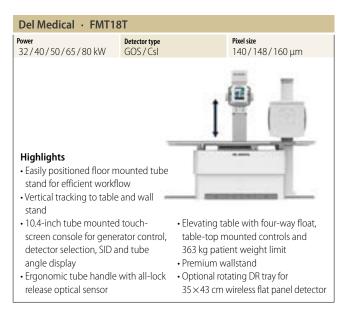

floating table top and 318 kg patient


weight limit - table base can house


DR

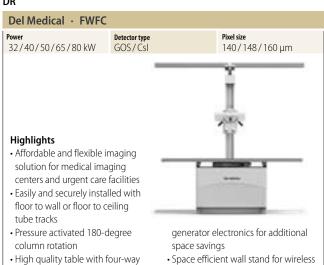

or fixed detectors with 148 cm

vertical travel range

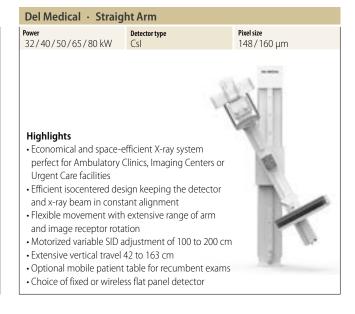

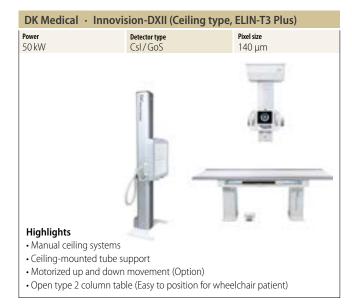





• User convenience with APR






- Conbined system consisting of tube support, table and generator (space saving)
- \bullet Bucky tray following in same direction with tube stand movement
- · Collimator turining on automatically
- User convenience with APR
- Standing knee position (Enable users to take images more convenintly without any hopital tool.


DRgem · Essential Ceiling System (GXR-SD Series) Detector type GOS/CsI Pixel size 32/40/52/68/82kW Highlights · System concept: Highly customizable DR ceiling system · Fundamental solution at affordable price • Ergonomic design for smooth movements and optimized workflow • Intuitive direction movement indicator and user-friendly interface · Options include AEC, carbon tabletop, · Elevating or floating table with high dual speed rotor and premium upgrade patient load up to 300 kg · Detector format: 17×17 " / 17×14 " wired/wireless • Integrated lock function








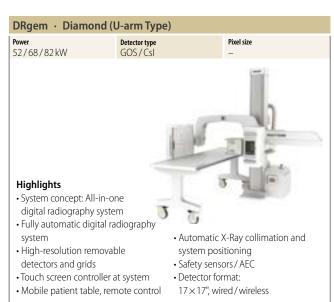

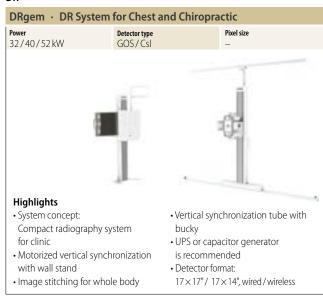








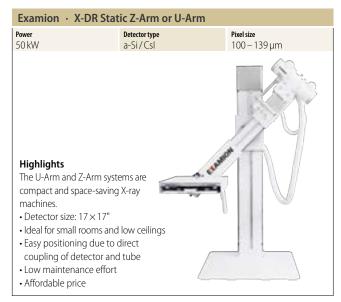





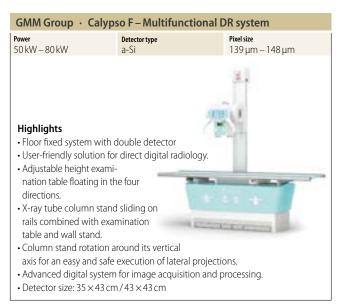

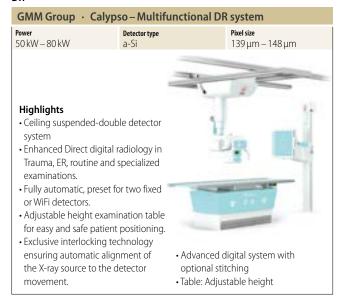

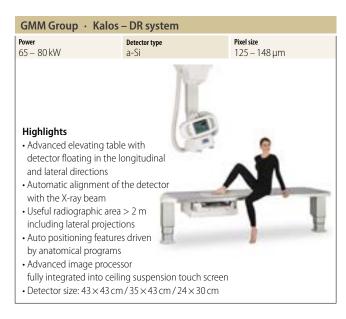


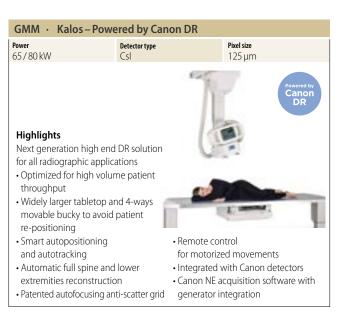


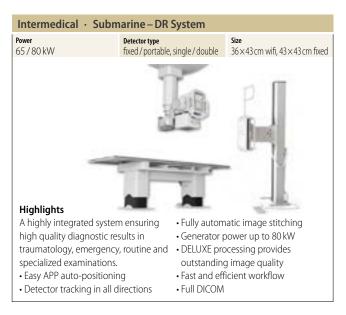







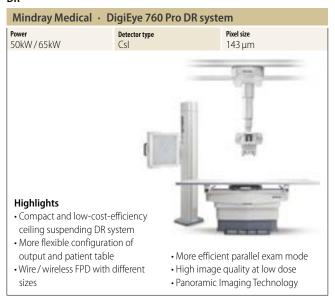

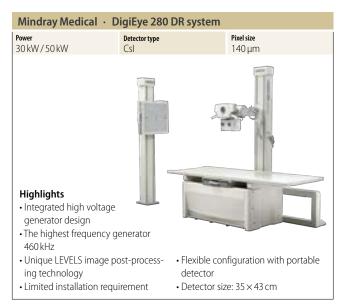



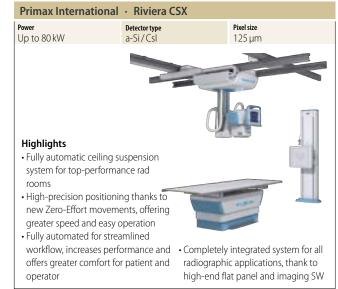

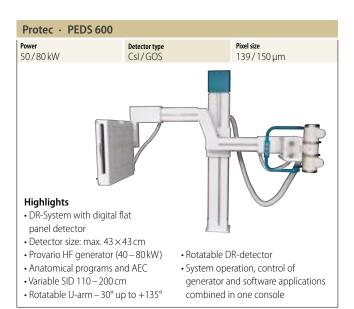







- Automatic coupling device to center tube and bucky
- Including wall bucky stand; stitching as optional solution
- Floating carbon fibre table top
- Adjustable height combined with undertable generator
- Fully digital DR-System with flat panel detector technology, different configurations from single to dual detector systems


touch-screen



system directly

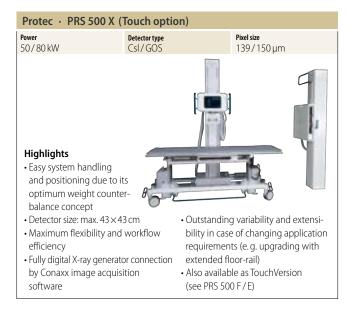

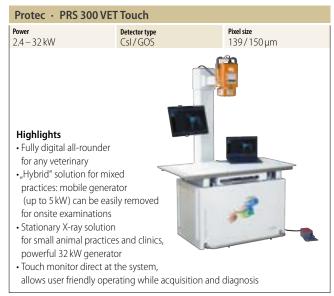


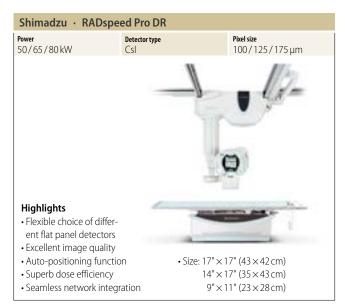








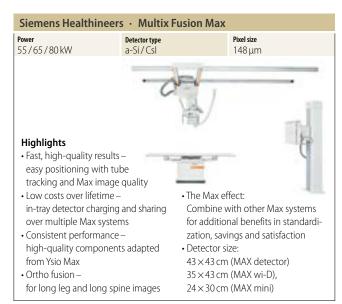





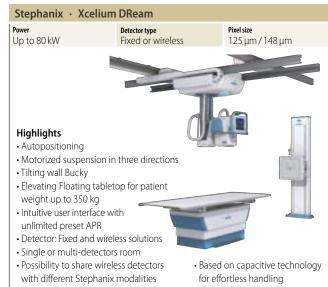







· Wireless remote





tabletop, on wheels









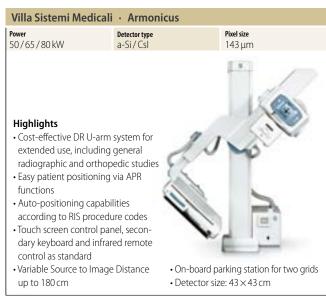

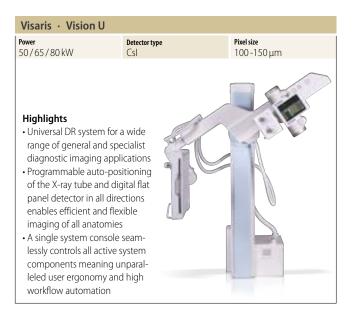



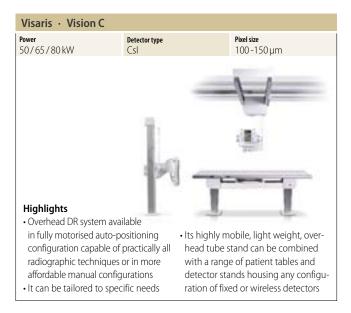

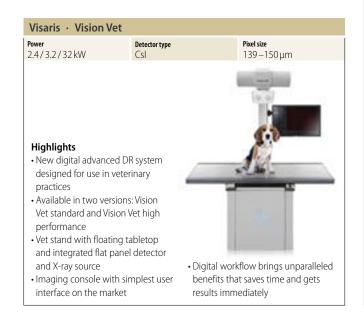





Portable DR



Design: DR upgrade mobile in 2 minutesSystem concept: 2 components

Examion · X-DR Portable Size Detector type a-Si / Csl Pixel size 150 μm

Highlights

Portable case solutions for emergency X-ray. All functions for acquisition, diagnosis and archiving on a single mobile PC.

- Wireless Digital X-Ray
- · Excellent image quality
- Patient administration with mini-PACS
- Radiological viewer
- Synchronization with stationary image archives (optonal)
- Detector size: 14 × 17"
- Pixel size: 150 µm

SIUI · SR-1000 Touch		
	Detector type a-Si	Pixel size 150 μm
Highlights - SR-1000 Touch is the port DR system integrated with X-ray generator, flat panel and the workstation	n the	
Eco-friendly & Safe:The gecaters to the clinical demalower power and low radiintensity Battery power supply syst A full charge enables over exposure at full power and 300 times at average power.	ands with Wation getween: as 200 times dover Wation Wation Wation Graph was as 200 times dover with the wation with the wation was as a second with the wation was as a second with the wation was a second was a second with the wation was a second was a second with the wation was a second was a s	IFI connection technology: IFI connection among the X-ray enerator, flat panel detector and orkstation to control exposure as well data transfer, which largely simplify e operation and improve efficiency ifferent stands facilitate wide oplications

Villa Sistemi Medicali · ArtPix EZ2GO

Highlights

- Plug-and-play solution for immediate upgrade to digital radiography
- Lightweight and portable acquisition system based on Wi-Fi Flat Panel detector and tablet
- \bullet Extreme flexibility and ease of use thanks to wireless connections
- Multi-use solution for shared use with general radiographic systems and mobile units
- Powerful acquisition software complete with post-processing tools and DICOM functions





- "Musica processing" for excellent contrast detail and exam-independent, consistent image quality
- Convenient size and light weight provide optimal convenience and portability, for all Genrad examinations
- Easier patient positioning enhances workflow
- High DQE and optimal pixel size, for low dose examinations
- Improved workflow and examination speed

Agfa HealthCare · DR 14s detector			
Size 14×17"	Detector type Csl		Pixel size 125 µm
			- 7
U:abliabte	1		
• Extremely long battery a	utonomy of	***	-
up to eight hours "Musica processing" for e contrast detail & exam-in consistent image quality	ndependent,	9	t, small, high resolution Exposure Detection (AED)
 Choice of Cesium lodide Gadolinium Oxy-Sulphic detector scintillator 	de (GOS)	• High DQE & low dose ex	
 Improved workflow & exspeed 	amination	Easy to clearSize: 46 × 38	



• High DQE & optimal pixel size, for low dose examinations • Extremely long battery autonomy of up to eight hours • "Musica processing" for excellent contrast detail & exam-independent, consistent image quality • Size: 251.0 × 314.5 mm (10 × 12 inch) (Effective area)









- Our proven advanced fine CsI/Tl and direct deposition technologies provide high DQE and better resolution
- Unique moisture-proof sealing method provides an extremely reliable CsI/TI screen that is protected from degradation
- High speed & low-noise ROIC prvide low-noise and real time image

Pixel size Detector type $35 \times 43 \,\mathrm{cm}$ 140 µm



- · Portable flat panel detector
- Our proven advanced fine CsI/TI and direct deposition technologies provide high MTF and excellent resolution
- Unique moisture-proof sealing method provides an extremely reliable CsI/Tl screen that is protected from degradation
- Standard cassette size
- Prompt display of preview/full images and the short cycle time
- enable fast image acquisition • Compact and lightweight for easy
- DC power input type is selectable

Canon · One Shot Long-Length (optional) Software

Highlights

One shot long-length exams enhance efficiency compared to conventional stitch exams; shorter examination time lower risk on patient movement, reduced dose and increased image quality.

- · Patient positioning stand with motorised height adjustment for easy positioning
- · Fixed installation or mobile for convenient reallocation
- · Large, ergonomic grip rails for confident patient positioning





- Optional grid
- Abiliy to use three existing detectors for cost-effective one shot longlength imaging
- · Versatile configuration; use either $3 \times 43 \times 42$ cm (410 CW) or $3 \times 35 \times 43$ cm (710 W) wireless detectors


Canon · CXDI Control Software NE

Highlights

CXDI control software NE is made exclusively for use with Canon digital radiography systems. This software helps to

optimise workflow and reduce the procedure steps needed to complete exams

- Instant viewing of high quality images
- Optimised workflow with minimum operation steps
- Interactive GUI for intuitive operation
- Single and prepacked protocols
- Emergency study capability
- Suspend exam / Reject analysis
- Automatic forwarding rejected images to a designated analysis workstation
 - · Automatic image stitching included
 - Scatter correction software (optional)
 - Advanced edge enhancement software (optional)

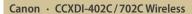

Pixel size Detector type $35 \times 43 / 43 \times 42 \text{ cm}$ 125 µm



Highlights

Wireless flat panel detector range

- Durable and ergonomic shaped wireless detectors
- Ergonomic detector design for easy hold, easy handle and easy position
- Dust- and water proof (IP55)
- · Optional docking station for detector check-in, detector battery charging and image transfer
- · Equipped with last image hold for secured image transfer

Canon · CXDI-410C/710C/810C Wireless


Pixel size Detector type 27.4×35/35×43/43×42 cm 125 µm



Highlights

Wireless flat panel detector range

- · Ultralight wireless detectors
- Increased durability by strong carbon fiber construction techniques
- Ergonomic detector design for easy hold, easy handle and easy position
- Dust- and water proof (IP57)
- · Docking station for detector checkin, detector battery charging and image transfer
- Equipped with on-board memory where 99 images can be stored (in stand-alone-mode)

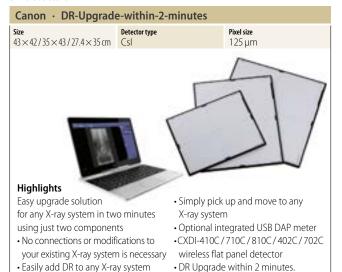

Freedom within reach

· Vet software available

 $17 \times 17'' / 17 \times 14''$, wired/wireless

Detector format:

Del Medical · Delworks EDR		
Size 43×108/43×43/35×43/24×30 cm	Detector type GOS/CsI	Pixel size 140 / 148 / 160 μm
Highlights • Powerful and user-friendly retrofit DR system • Choice of wireless detectors E17C fixed: 43×43 cm E14C wireless, premium: 35×43 cm E14CE economical: 35×43 cm E2430: 24×30 cm Delworks LLI monolithic: 43×108 cm for long length imaging	• Delworks Fl workstation	IT mobile tablet-based n option for ultimate
Single workstation for image display and processing and integrated generator control	• Cyber secu	re and HIPAA compliant 0 operating system



using just two lightweight components

a fully digital radiography system

(AED function included)

• Easy to apply to any X-Ray generator









• Robust, IPX6 waterproof, carbon monocoque housing


• Full image acquisition within four seconds only


• Charging time of only 13 minutes

AeroSync





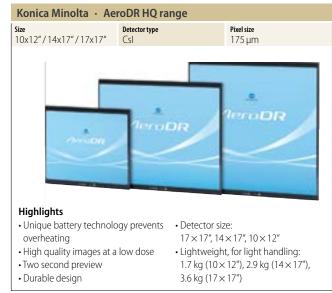

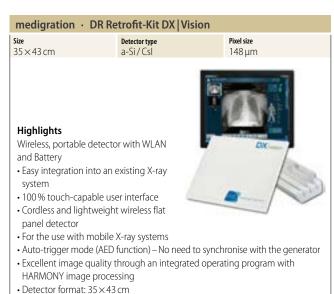



• Lightweight: < 3.0 kg

once the X-ray source starts the emission

• Enhanced productivity with DICOM functions



• Detectors are ISO 4090 compliant, existing Bucky can be used for DR Retrofit

• 16-bit dynamic range and high DQE for excellent image quality in three seconds

• Just one flatpanel required for integration into bucky table + wall stand



DICOM connectivity

• Shareable solution with other Stephanix modalities



• System equipped with battery charger and two batteries as standard










Del Medical	· MDR		
Power 40 kW	Width 57.6 cm	Weight 435 kg	
radiography so • Compact design column • Fully integrated station with characters	gn with fully collapsible d DELWORKS DR work- noice of detectors		
 Effortless maneuverability, allowing navigation through tight spaces Motor assisted inching from the tube Convenient storage for wireless detector, grids batteries, wipes, and lead apron 			
head • 35 × 43 cm / 24 × 30 cm detectors • On board detector charging • Pixel size: 148 µm			









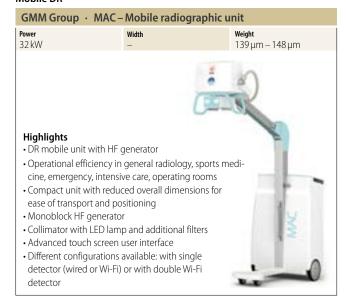

Intermedical · Com	oact DR Plus	
Power 32 kW	Width 57.6 cm	Weight 412 kg
Highlights Motorized mobile unit, ba powered, easy to handle a Telescopic arm • Wide choice of available • Full DICOM connectivity • 19" touchscreen user frie • Available in analogue ve	nd operate detectors ndly interface	







Primax International · Raybow XE				
Power 20/32/40/50 kW	Width 54 cm		Weight 520 kg	
Highlights Raybow XE redefines the of-care imaging Design centered on User delivering great improve reliability, image quality Effortless maneuverability	e concept of point r and Patient need: ments in power, and ergonomics ty and energy		JZU NG	
autonomy thanks to its exclusive patented Easy Moving System with Energy Recovery technology				
Telescopic Column for fu High-end WiFi portable compared to the control of the control				


RADBOOK 2020 103


and profitability








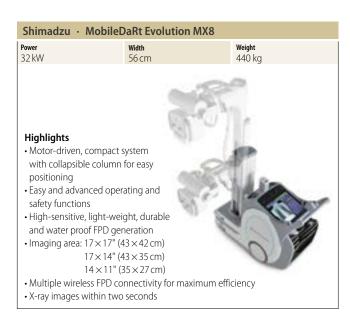

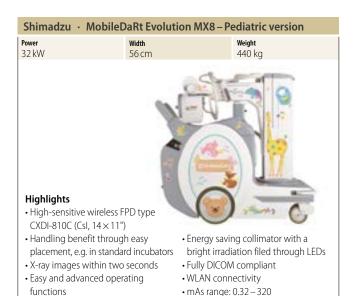

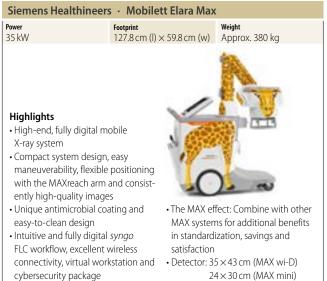



• Battery operating time up to 9 hours

• 10 min charging -

1 hour operation time

· 8 years battery warranty

only for batteries loading

(option)

• Colour LCD touch screen 19"

• Login /identification by badge

• Same interface as Stephanix RAD

rooms, intuitive with unlimited APR

Stephanix · Movix Series DReam Width Weight 20/32/40/50kW 580 kg 67 cm Highlights Compact and light design Motorized up to 5 km/h · Independent from mains, only for batteries loading · Telescopic column and arm, offering wide range of movements for easy positioning • X-ray tube with rotating anode, thin dual focal spots and high heat capacity • Color LCD touch screen 17" • Shareable solution

• kV Range: Up to 150 kVp

• mAs Range: Up to 500 mAs

• Same interface as Stephanix RAD

rooms, intuitive with unlimited APR

Stephanix · Movix DReamy Power 20/32/40/50kW Weight 520 kg Highlights · New ultra-compact and streamlined design • Motorized up to 5.5 km/h • Telescopic column and arm, offering wide range of movements for easy positioning • X-ray tube with rotating anode Up to 150 kV, up to 500 mAs Independent from mains,



• Ready to use within 10 seconds

• Remote diagnostic handle

• Height and reach adjustable drive handle

• Motorized collapsible column support

- Possibility to share wireless detectors with different Stephanix modalities
- Based on capacitive technology for effortless handling

Stephanix · Movix 4/8 DReam Weight 4/8 kW 94 kg Highlights · Lightweight, less than 100 kg • Design for in /outdoor operation • Well-suited for applications at patient bedside, traumatology, paediatrics • Foldable system easy to store and to transport on field · Same interface as Stephanix RAD rooms, intuitive with unlimited APR • Secondary generator control console on monoblock tube head Shareable solution • kV Range: Up to 125 kVp

SternMed · Xenox M	1100	
Power 32 kW	Operation Mains	Motorized No
Highlights Light weight 32 KW rated power suitable for wide clinical applications KV range: from 40 to 125 kV mas range: 0.5 to 200 mAs Exposure time range: 1ms – 1.3s Articulated arm suitable for different		

X-Ray positions • User-friendly workstation with 15" LCD touch screen · Advanced APR function · Advanced flat panel technologyamorphous silicon · Multiple peripheries including USB, · Automatic image processing DICOM and burning CD

Technix · TMB 400 / TMB 400 DR				
Power 40 kW	Width 57.6 cm	Weight 435 kg		
Highlights Battery-motorized unit for easy maneuvering and bedside positioning Freeview technology thanks to telescopic column Battery powered X-ray exposures Two different versions: analogue and digital X-ray Housing Compact design Anatomical programs Telescopic arm 19" touch screen user interface				
Swiveling column Full DICOM connectivity+WLAN				

• Multiple detectors can be interfaced

105 RADBOOK 2020

Integrated generator














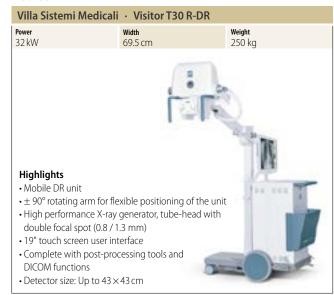



Flatpanel Fluoro










Flatpanel Fluoro









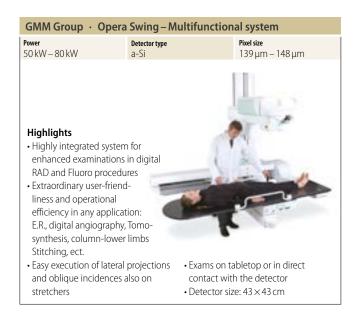



Flatpanel Fluoro

the tube head



- Single end suspended carbon-fibre patient tabletop for total accessibility from any side
- Elevating tabletop with 50 cm minimum distance from the floor
- Full-length patient examination in both vertical and horizontal position
- Full integration with optional ceiling suspension and Wi-Fi detector
- Detector size: $43 \times 43 \,\text{cm}/35 \times 43 \,\text{cm}/24 \times 30 \,\text{cm}$






• Table: Motorized, carbon fiber, floating top with 340° rotation





Shimadzu · FLUOROspeed X1 edition* Pixel size Detector type 80 kW 160 µm Highlights • Premium patient-side R/F system with dynamic flat panel detector • FPD size: 17 × 17" (43 × 43 cm) • 2nd tube for multipurpose room • Real-time image enhancement solution processing technologies • Comprehensive dose management • Superb table operability for easy operation and patient convenience package Bariatric functionality * Product is not available in all countries

Flatpanel Fluoro

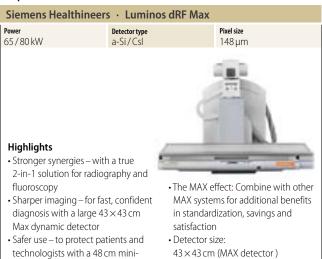

35 × 43 cm (MAX wi-D)

24×30 cm (MAX mini)





• Safer use – to protect patients and technologists

mum table height, full patient ac-

cess from all sides and SmartTouch

- The MAX effect: Combine with other MAX systems for additional benefits in standardization, savings and satisfaction
- Detector size: 43×43 cm (MAX detector), 35×43 cm (MAX wi-D), 24×30 cm (MAX mini)

Ysio Max options: • Fully integrated ceiling-suspended tube with bucky tracking

- MAX wi-D and MAX mini detectors
- · SmartOrtho: long leg and full spine imaging








Flatpanel Fluoro





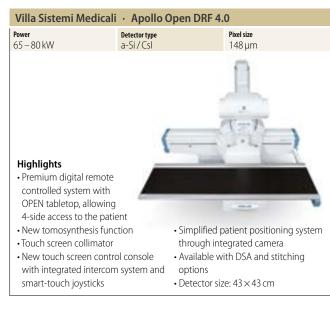



DXA

orthopedic software

60 seconds per site

• Exams can be performed in only

• Powerful easy-to-use software



reference curves.

and lean tissue

· Body composition application for

weight management, tracking fat

DMS Imaging · Stratos DR Detector Scan method System
HF Monoblock Generator 256 elements Pencil beam Highlights • Technology: 2D-Fan Beam • Detector: 256 elements, highest 30 seconds in routine mode image resolution • Compatible with new 3D-DXA technol- Complete solution for an optimal fracture risk diagnosis

- Full options including peadiatric and orthopedic software
- Powerful easy-to-use software
- Exams can be performed in only
- ogy to provide trabecular and cortical bone information with reference curves.
- · Body composition application for weight management, tracking fat and lean tissue

Accessories / Complemantary Systems

Canon Electron Tubes & Devices · XRR-3332 X



Highlights

- 3 inch ROTANODE X-ray tube assembly Adopt large capacity anode target for Mobile systems
- 20 % smaller size / 22 % lighter weight housing than previous model • Size: 1.2/0.6
- XRR-3332X is useful for designing smaller and excellent mobile system.
- to support multipurpose diagnostic application.
- High power input: 46 kW/20 kW (0.1 s) Power: 46 kW/20 kW
 - · Capacity: 300 kHU (Anode heat capacity) 870W (Anode heat dissipation)

Canon Electron Tubes & Devices · XRR-6652 X



- Highlights
- 4 inch ROTANODE X-ray tube assembly for RF systems
- · 20% smaller housing than previous model
- · Can be used as a replacement part for similar models
- Size: 0.8 / 0.3 (Focal Spot) • Power: 52 kW / 12 kW (Input Power)
- Capacity: 600 kHU (Anode heat capacity) 1670 W (Anode heat dissipation)
- High throughput
- (1000W continuous input power)
- Enhanced heat transfer performance
- of heat exchanger
- High resolution image with small focal spot size

Canon Electron Tubes & Devices · XRR-4631G



Highlights

- 4 inch ROTANODE X-ray tube assembly High power input: 100 kW / 40 kW (0.1 s) for DR systems
- 20 % smaller housing than previous model
- · Can be used as a replacement part for similar models
- High cooling rate provided by housing
- · Size: 1.2/0.6 (Focal Spot)
- Power: 100 kW / 40 kW (Max Rating)
- · Capacity: 400 kHU (Anode heat capacity) 1,200W (Anode heat dissipation)

Canon · OrthoMod3D



DRgem · Mobile DR Imaging System for Chest and Chiropractic

Power 32/40/52 kW

 $17 \times 17'' / 17 \times 14''$, wired/wireless



Highlights

- System Concept: Compact Radiography System for Mobile RAD Room
- Mobile imaging radiography system
- Motorized vertical synchronization with wall stand
- Image stitching for whole body
- · Auto numbering function with barcode scanner available
- UPS or Capacitor Generator is recommended

DRgem · GXR series – X-Ray generator

32/40/52/68/82kW

32/40/52kW

32/40kW

Highlights

- · High-frequency generator, perfect for general radiography
- Excellent reproducibility, accuracy, and linearity
- Smaller, lighter modular design
- 1,280 APR conditions with APR utility software
- Tube overloading and housing overheating protection
- Real-time monitoring and self-diagnosis
- Remote diagnosis and automatic calibration
- Adaptive calibration for long-term
- Capacitor Type: compatible with standard wall outlet



- UPS Type: 800 W, free-voltage (100 ~ 240 VAC) line power
- UPS Type: Operation time of up to 12 hours and 3,500 x-ray shots during a power failure



Innovative X-ray technology made in Germany

PROTEC is the synonym for state-of-the-art X-ray technology, designed to offer the highest treatment comfort for patients and medical staff. Sharper, faster and more accurate images enable doctors all over the world to obtain an extensive overview for a clear diagnosis. PROTEC – made in Germany.

The perfect alternative to ceiling mounted systems

The PRS 500 B is a motorized X-ray system with automatic tracking for versatile application and high patient numbers. Thanks to the integrated, intelligent auto tracking control it offers a highly comfortable, fast and efficient daily workflow. Supporting any exposure technique the system provides you with excellent image quality for an extensive and clear diagnosis.

A large touch screen monitor in combination with the CONAXX 2 acquisition software enables users to control the system simple and intuitively. 360° rotation of the complete tube column allows maximum application flexibility, ideal for immobile patients that have to be X-rayed in bedside. Therefore, the PRS 500 B is the perfect cost-saving alternative to ceiling mounted systems.

Symbiosis of manual and motorized

In 35 years of manufacturing X-ray technology we learned that in the daily work flow nothing is more important than performing fast, accurate and patient comfortable X-ray examinations. A drawback of purely motorized systems in aligning as fast as a manual movement can be overcome:

the PRS 500 B combines manual and motorized movements in a perfect symbiosis for the benefit of the user.

While manual movements are supported by precise distance and angle information, the auto tracking system aligns various movements e.g. the SID motorized and automated. As the result users and patients benefit form immense time savings and less examination strain.

Precision, ease of operation and highest treatment comfort – that's what PRS 500 B stands for.

Auto Stitching

The fully digital auto-stitching function of the PRS 500 B is an outstanding solution for spine and full leg examinations. After selecting a stitching examination, the system is automatically positioning itself. Users are guided through all necessary preparation steps directly on the touch screen monitor.

To ensure maximum patient comfort the positioning is carried out in the last step before the actual exposure. With the start of the imaging process, collimator and wall stand automatically position themselves for the relevant partial images. A preview allows you to accept each exposure immediately. After final exposure the complete image is stitched together and ready for diagnosis within seconds.

To prevent projections resulting in blurred stitching images, the PRS 500 B uses a new innovative collimator technology. It automatically fades in all required exposures from a single position and a distance of three meters. This results in razor sharp edges and stitching images for a clear and safe diagnosis.

www.protec-med.com



Accessories / Complementary Systems

Examion · X-Emergency



Highlights

Customized container for digital X-ray.

- U-Arm or Z-Arm design. Z-Arm allows lateral exposures on lying patients
- · Low maintenance effort
- Excellent image quality
- Patient administration
- Mini-PACS or connection to central archives
- Radiological viewer
- Power: 50 kW
- System concept: Wireless or Wired
- Detector size: 14 x 17" / 17 x 17"
- Pixel size: 100 150 μm

I.A.E. · RTC 600



Highlights

- Rotating anode graphite XRay tube, specifically designed for remote controlled table and digital systems
- Enhanced anode heat dissipation, provided by high emittance coating and target design
- Severe tests during conditioning assure reliable performances
- High anode heat storage for repeated loading
- Ground glass window for consistent HVL
- Variety of housings allows flexible systems configurations

I.A.E. · C31-RTM 72

 Size
 0.6/1.2

 Power
 30 kW/75 kW

 Capacity
 300 kHU (Anode 500 W (Anode



Highlights

- Rotating anode X-ray tube unit for mobile x-ray equipment with film and digital detectors
- Lead lined aluminium body
- H.T. cable sockets: type MINI75 4 pin
- Storage and shipment temperature range –10°C/+80° C
- Optional mounting plate for tilting brackets


PTW · Diamentor RS / Diamentor C-RS – Wireless DAP systems



Highlights

- Wireless DAP System
- Integrated DAP chamber and DAP electronics housing (Diamentor RS)
- Automatic air density correction of DAP measurements (Diamentor RS)
- Highly integrated electronics (Diamentor C-RS)
- Available with RS232 or RS485 interface
- Optional Diamentor BT interface for wireless data transfer from DAP Meter to an X-ray unit

PTW · Normi RAD / FLU – X-Ray test object

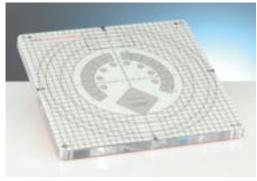

Highlights

- Checks all relevant parameters of analogue and digital fluoroscopic and radiographic X-ray units
- Suitable for routine quality checks on over/under couch tubes and C arms
- Includes an attenuation plate for patient simulation
- Complies with DIN 6868-4 and 6868-150
- Available with the outer format of 300×300 mm or 200×200 mm


PTW · Nomex Multimeter




Highlights

- Dosimetry system, CE marked, class Ilb certified, acc. to IEC 61674
- Single exposure captures all dose values, kVp, time, TF, HVL, frequency, pulses and waveforms
- Angular independent for positioning within the beam
- Fully automatic adjustment
- \bullet Data and waveform export to Excel
- Accessories: Normi RAD/FLU test object

Molecular Imaging

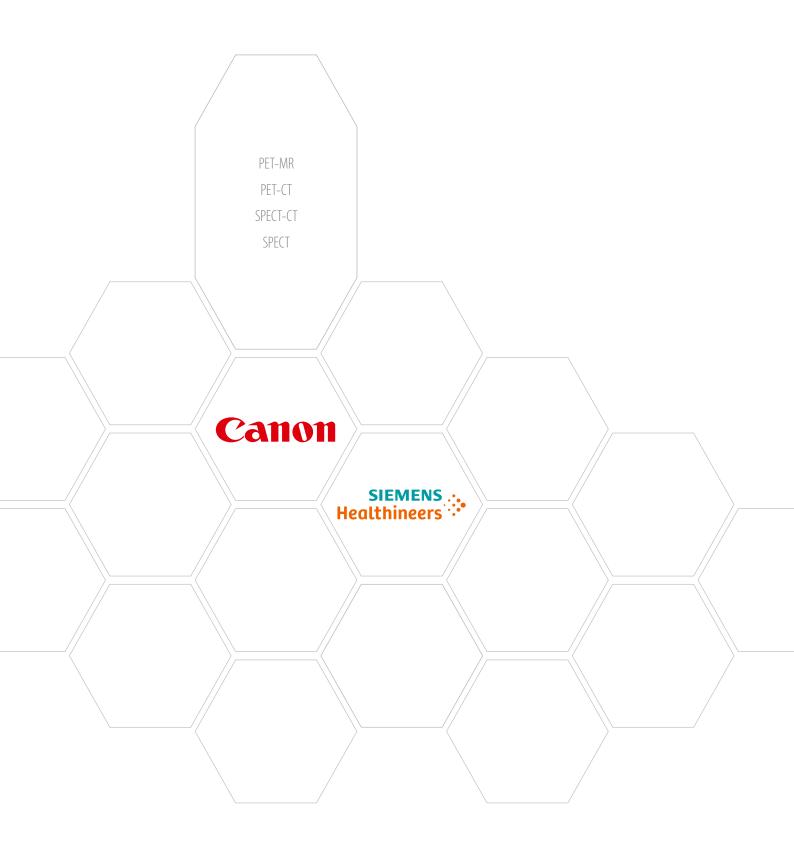

PET-MR

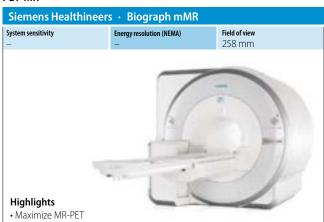


- Benefit from motion-free PET images with MR-based motion compensation beyond gating
- Advance PET attenuation correction with whole-body 5-compartment model including bones and HUGE
- Deliver exceptional quality and speed in MR-PET with the latest MR innovations

PET-CT

· Respiratory gating (option)



· Registration accuracy and repeatability

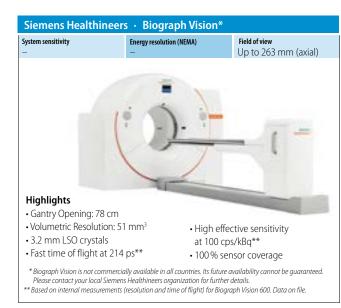





resolution that BGO crystals

• FlowMotion continuous bed motion

• Large 78 cm bore and table capacity of 227 kg (500 lb)

RADBOOK 2020

Please visit us at

www.healthcare-in-europe.com

SPECT-CT




- to distinguish between degenerative disease and cancer
- The first and only system offering accurate and reproducible SPECT quantification
- Up to 68 % lower CT dose* with CARE Dose4D and up to 80 % lower injected dose* with IQ-SPECT to reduce patient radiation risk
- Productivity tools and IQ•SPECT save time and can double patient throughput
- * Based on competitive literature available at time of publication. Data on file.





- Reduction to see more details by reducing metal artifacts. iMAR lets you overcome the effects of metal artifacts in challenging exams
- SAFIRE Sinogram Affirmed Iterative Reconstruction - reduces radiation dose while maintaining image quality
- IVR Interleaved Volume Reconstruction – reconstructs up to 32 slices to evaluate small structures
- · Dual Energy Scan improves image quality with two sequential spiral scans at different energies

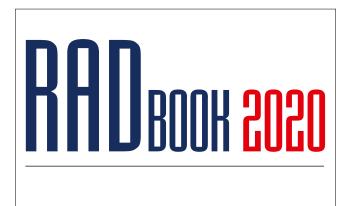


Please visit us at

www.healthcare-in-europe.com

SPECT



- high-capacity, low-height patient bed; and hospital bed imaging capabilities
- Industry-leading* image quality delivers accurate and reproducible clinical information to support physicians' diagnostic confidence
- *Based on competitive literature available at time of publication. Data on file. **Patients up to 227 kg

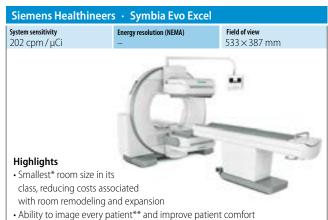

- with a larger bore; a high-capacity, low-height patient bed; and hospital bed imaging capabilities • Industry-leading* image quality delivers accurate and reproducible clinical
- information to support physicians' diagnostic confidence
- * Based on competitive literature available at time of publication. Data on file. ** Patients up to 227 kg

Displays / Printers



Displays - Mammo

JVC • MS-S500 Panel size 21.3" Resolution 2,048×2,560/2,048×7,680 Max. Luminance 3,000 cd/m²



Highlights

- LED Backlight
- 2,000 : 1 contrast ratio
- True 11 Bit grayscale
- ISD Support
- Front and ambient light sensor
- Remote management and calibration
- Integrated power supply
- DVI and DisplayPort interface

JVC · CL-S500		
Panel size 21.3"	Resolution 2,048 × 2,560	Max. Luminance 1,150 cd/m ²
Highlights • Panel technology: IPS • 2,000: 1 contrast ratio • Auto Text Mode • Dynamic Gamma • Front and ambient light • Remote management ar • Integrated power supply • DVI and DisplayPort inter	nd calibration	

Displays - Color

• Optional AR coating

• Auto Text mode and Dynamic Gamma








Panel size 21.3" Resolution 2,048×1,536 Highlights Panel technology: IPS 1,500:1 contrast ratio Front and ambient light sensor Remote management and calibration Integrated power supply DVI and DisplayPort interface Optional AR coating Auto Text mode and Dynamic Gamma

Displays - Color

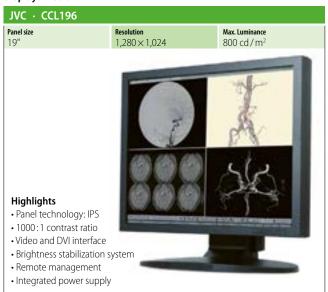

Displays - Grayscale











DVD Import

Nexus/Chili · Import Robot



Highlights

- Automatic import robot
- Import of patient CD/DVD
- 2, 5 or 10 drives
- 2 import trays (regular / express)
- 2 output trays (ok, failed)
- · Optional virus scan
- Correction of foreign data
- Automatic DICOM transfer
- · Works with any PACS

Printers

Nexus/Chili · Burn Gateway



Highlights

- Receives data by DICOM C-Store
- Burns data on one or more CD/DVDs
- Optional reports
- Individual label printing
- Client enabled (different logos)
- · CHILI viewer in report quality
- Alternative presentation as HTML and JPEG
- · Certified by OFFIS and DRG
- · Works with any PACS
- · External output tray

DVD Burner

medigration · CD-Imager



Hiahliahts

- Fully automatic compact system for creating DICOM patient CDs or DVDs
- Highly compatible with all digital DICOM modalities (multimodality)
- Individual labeling (practice / clinic logo)
- Easy integration of DICOM patient data
- Extremely cost effective due to quick printing times and low link consumption
- Format: CD-R, DVD-R, DVD+R, DVD-R DL, DVD+R DL
- · Capacity: 30 CDs/h or 15 DVDs/h (burn and print)
- Magazine size: 2 × 50 pcs



- Multi-application hardcopy solution, including digital mammography
- Integrated A#Sharp technology for optimized image quality
- Two multi-format trays, each supporting different film sizes and types
- Very short access time for extremely fast delivery of first four prints

Agfa HealthCare · Drystar 5503

Resolution Technology **Capacity** 100 films/h (14×17) Direct digital imaging 508 dpi

Highlights

- · Multi-modality, high throughput imager with film sorter
- Ideal for centralized workflow, can easily be connected to the network
- Integrated A#Sharp technology for optimized image quality
- Three multi-format trays, each supporting different film sizes and types
- Suitable for CT, MRI, DSA, digital R/F, CR, DR and optional mammography applications





Printers



· Very short access time ensures fast printing of small print jobs














Frequency range 1 –33 MHz	Display mode 2D/3D/4D	Display size 23"
Imaging, ApliPure ADF, Intelligent Su Imaging, Ultra Hig • CEUS; Advanced C MicroFlow imagin • 4D (surface, MPR, Shadow Glass) • FlyThru virtual end Fusion, Strain and wave elastograph	EUS incl. VRI, g, Quad View MultiView, Luminance, doscopy, Smart Quad View Shear- y, Dispersion croPure, Auto IMT,	ion



IMT, AUTO NT, Wall Motion Tracking



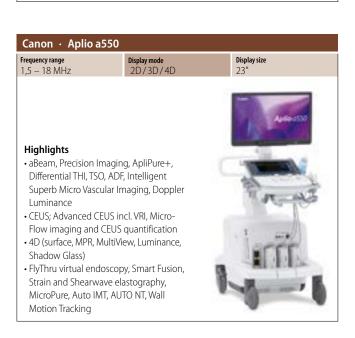



Frequency range Display mode Display size 1.8 - 18 MHz 2D/3D/4D 21.5"



Highlights

- aBeam, Precision Imaging, ApliPure+, Differential THI, TSO, ADF, Superb Micro Vascular Imaging, Doppler Luiminance
- CEUS; Advanced CEUS incl. VRI, Micro-Flow imaging
- 4D (surface, MPR, MultiView, Luminance)
- Smart Fusion, Strain and Shearwave elastography, MicroPure, Auto IMT, AUTO NT, Wall Motion Tracking

· High Density Beamformer, Precision Imag-

ing, ApliPure+, Differential THI, TSO, ADF,

• CEUS; Advanced CEUS incl. VRI, MicroFlow

• 4D (surface, MPR, MultiView, Luminance)

• Smart Fusion, Strain and Shearwave elasto-

graphy, MicroPure, Auto IMT, AUTO NT, Wall

Superb Micro Vascular Imaging

imaging and CEUS quantification

Canon · Aplio 400		
Frequency range 1.8 – 18 MHz	Display mode 2D/3D/4D	Display size 21.5"

Highlights

- High Density Beamformer, Precision Imaging, ApliPure+, Differential THI, TSO, ADF, Superb Micro Vascular Imaging
- CEUS; Advanced CEUS incl. VRI, MicroFlow imaging and CEUS quantification
- 4D (surface, MPR, MultiView, Luminance)
- Smart Fusion, Strain and Shearwave elastography, MicroPure, Auto IMT, AUTO NT, Wall Motion Tracking

Canon · Xario 200 Platinum Edition



Canon · Aplio 300

Highlights

 Frequency range
 Display mode
 Display size

 1.8 – 18 MHz
 2D/3D/4D
 21.5"



Highlights

Frequency range

1.8-18 MHz

- High Density Beamformer, Precision Imaging, ApliPure+, Differential THI, Tissue
- Enhancement, Advanced Dynamic Flow, Superb Microvascular Imaging (SMI)
- 4D-imaging; surface rendering, MPR, MultiView, Freehand 3D, Luminance
- Realtime Shearwave elasto with propagation curves, strain elasto, Auto IMT, Stress Echo, 2D Wall Motion Tracking, CEUS contrast imaging
- iStyle+ productivity suite with fully customizable panel, agile housing, height adjustable console, panel swivel, Quick Start, Quick Scan & Quick Assist



Display size

Canon · Xario 200G

Motion Tracking

 Frequency range
 Display mode
 Display size

 1.8 – 18 MHz
 2D/3D/4D
 21.5"

Highlights

- Up to eight hours battery autonomy, two seconds startup from standby, 21.5 inch wide screen display
- High Density Beamformer, Precision Imaging, ApliPure+, Differential THI, Tissue Enhancement, Advanced Dynamic Flow, Superb Microvascular Imaging (SMI)
- ${\color{red} \bullet} \, 4D\text{-imaging; SR, MPR, MultiView, Freehand 3D, Luminance } \\$
- Shearwave elastography with propagation map, strain elasto, Auto IMT, Stress Echo, 2D Wall Motion Tracking, CEUS contrast imaging
- iStyle+ productivity suite with fully customizable panel, agile housing, height adjustable console, panel swivel, Quick Start, Quick Scan & Quick Assist, Extensive line-up of transducers



Canon · Xario 100 Platinum Edition

 $\begin{array}{cccc} \textbf{Frequency range} & \textbf{Display mode} & \textbf{Display mode} \\ 1.8-18 \ \text{MHz} & 2 \ \text{D/} \ 3 \ \text{D/} \ 4 \ \text{D} & 19 \ \text{T} \end{array}$

Highlights

- High Density Beamformer, Precision Imaging, ApliPure+, Differential THI, Tissue Enhancement, Advanced Dynamic Flow
- 4D-imaging; surface rendering, MPR, MultiView, Freehand 3D, Luminance
- Realtime elastography, Auto IMT, Stress Echo, EUS contrast imaging
- iStyle+ productivity suite with fully customizable panel, agile housing, height selectable console, Quick Start, Quick Scan & Quick Assist



Canon · Xario 100G Display mode 2D/3D/4D Display size 21.5" Frequency range 1.8 – 18 MHz

Highlights

- Up to 4 hours battery autonomy, two seconds startup from standby, 21.5 inch wide screen display
- · High Density Beamformer, Precision Imaging, Apli-Pure+, Differential THI, Tissue Enhancement, Advanced Dynamic Flow
- 4D-imaging; surface rendering, MPR, MultiView, Freehand 3D, Luminance
- Realtime elastography, Auto IMT, Stress Echo,
- iStyle+ productivity suite with fully customizable panel, agile housing, height selectable console, Quick Start, Quick Scan & Quick Assist, Extensive line-up of transducers



Canon · Viamo c100		
Frequency range 1.5 – 18 MHz	Display mode 2D/3D/4D	Display size
Highlights		
Up to two hours15 inch screen dieHigh Density Bea	, ,	

- 3D-imaging: surface rendering, Virtual HD, depth cueing
- Auto IMT, needle enhancement, Quick Scan
- Magnesium alloy housing, height adjustable on cart, Quick Start (eight seconds from standby).



Canon · Viamo sv7 Frequency range Display size 1.5-12 MHz Highlights • Portable ultrasound system • Multi touch screen, Tablet mode possible • Single transducer input, expandable to three transducers • Battery and AC operation, very fast boot time

Highlights

Frequency range

1-22 MHz

· Multi-disciplinary Premium platform, ergonomic design

Hitachi · Arietta 850SE

- Pure Image Symphonic Architecture
- 23" LCD monitor for highest contrast
- Wide range of transducers for GI, interventional guidance, urology and TEE applications
- Advanced modalities: SWM, ATT, Real-time Elastography, CEUS, RVS Fusion, Needle and Body Motion tracking, 3D/4D
- Advanced analysis: TIC, eTracking, WI, 2DTT, Protocol assistant, Auto Measurements

Protocol assistant, Auto Measurements



Display size

Hitachi · Arietta 850		
Frequency range 1 – 22 MHz	Display mode 3D/4D	Display size 22"

(< 10 s from standby to scanning), three hours battery support

• Highly programmable Touch Screen, few buttons, easy to operate,

protocol assist. Precision Plus, Aplipure+, Cardiac measurement kit and

• High color sensitivity, exceptional image quality

OB measurement package kit are pre-installed.

Highlights

- Multi-disciplinary Premium platform, ergonomic design
- Pure Image Symphonic Architecture
- 22" OLED monitor for highest contrast
- Wide range of transducers for GI, interventional guidance, urology and TEE applications
- Advanced modalities: SWM, ATT, Real-time Elastography, Combi-Elasto, CEUS, RVS Fusion, 3D SIM navigator, E-field Simulator, Needle and Body Motion tracking, 3D/4D
- · Advanced analysis: TIC, eTracking, WI, 2DTT, Protocol assistant, Auto Measurements



Hitachi · Arietta 750				
Frequency range 1 — 18 MHz	Display mode 3D/4D		Display size 23"	
Highlights • Multi-disciplinary Premite ergonomic design • Pure Image Symphonic and a substitution of the properties of transduced guidance, urology and Total and a substitution of the properties of the prope	Architecture contrast ors for GI, intervention EE applications /M, ATT, Real-time oi-Elasto, RVS Fusion, atracking, 3D/4D	nal		

Frequency range Display mode Display size 1 – 18 MHz 3D/4D 21.5"

Highlights


- Compact multi-disciplinary platform with comfortable workflow, high definition imaging and useful application from premium platform
- Automated process features: Protocol Assistant, Auto-Optimizer, Auto Measurement
- Battery
- Unique image processing technology underpin outstanding image quality
- 21.5" LCD Widescreen Monitor, fully articulating monitor arm
- Wide range of transducers for GI, urology and TEE applications
- Advanced modalities & analysis: Dual gate Doppler, Strain Elastography, CEUS, 2DTT



Hitachi · Arietta 50 Frequency range 1 – 15 MHz 2D Display mode 21" Highlights · Lightweight compact multi-disciplinary platform with intuitive workflow, clear imaging, and easy-to-use applications · Automated process features: Auto Optimizer, Doppler auto trace, Automated Measurement (IMT, NT) · Simplified and user-friendly interface · Battery · Unique image processing technology underpin

Hitachi · Arietta 50L	E	
Frequency range 1 – 15 MHz	Display mode 2D	Display size 21"
Highlights Lightweight compact muplatform with intuitive wimaging, and easy-to-use Automated process featu Doppler auto trace, Automent (IMT, NT) Simplified and user-friend Battery Unique image processing outstanding image quali Widescreen 21" LCD mor	ulti-disciplinary orkflow, clear e applications ures: Auto Optimizer, mated Measure- dly interface g technology underpin ty nitor with fully	
articulating monitor armTransducers for GIAdvanced modalities: CE		



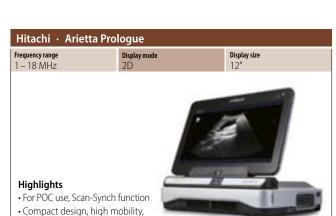

• Simple and intuitive to use, tablet-style with touch screen control

• Option of 14 transducers, offers high quality imaging for a broad range

of applications including MSK, rheumatology, emergency medicine,

in-built battery for portable use

anaesthesiology

• Hand carry, can be used with probe tray or cart

• Ethernet, Wi-Fi, Bluetooth network connections

Konica Minolta · Sonimage HS1

outstanding image quality

• Widescreen 21" LCD monitor

· Advanced modalities: CEUS

• Transducers for GI

requency range	Display mode	Display size
Jp to 18 MHz	2D	15"

Highlights

- Triad Tissue Harmonic Imaging (3THI)
- $\bullet\,\mathsf{SNV}-\mathsf{Simple}\,\mathsf{Needle}\,\mathsf{Visualization}$
- Newly developed multi-frequency probes up to 18 Mhz
- Portable system with built-in battery
- Start-up from standby within 15 seconds
- Excellent for MSK/orthopaedic, nerve, vascular and anaesthesia
- Rotatable and tiltable 15 inch touchscreen
- Weight: 7.8 kg



Konica Minolta · Sonimage MX1 Frequency range Up to 14 MHz Display mode 2D Display size 12.1"



Highlights

- Portable Point-of-Care ultrasound
- Dual Sonic Technology
- Tailored for MSK & Anesthesia
- One-touch image optimization for quick diagnosis
- Multi Parameter Adjuster function
- Simple Needle Visualization (SNV) for guided injections
- High quality Made in Japan
- Weight: 4,5 kg

Mindray Medical · Resona 7		
Frequency range 1 – 20 MHz	Display mode 3D/4D	Display size 21.5"
		'氖!
 Highlights Powered by ZST+ platt generation ZONE Son- nology based on Char 	ography Tech-	
Software processing. • A premium ultrasound helps customers to se • Faster and more accur • Complete set of tools Imaging, RIMT and CE	e more. ate images. –V Flow, Fusion	d

Mindray Medical · Resona 6 Frequency range 1 – 20 MHz Display mode 3D/4D Display size 21.5"

Highlights

- Powered by ZST+ platform, the next generation ZONE Sonography Technology based on Channel Domain Software processing.
- A premium ultrasound system that helps customers to see more.
- Faster and more accurate images.
- Complete functionality for Radiology and clinical research



Mindray Medical · DC-80 with X-Insight

quency range – 20 MHz	Display mode 3D/4D	Display size 21.5" / 23.8"

Highlights

clinical research

- Single Crystal with 3T technology transducer, ComboWave transducer, interventional transducer
- Best in class shear wave(STE & STQ), NTE (shell), UWN+ CEUS, ART Flow, TT QA, LVO
- Dual-wing floating arm, powerful and intuitive gestures, MedTouch



Mindray Medical · DC-80A with X-Insight

 Frequency range
 Display mode
 Display size

 1 - 20 MHz
 3D/4D
 23.8"



Highlights

- Superb 3D/4D with single crystal volume and Hyaline
- Outstanding ABD image in both penetration and resolution
- Most intelligent Smart Planes CNS and Smart Face
- Large touch screen (13.3") & Full HD monitor (23.8"), five active sockets
- Best in class shear wave (STE & STQ)
- Built-in battery for continuous scanning

Frequency range 1 – 20 MHz Display mode 3D / 4D Display size 13.3" / 21.5" / 23.8" Highlights Top in class 3D/4D with single crystal volume and Hyaline Best in class ABD image in both penetration and resolution Most intelligent Smart Planes CNS and Smart Face Largest Full HD monitor (21.5" / 23.8") and ultra-slim touch screen (13.3")








Visit us at ECR2020, July 15-19, Expo X5/543




Mindray Medical · MX7				
Frequency range 1 – 20 MHz	Display mode 3D/4D	Display size 15.6"		
Highlights				
• 15.6" IPS monitor,				
12.3" IPS touch screen		lastography imaging		
 Cutting-edge ZST⁺ platf Eight hours continuous 		tress echo DI and QA		
Magnetic power socket		VO		
Contrast imaging	• il	Needle+		





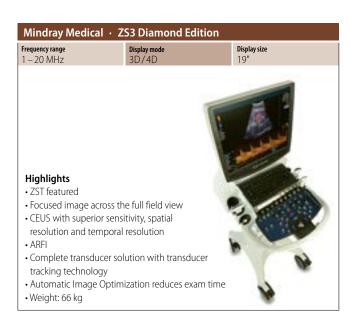

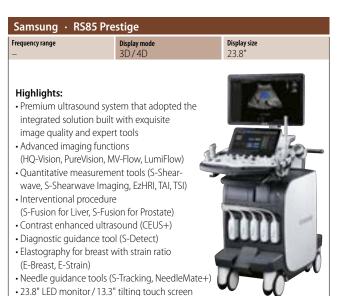

Frequency range Display mode Display size 3D/4D 23"

Highlights:

- Premium ultrasound system that offers superior imaging performance for radiology
- Interventional procedure (S-Fusion)
- Contrast enhanced ultrasound (CEUS+)
- Quantitative measurement tools (S-Shearwave, S-3D Arterial Analysis)
- Advanced imaging function (HQ Vision)
- Elastography for breast with strain ratio (E-Breast, E-Strain)
- Diagnostic guidance tool (S-Detect)
- Needle guidance tools (S-Tracking, NeedleMate+)
- 23" LED monitor / 13.3" tilting touch screen

(3 transducer ports/extended battery)

User-oriented features

(Quick Preset, EZ-Exam+, QuickScan)
• 21.5" LED monitor / 10.1" touch screen



Samsung · HS70A with Prime		
requency range –	Display mode 3D/4D/5D	Display size 23"
Highlights: - High-end ultrasoun high efficiency in pr - Diagnostic guidance - Contrast enhanced - Quantitative measu (S-Shearwave, Arter - Elastography for bre	actice e tool (S-Detect) ultrasound (CEUS+) rement tools ial Analysis)	

(E-Breast, E-Strain)

(NeedleMate+)

• Cardiac solutions (Strain+, Stress Echo)

• 23" LED monitor / 10.1" touch screen

• Needle guidance technology

Samsung · HM70A with Plus			
Frequency range —	Display mode 3D/4D	Display size	
Highlights: Laptop design to suit varenvironments Advanced imaging funct (ClearVision, HDVI, S-Flowns) Simple and accurate intime assurement (Auto IMT) Highly sensitive elastogrous Needle guidance technology User-oriented features (ADVR, EZ-Exam+, Quicknowns) Fast booting within 20 seconds of the second of the seconds of the second of the seconds of the second of the seconds of the second of the seconds of the seconds of the seconds of the seconds of t	rious diagnostic tions v) ma-media thickness) aphy (ElastoScan) ology (NeedleMate) Scan)		

_	Display mode 3D/4D/5D	Display size 21,5"
Highlights: • High-end ultrasour with versatile diagr • Advanced imaging (S-Harmonic, Clear • Simple and accurar measurement (Aut	nostic solutions of functions Vision, S-Flow) te intima-media thickness	

Samsung · HS5 Frequency range	Display mode	Display size
-	3D/4D/5D	21.5"
measurement (Aut	ns range g functions Vision, S-Flow) te intima-media thickness	

Frequency range —	Display mode 3D/4D/5D	Display size 21.5"
Highlights: Fully equipped ultrasc for everyday efficiency Advanced imaging fu (S-Harmonic, ClearVisi Simple and accurate i measurement (Auto II Cardiac solution (Strai Needle guidance tech User-oriented feature: (Quick Preset, EZ-Exan 21.5" LED monitor / 10 Gas lift / Side pocket	ound system nctions on, MultiVision) ntima-media thickness WT+) n+) nology (NeedleMate+) s n+, QuickScan)	



Siemens Healthineers · Acuson Sequoia Ultrasound System			
Frequency range 1 – 17.8 MHz	Display mode 2D		olay size .6" / 22"
Highlights • Unique BioAcoustic highest signal fideli penetration withou quality of up to 40 c • Ultra-fast, fully focu without degradatio	technology delivers ty for deep tissue t impacting diagnostic	300	
without flashing or	,	m //	-
Advanced applications, such as contrast-enhanced ultrasound and elastography allow a personalized			
approach to provid	e accurate and reprodu	cible results	





Siemens Healthineers · Acuson S Family Ultrasound Systems			
Frequency range 1 — 18 MHz	Display mode 2D/3D/4D	Display size 12.1 / 21.5"	
Highlights • Superior imaging perfori imaging, interventional i women's health with next HD transducer technolog. • Advanced applications to capabilities: eSieFusion in imaging, ARFI shear wav tography, contrast-enhal. • Intuitive, user-centric wowith simplified control punnecessary keystrokes	radiology and ext generation gy o expand clinical multi-modality ee & manual elasnced ultrasound orkflow design		

Siemens Healthineers · Acuson Juniper Ultrasound System				
Frequency range 1.1 — 18 MHz	Display mode 2D/3D/4D	Display size 13.3 / 21.5"		
Highlights - High-performance, share virtually every patient wis smallest footprint - Five active transducer posupport 16 transducers for capabilities – from radiol radiology, cardiology, urgand OB/GYN - High-fidelity acoustic signoise and offer premium	th one of the industry's orts and one CW port or a wide variety of ogy, interventional ology to orthopedics nals greatly reduce			

132 RADBOOK 2020

industry-leading elasticity solutions

Frequency range Display mode Display size 1.3 – 16 MHz 2D/3D/4D 10.4"/21.5"



HighlightsPowerful platform driven by efficiency and built for performance.

- Intuitive user interface with up to 28% fewer keystrokes and 3x more user-defined keys
 21.5" HD display and 220° endo-cavity
- transducer provides expanded field of view
- 10.4 inch touch display with swipe motion
- Transducer compatibility with existing and legacy Siemens Healthineers systems

Siemens Healthineers · Acuson NX3			
Frequency range 1.3 – 12 MHz	Display mode 2D/3D/4D	Display size 10.4"/21.5"	

Highlights

- Powerful platform driven by efficiency and built for performance
- Intuitive user interface with up to 28 % fewer keystrokes and 3 x more user-defined keys
- 21.5" HD display provides expanded field of view
- 10.4 inch touch display with swipe motion
- Transducer compatibility with existing and legacy Siemens Healthineers systems

Siemens Healthineers · Acuson NX2



Siemens Healthineers · Acuson NX2 Elite

 Frequency range
 Display mode
 Display size

 2 - 10 MHz
 2D
 21.5"



- Provides premium imaging performance using a cost-efficient, ten-transducer set to perform a wide range of exam types at a sustainable value
- Intuitive control panel design combined with up to four front-facing transducer ports optimize workflow efficiency
- Large 21.5" 1,080 p HD display; Twice the pixel density
- Migrated optional advanced clinical applications such as DTI, eSie Touch elasticity & advanced foursight technology



Frequency range

2 – 10 MHz

Frequency range 1 – 18 MHz

- Highlights
 Provides premium imaging performance using a cost-efficient, eight-transducer set to perform a wide range of exam types at a sustainable value
- Intuitive control panel design combined with up to four front-facing transducer ports optimize workflow efficiency
- Large 21.5" 1,080 p HD display; Twice the pixel density
- Simplified control panel designed to enable operator efficiency and speed-up completion of essential tasks



Display size

Display size

Siemens Healthineers · Acuson SC2000 PRIME

 Frequency range
 Display mode
 Display size

 1.25 – 10 MHz
 2D/3D/4D
 21.5"

Highlights

- One-click automated aortic and mitral valve modeling and measurements within seconds with eSie Valves
- Capable of 2D and 4D TTE, TEE, and ICE and TrueFusion on one system
- Comprehensive 1-click knowledge-based advanced clinical applications: eSie Measure, eSie LVA, eSie Left Heart and eSie Valves
- Advanced applications to support routine echo and interventional guidance with eSie PISA, eSie VVI, Volume Right Ventricular Analysis (RVA), Septal Guide, TrueFusion and more



Siemens Healthineers · Acuson S2000 Breast Volume Scanner

Display mode 2D/3D/4D

Highlights

- Automated volume acquisition for operator-independent, standardized 3D imaging to enable consistent, reproducible

• Excellent 2D imaging capabilities using hand-held high-frequency HD transducers

breast imaging

results improving the quality of

- Advanced technologies to expand clinical capabilities: manual and shear wave elastography, multi-modality review
- Read anytime. Anywhere. syngo.Ultrasound Breast Analysis reading software

Siemens Healthineers · Acuson Freestyle Elite Frequency range Display mode 2 – 15 MHz



Highlights

- With cable-free technology to offer unrestricted access to practitioners at the point of care, allowing quicker turnaround time
- Enhanced needle visualization and Pixelformer image processing architecture on an expanded image display may improve procedural confidence in interventional settings
- Automatically populate patient registration data between systems with Artis Patient Synchronization using Artis Access

Sigmons Hoalthingars Agusan Franctula			
Siemens Healthineers · Acuson Freestyle			
Frequency range 2 – 15 MHz	Display mode 2D	Display size 15"	

Highlights

- With cable-free technology to offer unrestricted access to practitioners at the point of care, allowing quicker turnaround time
- Enhanced needle visualization and Pixelformer image processing architecture on an expanded image display improve procedural confidence in interventional settings
- Empowered workflow with zero cable-drag and single-user operation via integrated scanning controls

KHUBOOK 2020

Please visit us at


www.healthcare-in-europe.com

SIUI · Apogee 6500

ransducer inputs	Display mode	Display size
-	3D/4D	23"/13.3"

Highlights

- · Radiology: SEQ (Shearwave Elastography Quantification) with convex probe is a quantitative assessment tool for easy measurement of tissue hardness; CHI empowers the detection and qualitative diagnosis of tumors in liver, kidney and other organs
- · Cardiology: Auto EF, Auto SG, AMM (Anatomic M-mode)
- OB/GYN: Auto OB supports auto measurements of HC, C, FL, BPD; Lumi 4D is an image rendering processing technology to present more lifelike 4D image with higher resolution

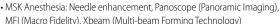

MFI (Macro Fidelity), Xbeam (Multi-beam Forming Technology) Urology & Anorectal: Convenient diagnosis without discomfort with convex probe, linear probe, transrectal linear probe and transrectal bi-plane probe

SIUI · Apogee 6300

Display size Transducer inputs

Highlights

- Radiology: Contrast Imaging empowers the detection and qualitative diagnosis of tumors in liver, kidney and other organs
- · Cardiology: Auto EF automatically measures cardiac volumes based on speckle tracking technology; Auto SG is the automated quantitative assessment of both global and regional myocardial functions; TDI helps to assess the myocardium movement ability, by TDI-PW to calculate the velocity and the ration E/E', A/A'
- OB/GYN: Auto OB supports auto measurements of HC, AC, FL, BPD; Lumi 4D
- SIUI MAI: Available for remote pre-sales demonstration, remote application and technical instructions, online public classes etc.



SIUI · Apogee 5800 Genius

Display size

Highlights

- SEQ (Shearwave Elastography Quantification) offers a much more accurate way to detect disease of some deeper organs
- Ultrasound Contrast Imaging applies in clinical difficult parts such as 2D liver, 4D oviduct contract imaging etc.
- Professional functions like Auto EF, Auto Strain, Stress Echo offers ideal clinical solutions for cardiology
- · SIUI MAI: Available for remote pre-sales demonstration, remote application and technical instructions, online public classes etc.









Ultrasound			
SIUI · Apogee 5300V Neo			
Transducer inputs 4	Display mode 3D/4D	Display size 10.4"	
Advanced technology FusionFreq, Fusion Innovative diagnology Flow, Smarchive Complete application Systems, cardiology Powerful data mand hard disk, DICOM	ostic tools: ECG,TDI, CW, VS tion: Abdomen, reproducti	ve .	

SIUI · Apogee 3500 Elite			
Transducer inputs 4	Display mode 3D/4D	Display size 18.5" / 8.4"	
Highlights • 18.5" LED monitor / 8.4" independent probe and • OB / GYN: The system su transvaginal volumetric pimaging in OBG, and line graphy image • Cardiology: The system is complete functions and cardiac measurement pa 800 items • Radiology: The system b for radiology including V Needle enhancement et	cable management pports both abdominal and probes to realize 3D/4D ear probes with Elastose endowed with professional ackages with over poasts all functions S Flow, Panoscope,		











ment package, DICOM 3.0

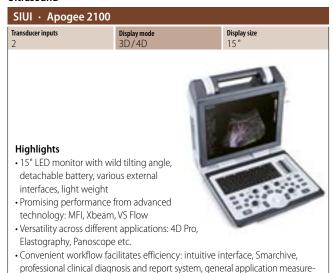





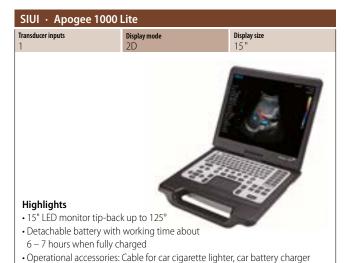


• A variety of professional probes cover all applications including OBG,

Cardiology, Urology, abdomen, small parts etc.

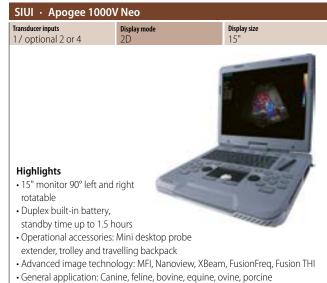

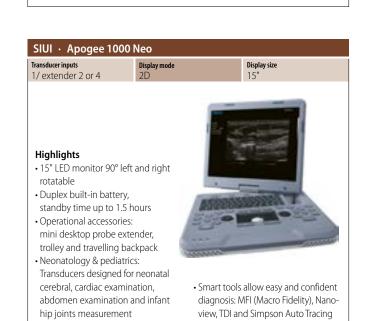



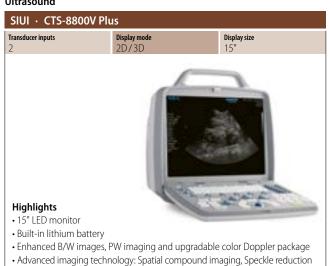

Transducer inputs 1	Display mode 2D	Display size 7"	
Highlights • Hand held ultrasoun: • 7" WVGA LCD monitc • Environmental rating and IP67 (probe head • Three hours operatin • Software and report system • Gravity sensor for lay (trans- verse/vertical) • Measurement for dis	g: IP54 (main unit) d) g time for reproductive out change	als	

SIUI · IBUS 60			
Transducer inputs 4	Display mode 3D	Display size 21.3"/15"/10.4"	
Highlights Intelligent breast volume system, third of its kind ireduces missed diagnosi breast and dense breast Acquires data from multimulti-angle and multide coronal scanning intuitive anatomical information Exceptional diagnosis ex 3D-positioning, Nanopu system, mute scanning	n the world age and s, ideal for exams -section, rection ely shows f breast tissue perience:	uction technology, breast health report	



Accessories / Complementary Systems

technology, Trapezoidal imaging etc.

· General application: Canine, feline, bovine, equine

GCTechnology · CIRS Phantoms





Fetal ultrasound phantom family

Multimodality abdominal phantom

Highlights

phantoms

- Doppler flow phantom
- Quality assurance test phantoms
- Ultrasound accreditation phantoms
- Male and female ultrasound pelvic
- Thyroid ultrasound training phantom Elastography phantoms
- Prostate phantom family -Breast phantom family
- Kidney training phantom
- · Vascular access training phantom kit
- Shear wave liver fibrosis phantoms




Please visit us at

www.healthcare-in-europe.com



IBA Dosimetry · DIGI-13

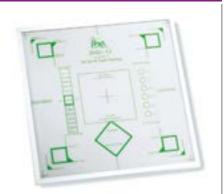

Highlights

For quality checks at digital radiographic systems (CR/DR) according DIN 6868-13.

Test parameter:

- Uniformity
- Spatial resolution
- · Alignment of light and beam field
- Image scale
- Artifacts
- Geometrical distortion

IBA Dosimetry · ETR1 incl. Centering Tube



Highlights

For quality checks of conventional radiography systems; according DIN 6868-3; including holder for chest wall stand.

Test parameter:
• Spatial resolution

Low contrast

- Alignment of light and beam field
- Geometrical distortion
- · Measuring areas for optical density

IBA Dosimetry · Test Device Primus A



Highlights

Test device Primus A is designed according DIN 6868-150 & DIN 6868-4 for Quality assurance at radiography and fluoroscopy systems.

- 17 steps for dynamic verification
- 8 low contrast sensitivity circles
- \bullet Grid for easy and efficient determination of light- & beam field alignment as well as geometrical distortions

IBA Dosimetry · Mammo-14



Hignlights

For quality assurance / constancy test at digital mammography systems according DIN 6868-14.

- 40 mm base plate with integrated Al step wedge and 2 rows of steel balls, for checking the image limitation towards the thorax side.
- 6 mm structural plate with recess for test inserts
- Test insert: PMMA, SDNR & High Contrast
- $3 \times 20 \text{ mm} / 1 \times 10 \text{ mm} / 1 \times 4 \text{ mm}$ PMMA attenuation plates
- 2×20 mm PMMA full field attenuation plate (260×320 mm)

IBA Dosimetry · Mammo-152



Highlights

For quality assurance / acceptance and constancy tests according DIN 6868-152, DIN 6868-7, IEC 61223-3-2 and EPQC (EUREF) in conventional mammography.

Test parameter:

- Object thickness and tube voltage compensation resp. AEC reproducibility
- Spatial and contrast resolution
- Artifacts / Geometry
- Check of the image limitation towards the thorax side

IBA Dosimetry · Mammo-162



Higninginus

For quality assurance / acceptance test of digital Mammography Systems, according DIN 6868-162.

- 40 mm base plate with integrated AI step wedge and 2 rows of steel balls, for checking the image limitation towards the thorax side.
- 6 mm structural plate with recess for test inserts
- Test insert: PMMA, SDNR & High Contrast
- $\cdot 3 \times 20 \text{ mm} / 1 \times 10 \text{ mm} / 1 \times 4 \text{ mm}$ PMMA attenuation plates
- 1 \times 20 mm PMMA full field attenuation plate (260 \times 320 mm)

IBA Dosimetry · 2-part PMMA CT-Phantom



Highlights

Phantom for measurements of CTDI according IEC 60601-2-44, IEC 61223-3-5,

- 1 Adult Head-Phantom, 16 cm diameter, 5 holes
- 1 Adult Body anulus, 32 cm diameter, 4 holes
- 9 Acrylic rods for plugging in all phantom holes

IBA Dosimetry · 3-part PMMA CT-Phantom



Highlights

Phantom for CTDI measurements, according IEC 60601-2-44, IEC 61223-3-5, IEC 61223-2-6.

- Innovative 3-part nested phantom according FDA 21 CFR 1020.33.
- 1 Pediatric Phantom, 10 cm diameter, 5 holes
- 1 Adult Head anulus, 16 cm diameter,
- 4 holes
- 1 Adult Body anulus, 32 cm diameter, 4 holes
- 13 Acrylic rods for plugging in all phantom holes

IBA Dosimetry · DSA Test Device



Highlights

For Quality Assurance of "Digital Subtraction Angiography" (according DIN 6868-150, DIN 6868-4, IEC 61223-3-3)

Test parameter:

- Cupper dynamic step wedge with logarithmic check
- DSA contrast sensitivity
- Artefacts

IBA Dosimetry · Multimeter MagicMaX Universal



- XM Mammography
- DCT10-MM Ionization Chamber for CT

Measurement parameter:

Dose/dose rate - dose per pulse - kVp/PPV -time -total filtration - HVL wave form - dose, dose rate length product for CT

IBA Dosimetry · Dosimax plus A (HV)*

*Dosimax plus A HV with integrated high voltage for meas-urements at CTs with ionization chamber DCT10-RS



Highlights

Single channel dose meter according IEC 61674 for acceptance tests at Radiography-, Fluoroscopy-, Dentaland Mammography systems. Available with RQA/RQM/DCT10-RS* • Time: 1 ms - 9,999 s

Measurement parameter (RQA):

- Dose: 200 nGy 9,999 mGy
- Dose rate: 80 nGy/s 70 mGy/s

IBA Dosimetry · Dosimax plus I



Highlights

Single channel dose meter according IEC 61674 for quality assurance at Radiography-, Fluoroscopy-, Dentaland Mammography systems. Available with ROA/ ROM / DEDX

Measurement parameter (DEDX):

- Dose: 20 μGy 9,999 mGy
- \bullet Dose rate: 20 μ Gy/s 400 mGy/s
- Time: 1 ms 9,999 s

IBA Dosimetry · DVT-3D



Highlights

without artifacts.

Test of 3D image quality of "Digital Volume Tomography" (DVT) systems, according DIN 6868-150 / DIN 6868-4. Optional Carbon adapter for easy and precise positioning in the beam

Spatial parameter:

- Detail resolution
- Uniformity and noise
- Laser marks for convenient positioning in iso-center

IBA Dosimetry · Spot-Luminance meter LXcan



Highlights

For luminance measurements at image display devices according DIN 6868-157, DIN V 6868-57, IEC 61223-2-5 and AAPM TG18.

- Distance and contact measurement
- Easy targeting with a built-in camera and display



Optional photometric detector LX-LS to measure the Illuminace in combination with I Xcan

IBA Dosimetry · KermaX plus SDP



Highlights

Easy to install standard dosimeter dedicated to measure DAP and DAP rate for patient dose monitoring. Rectangular, transparent ionization chamber and separate 10-digit background lighting LCD "Single Line Display".

Measurement parameter:

- DAP rate:
- $0.01 \, \mu \text{Gym}^2/\text{s} 3,000 \, \mu \text{Gym}^2/\text{s}$
- DAP resolution: 0.01 μGym²
- Interface: 1 × RS232 (RIS/HIS or printer)

IBA Dosimetry · KermaX plus IDP



Highlights

Ideal solution for a quick and convenient retrofit installation to measure DAP and DAP rate for patient dose monitoring.

Rectangular, transparent ionization chamber with integrated 10-digit internal background lighting LCD display.

Measurement parameter:

- DAP rate:
- $0.01 \mu \text{Gym}^2/\text{s} 3,000 \mu \text{Gym}^2/\text{s}$
- DAP resolution: 0.01 μGym²
- Interface (optional): RS232, RS485

IBA Dosimetry · KermaX plus DDP "Duo"



Highlights

Multifunctional duo-channel dosimeter dedicated to measure DAP, DAP rate and exposure time in patient dose monitoring.

Two Rectangular, transparent ionization chamber with integrated electronics and one separate "Dual Line Display" with two very bright LED display lines.

Measurement parameter:

- DAP rate: 0.01 μGym²/s 3,000 μGym²/s
- DAP resolution: 0.01 μGym²
- Interface: $2 \times RS$ 232 (RIS/HIS and printer)

IBA Dosimetry · KermaX plus TinO IDP



Highlights

Two in One – Dose Area Product and dose measurements in one Chamber. Rectangular, transparent ionization chamber with integrated 10-digit internal background lighting LCD display for easy and smart installation at collimator rails.

Measurement parameter:

- DAP rate:
- $0.01 \, \mu \text{Gym}^2/\text{s} 3,000 \, \mu \text{Gym}^2/\text{s}$
- DAP resolution: 0.01 μGym²
- Interface (optional): RS232, RS485, CAN

Quart · dent/digitest Dental QA/QC Test Phantom



Highlights

- Quart dent/digitest 2D dental test phantoms are designed to assess x-ray imaging parameters according DIN and IEC QA/QC requirements.
- Features patient equivalent filtration and test objects to perform full-scale x-ray image quality analyses.

Parameters:

- Spatial resolution
- High-contrast resolution
- Low-contrast resolution
- Homogeneity/artefacts
- Radiation field/tube alignment

Quart · DVT AP Cone-Beam CT Test Phantom



Highlights

- The Quart DVT AP phantom is designed for QA/QC at conebeam CT (CBCT), dental volume tomography (DVT) and 3D imaging equipment.
- It is to be used in dental 3D imaging (according DIN 6868-161 requirements) as well as angiography in C-arm x-ray applications (manufacturer-specific applications). Based on latest research, the solution can also be utilised for standard CT IQ tests.

Quart · DVT 150 CBCT IQ Test Phantom



Highlights

- The Quart DVT 150 phantom is designed to meet the requirements of the German DIN 6868-150 x-ray imaging acceptance test standard.
- Handling and positioning of the phantom is easy and straight-forward.
 It enables quick and simple contrast resolution tests for 3D, ENT and angiography x-ray applications.

Quart · mam/digi Mammography IQ Phantom



Highlights

- The Quart mam/digi phantom is designed to be used as universal tool for QA/QC routine testing in digital and analog mammography. The phantom creates a link between technical and clinical image quality. It can also be used as QA tool for digital tomosynthesis.
- The phantom incorporates Quart's unique Landolt ring objects. They serve to verify low-contrast and perceptibility limits.

Quart · SP dl R/F IQ Phantom



Highlights

- The Quart SP dl phantom enables assessment of digital x-ray equipment according to the german DIN 6868-150 and DIN 6868-4.
- The phantom is available with a unique kV test object to assess radiation quality and generator performance on a routinely basis.
- For ease of use, a frame / extension is provided as well as a wire-mount system for use with wall stand units.

Quart · Anthropomorphic X-Ray Phantoms

Highlights

- Our german-made anthropomorphic phantoms allow repeated x-ray imaging of specific body regions. They are used in x-ray trainings or for specific equipment tests under life-like conditions.
- The phantoms comprise of real human bones embedded in tissue-equivalent material.

Available phantom versions

- Full Body
- Head
- Hand / arm
- Hip/spine
- Foot/leg
- Special training phantoms



Quart · didoNEO R Diagnostic X-Ray Dosemeter



Highlights

The Quart didoNEO introduces a new approach to diagnostic x-ray meters: it features the most compact base unit and most compact detector in the x-ray meter industry. The didoNEO R is used for

QA and service in radiography, (Pulsed) fluoroscopy, DSA, dental, 3D (CBCT).

- Compact multi-functional state-of-the-art solid state detector
- Enables measurements in spots with limited space
- Measures behind scatter radiation grids
- Direct measurement of DLP/DWP in dental OPG

Quart · didoCT Pencil Chamber Meter



Highlights

The Quart didoCT pencil-shaped ion chamber meter is designed for easy and precise dose-width product measurements.

- \bullet The meter does not require any pre-setting procedure for direct reading of DWP, rate and time.
- As an optional feature, the Quart didoCT can be supplied with free-in-air direct HVL measurement capability. This device feature is unique and had only been introduced by Quart in a CTDI chamber.

Quart · didoEASY Diagnostic X-Ray Meters



Highlights

- The Quart didoEASY meters are designed for quick measurements of dose, dose rate and exposure time in X-ray QA/QC and service.
- didoEASY meters automatically compensate all radiation qualities in their area of application. Three meter versions are available: for R/F and dental (50 – 150 kV), for mammography (25 – 40 kV), and one for the full diagnostic range (25 – 150 kV).

Quart · didoSVM Precision Survey Meter



Highlights

- The Quart didoSVM medical survey meter is designed to detect beta, gamma and x-ray sources of very low intensity around diagnostic x-ray equipment as well as in radiation therapy environments. Excellent energy response to measure radiation rate and dose.
- Its detection technology is based on solid-state components, enabling measurements with high sensitivity and very quick response.

Quart · nonius Digital X-Ray Ruler



Highlights

- The Quart nonius is a sophisticated, fully electronic x-ray ruler to verify size and geometrical properties of x-ray fields in radiography and mammography. It can also be used to analyse fanned CT or dental OPG x-ray beams.
- Its resolution capabilities and precision go down into to the nonius range of 0.1 mm!
- ullet Take only 3 steps to obtain the test result: Position Expose Evaluate.

Quart · DSA Test Phantom



Highlights

- The Quart DSA phantom features longitudinal sliding technique to minimise structural movement artefacts in the test image. It complies with DIN 6868-4, 6868-150 and IEC 61223-3-3.
- A special characteristic of the phantom is that it realistically reproduces the injection procedure of the contrast agent into vessels with different attenuation properties.

Testing Devices

Radcal · Accu-Gold Touch Professional



Highlights

- x-ray quality assurance instrument
- Includes all of the capabilities of the Accu-Gold Touch system
- Supports the full range of x-ray modalities incorporating the most extensive sensor selection in the
- Rechargeable Battery provides hours of continuous use
- · Embedded memory storage
- · The Touch Pro includes wired and wireless (wifi) computer interfaces
- · A computer interface allows automated report generation - when your measurements are done, your report is done
- Powerful waveform analysis tools provide a deeper understanding of the x-ray system

Radcal · Accu-Gold Touch



Highlights

- · A comprehensive stand alone diagnostic test instrument
- Supports the full range of x-ray modalities incorporating the most extensive sensor selection in the
- · Access to Radcal's full line of solidstate, ion chamber, and mA sensors
- Includes a large easy-to-read 5" touchscreen display
- · Simple, accurate, and reliable measurement of Dose, Dose Rate, kV, HVL, Filtration, Pulse parameters, and mA
- Provides the ability to simultaneously measure with solid-state and ion chamber sensor
- Rechargeable Lithium Ion Battery provides hours of continuous use
- · Embedded memory stores all measurement data for over 10,000 exposures

Radcal · Accu-Gold+

Highlights

- •The most versatile dose measure ment instrument available
- · An excellent solution for Radiography, Fluoroscopy, Mammography, CT, and Dental applications
- Supports Radcal's full line of ion chambers and solid-state dose sensors
- Complimentary Accu-Gold software incorporates an easy-to-use interface customizable for your individual needs
- A streamlined Excel interface allows automated report generation -
- when your measurements are done, vour report is done.
- Powerful waveform analysis tools provide a deeper understanding of the x-ray system
- Secure WiFi connectivity is available with the (optional) Nugget adapter

Radcal · Accu-Dose+



Highlights

- · The most versatile dose measurement instrument available
- · An excellent solution for all X-ray modalities
- Supports Radcal's full line of ion chambers and solid-state close
- Complimentary Accu-Gold software incorporates an interface customizable for your individual needs
- · A computer interface allows automated report generation
- Powerful waveform analysis tools provide a deeper understanding of the x-ray system
- Secure WiFi connectivity is available with the (optional) Nugget adapter

RTI Group · Ocean software



Highlights

Minimise your working time with the market-leading software for X-ray QA. Versatile, traceable and swift, Ocean provides simple-to-use features for routine inspection – using Quick Check – through to highly customisable QA programs for deep analysis and record keeping in Ocean professional. Ocean is central to RTI's X-ray solutions and is compliant with all our product lines.

RTI Group · CT Ion Chamber



Launched amid much fanfare at RSNA 2019, the new CT Ion Chamber 10 cm and CT Ion Chamber 30 cm are both pencil type ion chambers intended for measuring the exposure output level of CT scanners in a CTDI Phantom or freein-air. The Ion Chambers are compatible with the RTI Chamber Adapter for use with the Piranha and Cobia. They can also be used with the older Barracuda and Solidose 400

Testing Devices

RTI Group · Cobia



Highlights

The RTI Cobia family offers a straightforward instrument that everyone who needs to make reliable quality control can use. All Cobias are tailored for efficiently measuring a variety of radiography and fluoroscopy parameters. All four models are wireless, have bluetooth connection and are fully compliant with RTI's Ocean software. Run your Cobia in English or choose from several languages.


RTI Group · Piranha



Highlights

The Black Piranha is RTI's premium platform for advanced, reliable QA. The Piranha Multi can be used for X-ray QA of all modalities – R/F, Dental, Mammo and CT – whereas the other four Piranhas are dedicated to one specific modality. All Piranha models are compliant with RTI's Ocean software and come with bluetooth connection. Automatic connection to various RTI accessories - just plug 'n' play.

VacuTec · **VacuTec AEC Chamber**

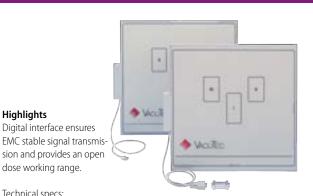

Technical specs:

dose working range.

Highlights

- Tube voltage: 40 kV ...150 kV
- Dose rate range: 0.5 . . . 1,000 μ Gy/s
- Aluminum equivalent: <0.75 mm Al
- · Analog interface: ramp voltage 0 – 10 V
- Digital interface: differential pulses (RS422)
- Resolution: 0.025 μGy
- Pulse width: 2 μs

VacuTec · VacuDAP-C/VacuDAP-C duo



Highlights

The VacuDAP-C systems for measurement of DAP and Dose are basically integrated in interventional devices with customized calibration settings.

Technical specs:

- Resolution DAP: 0.01 µGym²
- Resolution Dose: 0.005 mGy
- Interface: RS485, RS232, Bluetooth, CAN
- Active area: Ø (8 . . . 100) mm

VacuTec · VacuDAP Bluetooth



Highlights

- VacuDAP chamber is now available with Bluetooth technology.
- Perfect suitable for DR upgrades and mobile X-ray units.
- The battery ensures simplest installation ever.

Technical specs:

- Resolution DAP: 0.01 µGym²
- · Active area:
- 123×123 mm/147×147 mm
- Battery operation time: about 24 h

VacuTec · VacuDAP/VacuDAP duo



Highlights

The VacuDAP family provides a wide range of DAP and Dose measuring solutions for most of the diagnostic X-ray systems in the market.

Technical specs:

- Resolution DAP: 0.01 µGym²
- Resolution Dose: 0.003 mGy
- Interface: RS485, RS232, Bluetooth, CAN
- · Active area:

123×123 mm/147×147 mm

146

		■ Computed Tomography	■ MRI	Injectors	Interventional Systems	Artificial Intelligence	■ IT Systems	■ Mammography	R/F Film-Screen	R/F Digital	Molecular Imaging	■ Displays	■ Printers	■ Ultrasound	Testing Devices
AB-CT – Advanced Breast-CT GmbH Henkestrasse 91 91052 Erlangen, Germany Tel +49 9131 97310 0 ask.crm@ab-ct.com www.ab-ct.com	АВ-СТ	15						67							
Agfa HealthCare Septestraat 27 2640 Mortsel, Belgium tel +32 3 444 94 44 agfahealthcareinfo.be@agfa.com www.agfa.com	AGFA 490 Health Care						53 54 55 59 62 63			82 83 96 97 101 102 107			121 122		
AlgoMedica Europe GmbH Waldhofer Strasse 102 69123 Heidelberg, Germany tel +49 6221 321 8715 info@algomedica.eu www.algomedica.com	& Algo Medica					50	55								
allMRI GmbH Südstrasse 23 74226 Nordheim, Germany tel +49 71 33 2 37 02 20 mail@allmri.com www.allmri.com	alMR		30 31												
BMS Informationstechnologie GmbH Diesterweggasse 7/1 1140 Vienna, Austria tel +43 1 524 81 34 00 info@bms-austria.com www.easydose.eu	BMS INFORMATIONSTECHNOLOGIE® GMBH						55 63 64								
Bracco Injeneering S.A. Avenue de Sévelin 46 1004 Lausanne, Switzerland tel +41 21 621 74 00 info:injeneering@bracco.com www.imaging.bracco.com	BRACCO INJENEERING			33 36			55 64								
Canon Medical Systems Europe B.V. Zilverstraat 1 2718 RP Zoetermeer, The Netherlands tel +31 79 368 92 22 eu.medical.canon	Canon	7 8 12 15	25 26		38 39 40	50	54 55 60 61	67	73 75	84 85 107 108 112	116			124 125 126	
Canon Europe NV Medical Components Business Group Bovenkerkerweg 59 1185 XB Anstelveen, The Netherlands tel +31 205 45 89 26 medical.drsales@canon-europe.com www.canon-europe.com/medical	Canon						65			84 89 90 92 98 99 105 109					
Canon Electron Tubes & Devices Co., Ltd 1385 Shimoishigami Otawara-Shi, Tochigi 324-8550, Japan tel +81 287 26 66 66 https://etd.canon/eng	CANON ELECTRON TUBES & DEVICES CO., LTD.	16			48					97 98 112					
Cefla s.c. Via Selice Provinciale 23A 40026, Imola (BO), Italy tel +39 045 820 27 27 info@newtom.it www.newtom.it	Newton what's next	15													
DEL MEDICAL 28 Calvert Street, Harrison, NY 10528, USA Tel: +1800 261-9808 241 Covington Drive, Bloomingdale, IL 60108, USA Tel: +1800 800-6006 www.delmedical.com	DEL MEDICAL				44		53 54 55			85 86 99 102 108					

		■ Computed Tomography	■ MRI	■ Injectors	■ Interventional Systems	Artificial Intelligence	■ IT Systems	Mammography	R/F Film-Screen	R/F Digital	Molecular Imaging	■ Displays	Printers	■ Ultrasound	Testing Devices
DK Medical Systems Co., Ltd 18. Baumoe-ro 7gil, Seocho-gu, Seoul, 06762, Korea tel +82-2 529 6190 global@dk.co.kr www.dk.co.kr	DK *Medical Systems				45					86 87					
DMS Imaging 393 rue Charles Lindbergh 34130 Mauguio, France tel +33 467 50 49 00 www.dms-imaging.com	/ T DMS L / IMAGING							68		96 111					
DRGEM Corporation 7F, E-Bt/D Gwangmyeong Techno-Park, 60 Haan-ro, Gwangmyeong-si, Gyeonggi-do, Korea tel +82 2 869 85 66 drgem@drgem.co.kr www.drgem.co.kr	DRGEW Your Best Healthcare								75	87 88 99 102 112					
Philips Medical Systems DMC GmbH Röntgenstrasse 24 22335 Hamburg, Germany marketing.dunlee@philips.com www.dunlee.com	DUNLEE	17	29												
EXAMION GmbH Erich-Herion-Strasse 37 70736 Fellbach, Germany tel +49 711 12 00 02-0 info@examion.com www.examion.com	EXAMION * X-Ray Systems · Digital Imaging · Service						53 54 55	70	73 76	82 88 96 99 102 114					
Febromed GmbH & Co. KG Am Landhagen 52 59302 Oelde, Germany tel +49 25 22 30-532 fax +49 25 22 30-8489 info@febromed.de www.febromed.com	febromed	21	31												
FUJIFILM Medical Systems Europe Heesenstrasse 31 40549 Dusseldorf, Germany tel + 49 211 5 08 90 www.fujifilm.eu	FUJ <mark>i</mark> FILM	12 13				50	53 54 55 56 60 61 64	67		88 99 102					
GCTechnology GmbH Freidling 12 84172 Buch am Erlbach, Germany tel +49 87 06 94 15 00 info@ctech-gmbh.com www.gctech-gmbh.com	GCTechnology GmbH	21	31		48			70	78					138	
GMM GROUP Via Partigiani, 25 24068 Seriate (BG), Italy tel +39 035 452 53 11 info@gmmspa.com www.gmmspa.com	GMM				42 45				73	88 89 103 108 109					
Guerbet BP 57400 95943 Roissy CdG Cedex, France tel +33 145 91 50 00 LF@guerbet.com www.guerbet.com	Guerbet ##			33 34 36			55 64								
Hitachi Medical Systems Europe (Holding) AG Sumpfstrasse 13 6300 Zug, Switzerland tel 441 41748 63 33 welcome@hitachi-medical-systems.com www.hitachi-medical-systems.com	HITACHI Inspire the Next	8 13	26 28 29											126 127	

		■ Computed Tomography	■ MRI	Injectors	■ Interventional Systems	Artificial Intelligence	■ IT Systems	■ Mammography	R/F Film-Screen	R/F Digital	Molecular Imaging	■ Displays	Printers	■ Ultrasound	Testing Devices
Hologic bvba Da Vincilaan 5, Building Caprese 1930 Zaventem, Belgium tel +32 2711 46 80 be-info@hologic.com www.hologic.com	HOLOGIC						60								
LA.E. S.P.A. Via Fabio Filzi, 53 20032 Cormano (MI), Italy tel +39 02 66 30 32 55 iaexray@iae.it www.iae.it	iae))	21			48			71	78	114					
IBA Dosimetry GmbH Bahnhofstrasse 5 90592 Schwarzenbruck, Germany tel +49 91 28 607-14 dosimetry-info@iba-group.com www.iba-dosimetry.com	iba Dosimetry														140 141 142
IMAGE Information Systems Europe GmbH Lange Strasse 16 18055 Rostock, Germany tel +49 381 496 58 20 info@image-systems.biz www.image-systems.biz	IMAGE Information Systems						54 55 56 62 63					120			
IMS Giotto S.p.A. – GMM GROUP – Via Sagittario, 5 40037 Sasso Marconi (BO), Italy tel +3951 84 68 51 imscomm@imsgiotto.com www.imsgiotto.com	IMS Siotto							67 68 70							
INTERMEDICAL SRL Via E. Fermi, 26 24050 Grassobbio (BG), Italy tel +39 035 659 48 11 info@inter-med.it www.inter-med.it	INTERMEDICAL				42 45			67 69		89 103					
i-SOLUTIONS Health GmbH Am Exerzierplatz 14 68167 Mannheim, Germany tel +49 621 39 28-0 info@i-solutions.de www.i-solutions.de	-SOLUTIONS						54 55 56 60 64 65								
ITZ Medicom GmbH & Co. KG Siemensring 44 a 47877 Willich, Germany tel +49 21 54 49 79 60 info@itz-medi.com www.itz-medi.com	itz-medi.com						54 55 56 57 63								
JVCKENWOOD Deutschland GmbH Konrad-Adenauer-Allee 1-11 61118 Bad Vilbel, Germany tel +49 21 61 69 84-180 medical-display.e@jvckenwood.com healthcare.jvc.com	JVC											119 120			
Konica Minolta Business Solutions Europe GmbH Healthcare Division Capellalaan 65 2132 Jl. Hoofddorp, The Netherlands tel + 31 20 658 41 00 healthcare@konicaminolta.eu www.konicaminolta.eu/healthcare	KONICA MINOLTA									83 89 99 100 103			122	127 128	
LEONI Special Cables GmbH Business Unit Healthcare Eschstrasse 1 26169 Friesoythe, Germany tel +49 44 91 291-5040 healthcare@leoni.com www.leoni-healthcare.com	LEONI	22	31					71							

		■ Computed Tomography	■ MRI	Injectors	■ Interventional Systems	Artificial Intelligence	☐ IT Systems	Mammography	R/F Film-Screen	R/F Digital	Molecular Imaging	■ Displays	■ Printers	■ Ultrasound	■ Testing Devices
mediaire GmbH Möckernstrasse 63 10965 Berlin, Germany te +49 30 28 64 90 67 info@mediaire.de www.mediaire.de	medi ai re					50	55								
medigration GmbH DrRudolf-Eberle-Strasse 8-10 76534 Baden-Baden, Germany tel +49 72 23 966 98 60 info@medigration.de www.medigration.de	bender gruppe medigration						54 55 57 60 61 63 65			100		121			
MEDTRON AG Hauptstrasse 255 66128 Saarbruecken, Germany tel +49 681 970 17-0 info@medtron.com www.medtron.com	MED (TRON° AG			34 36											
Medtronic International Trading Sàrl Route du Molliau 31 1131 Tolochenaz, Switzerland tel +41 21 802 70 00 www.oarm.com www.medtronic.com	Medtronic Further, Together				42										
SHENZHEN MINDRAY BIO-MEDICAL ELECTRONICS CO., LTD. Mindray Building, Keji 12th Road South Nanshan, Shenzhen 518057, China tel +86 755 81 88 89 98 intl-market@mindray.com www.mindray.com	mindray		29			51	55			90 103				128 129 130	
NEXUS/CHILI GmbH Friedrich-Ebert-Strasse 2 69221 Dossenheim / Heidelberg, Germany tel +49 62 21 180 79 10 sales@nexus-chili.com www.nexus-chili.com	nexus chili Imaging & radiology solutions						54 55 57 58 61 63					121			
NORAS MRI products GmbH Leibnizstrasse 4 97204 Hoechberg, Germany tel +49 931 29 92 70 mri@noras.de www.noras.de	NORAS MRI products		30												
NRT X-RAY A/S Birkegaardsvej 16 8361 Hassealager, Denmark tel +45 86 28 35 00 nrt@nrtxray.com www.nrtxray.com	(B)NRT									90 109					
Planmed Oy Sorvaajankatu 7 00880 Helsinki, Finland tel +358 20 779 53 00 sales@planmed.com www.planmed.com	Planmed	15						67 69 70							
Primax International SrI Via A.Volta, 10 24060 Torre De' Roveri, Italy tel +39 035 450 00 02 sales@primaxint.com www.primaxint.com	Primak				42					90 103 104					
PROTEC GmbH & Co. KG In den Dorfwiesen 14 71720 Oberstenfeld, Germany tel +49 70 62 925 50 protec@protec-med.com www.protec-med.com	PROTEG TEAM SPERT ABILITY						54 55 58		73 78	90 91 100 101 104					

		■ Computed Tomography	■ MRI	Injectors	Interventional Systems	Artificial Intelligence	■ IT Systems	Mammography	R/F Film-Screen	R/F Digital	Molecular Imaging	■ Displays	Printers	■ Ultrasound	Testing Devices
PTW Freiburg Physikalisch-Technische Werkstätten Dr. Pychlau GmbH Lörracher Strasse 7 79115 Freiburg, Germany tel +49 761 490 55-0 info@ptwdosimetry.com ptwdosimetry.com	THE DOSIMETRY COMPANY	22						71		114					
QUART GmbH Kirchenweg 7 85604 Zorneding, Germany tel +49 81 06 24 91 18 info@quart.biz www.quart.de	QUART Guarry Attendages Technologies														143 144
Radcal Corporation 426 West Duarte Road Monrovia, CA 91016, USA tel +1 626 357 79 21 sales@radcal.com www.radcal.com	Radcal														145
RTI Group Flöjelbergsgatan 8C 43137 Mölndal, Sweden tel +46 31 746 36 27 sales@rtigroup.com www.rtigroup.com	R I INDEPENDENT X-BAY QUALITY ASSURANCE														145 145
Samsung Healthcare 2-dong Alpharium Tower, 145, Pangyoyeok-ro, Bundang-gu, Seongnam-si, Gyeonggi-do, 13530, Republic of Korea tel +82 2 2194 1000 www.samsunghealthcare.com	SAMSUNG													130 131 132	
ScreenPoint Medical Toernooiveld 300 6525 EC Nijmegen, The Netherlands tel +31 242 020 020 info@screenpointmed.com www.screenpoint-medical.com	SCREENPOINT Medical					51	55								
Shimadzu Europa GmbH Medical Systems Division Albert-Hahn-Strasse 6 – 10 47269 Duisburg, Germany tel +49 203 76 87-0 info@shimadzu.eu www.shimadzu.eu	SHIMADZU Excellence in Science				39 40 41 45				73 74 75 76	91 104 109					
Siemens Healthineers Headquarters Siemens Healthcare GmbH Henkestrasse 127 91052 Erlangen, Germany tel +49 800 188 188 5 siemens.com/healthineers	SIEMENS Healthineers	7 8 12 13 16	25 26 27 28		38 39 41 42 43 45 46	51	54 55 58 59 61 62 65	68 69 70	74 76	92 104 110	116 117			132 133 134	
#77, Jinsha Road 515041 Guangdong, China tel +86 754 88 25 01 50 siui@siui.com www.siui.com	See the future					51	55			96				134 135 136 137 138	
Solutions for tomorrow Saxagårdsvägen 5 36251 Väckelsång, Sweden tel +46 10 456 4500 info@solutionsfortomorrow.se www.solutionsfortomorrowse	Solutions for tomorrow									105					
STEPHANIX 10, Rue Jean Moulin 42150 La Ricamarie, France tel +33 4 77 47 81 60 contact@stephanix.com www.stephanix.com	- STEPHANIX MEDICAL IMAGING SOLUTIONS				43 46				74 76	92 93 94 101 105 110					

		■ Computed Tomography	■ MRI	Injectors	Interventional Systems	Artificial Intelligence	■ IT Systems	Mammography	R/F Film-Screen	R/F Digital	Molecular Imaging	■ Displays	Printers	■ Ultrasound	■ Testing Devices
SternMed GmbH Schubertstrasse 31 88214 Ravensburg, Germany tel +49 751 35 97 80 email@sternmed.de www.sternmed.de	STERNMED © We make it possible	13	28 29		43 46			69	77	105					
SuperSonic Imagine Les Jardins de la Duranne, Bât E & F 510, Rue René Descartes 3857 Aix-en-Provence, France tel +33 442 99 24 24 contactsFR@supersonicimagine.com www.supersonicimagine.com	HOLOGIC' SUPERSONIC Imagine													138	
Technix S.p.A. Via Fermi 45 24050 Grassobbio (BG), Italy tel +39 035 384 66 11 technixd@technix.it www.technix.it	TECHNIX				43					105 106					
THERAPIXEL Le Village By CA, Rue Claude Daunesse 06560 Valbonne, France Phone: +33 9 72 55 20 39 contact @ therapixel.com www.therapixel.com	THERAPIXEL					51	55								
Ultrasound Technologies LTD Lodge Way, Portskewett, Caldicot, South Wales, NP26 5PS, U.K. tel +44 12 91 42 54 25 ultratec@doppler.co.uk www.doppler.co.uk	ultrasound technologies	22													
VacuTec Meßtechnik GmbH Dornbluethstrasse 14a 01277 Dresden, Germany tel +49 351 317 24-0 info@vacutec-gmbh.de www.vacutec-gmbh.de	♦ VacuTec														
Varex Imaging Deutschland AG Otto-Brenner-Strass 10 47877 Willich, Germany tel +49 21 54 92 49 80 info@vareximaging.com www.vareximaging.com	WYAREX	22						71							
VILLA SISTEMI MEDICALI s.p.a. Via delle Azalee, 3 20090 Buccinasco (MI), Italy tel +39 02 48 85 91 sales@villasm.com www.villasm.com	VILLA				46			68 69 70	74 75 77 78	94 96 101 106 107 111					
VISARIS d.o.o. Batajnički drum 10., deo 1B 11 080 Zemun, Serbia tel +381 11 201 76 00 · fax +381 11 201 76 70 info@visaris.com https://visaris.com	VISARIS						54 55 59			94 95 96 101 107 111					
Ziehm Imaging GmbH Lina-Ammon-Strasse 10 90471 Nuremberg, Germany tel +49 911 21 72 - 0 info@ziehm.com www.ziehm.com	ziehmimaging				44 47 48										

Note

Notes	

SOMATOM X.cite with myExam Companion

Intelligent imaging. Excellence empowered.







The number and complexity of radiological procedures are increasing. This calls for new solutions: Intelligent navigation. Patient-friendly design. Personalized imaging. Consistent standards.

Our single source CT scanner SOMATOM® X.cite offers exactly that. Together with myExam Companion it launches the era of intelligent imaging. myExam Companion guides users intuitively through any procedure, adjusting key parameters to individual patients. This allows technologists of any skill level to unlock the full potential of this new

high-end system, including its Vectron $^{\text{TM}}$ X-ray tube and large, 82 cm bore.

No matter the patient or complexity of the procedure, SOMATOM X.cite and myExam Companion together enables users to reduce unwarranted variations – and generate consistent, comprehensive results.

SOMATOM X.cite empowers excellence in computed tomography.

siemens-healthineers.com/somatom-xcite







Deep Intelligence Innovate. Illuminate. Initiate.

Introducing the all-new Aquilion ONE / PRISM Edition.

Combining the power of Deep Learning Spectral Reconstruction with our Advanced intelligent Clear-IQ Engine (AiCE), we have created a CT system that is designed for deep intelligence.

By seamlessly integrating AI technologies, you can now access all the tools you need to automate your workflows, while capturing vital diagnostic information that will help define your patient's medical care journey.

With Aquilion ONE / PRISM Edition, the possibilities are endless – for you, your patients and your business. It's time to innovate with technology, illuminate clinical insights and initiate new business opportunities for a brighter, smarter future.

https://eu.medical.canon/prism